

Instruction Manual

Please do not apply any skincare products before using this product.

Please follow the instructions provided in this manual for product installation and usage.



This product is intended for use in a home care environment, regardless of whether it is used by inexperienced operators or trained medical personnel.

Thank you for purchasing this product. To facilitate its use, please read and keep this instruction manual before use. The images in this manual are for illustration purposes only, and the actual product may vary.

Contents

1. Package Contents	01	6. Cleaning	09
2. Safety Guide	02	7. Storage	10
2.1 Proper Use	02	8. Troubleshooting	10
2.2 Safety Instructions	02	9. Explanation of Symbols	11
3. Technical Data	04	10. Disposal	12
4. Charging Instructions	04	11. Warranty	12
5. How to Use	06		
5.1 Sensitivity Test	06		
5.2 Operating Instructions	06		
5.3 Mode Selection	08		
5.4 Skin Treatment Regimes	09		

1. Package Contents

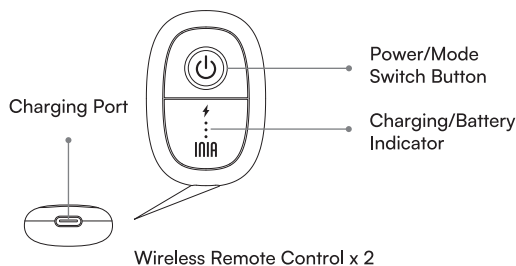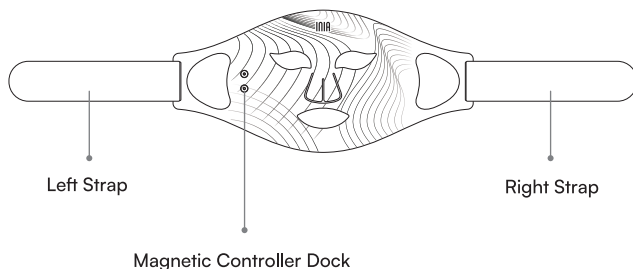

Protective Eye Pads



USB Charging Cable

2. Safety Guide

2.1 Proper Use

This product is specifically engineered to rejuvenate your skin, enhancing its radiance, elasticity, and texture. It is also effective in reducing the visibility of fine lines, wrinkles, and other skin imperfections.

Prior to using the product, we strongly recommend thoroughly reading the user manual, paying close attention to the safety instructions. Neglecting to do so could potentially lead to personal injury or damage to the product.

2.2 Safety Instructions



Contraindications

- **Pregnancy and Nursing:** Do not use the mask if you are pregnant, planning to become pregnant, or breastfeeding.
- **Skin Conditions:** Avoid using the mask if you have conditions such as Lupus erythematosus, photosensitive eczema, or Albinism.
- **Photosensitivity:** Refrain from using the mask if you suffer from any photosensitive disorder (sensitivity to light) or are taking medications that can cause photosensitivity.
- **Skin Health:** This product should not be used on skin with existing diseases, wound ulceration, trauma, or post-operative wounds.
- **Medical Conditions:** The product is not suitable for individuals with epilepsy, thyroid disorders, or those using thyroid medications.



Warnings

- **Preparation:** Prior to treatment, ensure that you have removed creams, lotions, or makeup from your skin. Using the mask on clean,

makeup-free skin is essential.

- **Start the Treatment:** It is recommended to first put on the face mask, and then magnetically attach the controller to the mask to power on.
- **Discomfort and Skin Reactions:** Do not use the mask if you experience any discomfort, abnormal reactions (such as itching or uncomfortable inflammation), or if your skin becomes excessively warm to the point of discomfort. Stop use immediately in such cases.
- **Avoid Direct Staring:** Never look directly at the light source part of the light therapy mask to protect your eyes from potential harm.
- **Session Duration:** Limit each session to a maximum of 20 minutes. Prolonged use beyond this duration may result in skin damage.
- **Keep Dry:** Do not immerse the product or power cable in water or any other liquids. Keep them away from moisture at all times.
- **Heat and Flames:** Avoid exposing the product to open flames, hot surfaces, or high temperatures, as this may damage the device.
- **No Modifications:** Do not attempt to modify or replace the power cable and battery by yourself. Seek professional assistance for any servicing or replacements.
- **Product Care:** Avoid folding or rolling up the product and refrain from stacking heavy objects on top of the face mask to prevent damage.



Dangers

- **Product Condition:** Do not use the product or power cable if they are damaged or defective. Using these components in a damaged state can be dangerous.
- **Charging Caution:** Avoid touching any part of the product with wet hands during the charging process to prevent electric shock.
- **Charging Specifics:** The included power cable is exclusively for charging the internal lithium battery within the controller and must not be used for any other purposes.
- **Unplugging:** Always remember to unplug the USB wall adapter when


you are not charging the product. This is essential for safety and energy conservation.

3. Technical Data


LED Color and Wavelength	Red Light: 620~630nm Blue Light: 460~470nm Orange Light: 590~610nm Near-infrared Light: 850nm
Number of LED Lights	68 LEDs (4-in-1 chips, 272 chips)
Product Dimensions	411 x 195 x 4mm
Product Weight	160g
Input Voltage	5V
Battery Capacity	1000mAhx2
Charging Time	2 hours
Working Temperature	0 C ~ +40 C
Storage Temperature	-10 C ~ +60 C

4. Charging Instructions


1. When the battery is low, pressing the power button will not turn the device on. All the indicator lights will flash five times with five quick beeps.

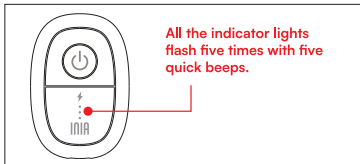

2. Fully charge the controller alone (ensure that the mask is disconnected) by inserting the provided USB-C cable into the controller, then connect the power supply to a power source, such as a computer or an adapter.

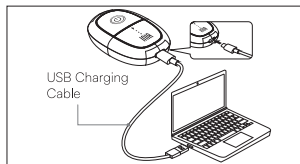

Using a Computer

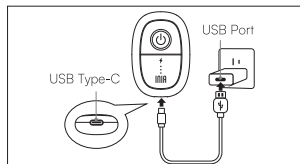

Using a Household Power Outlet

Note: Each controller cannot be used while charging.

3. The indicators on the controller will blink to indicate that charging has begun. Once charging is complete, the three indicators will change from blinking to steady. You can monitor the charging progress and available treatment duration based on the charging indicator:

Indicator Status	Charge Level	Available Treatment Duration
One Indicator Steady White	30% Charged	10 Minutes
Two Indicators Steady White	60% Charged	10-25 Minutes
All Indicators Steady White	Fully Charged	25-40 Minutes

5. How to Use

5.1 Sensitivity Test

Before carrying out your first treatment, you must perform a sensitization test. Follow the step-by-step guide below:

1. Place the light therapy mask on your inner arm.
2. Attach the controller to the mask magnetically to enable the device to complete a 10-minute treatment.
3. If any abnormal symptoms appear, such as redness, itching or allergies appear during use. Stop using it immediately.
4. Wait at least 6 hours after treatment. If the skin appears normal, begin your first treatment.

5.2 Operating Instructions




Scan this code to watch the tutorial video.




STEP 1

Thoroughly clean and dry your face and hands, ensuring they are free from any makeup or dirt.

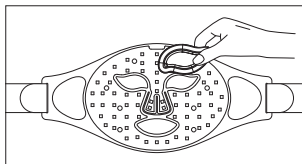

STEP 2

Insert the protective eye pads into the eye-holes of the mask.

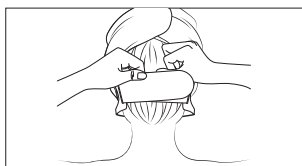


STEP 3

Place the mask comfortably on your face and fasten the adjustable Velcro straps around the back of your head.



STEP 4

Magnetically attach the controller to the mask. The device will automatically turn on in anti-aging mode.

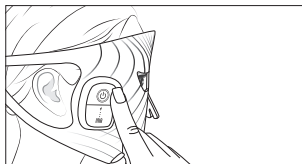

STEP 5


Briefly press the power/mode switch button to switch modes as needed, each accompanied by a quick beep.



STEP 6

Removing the remote during treatment puts the device in standby mode with the controller indicator light breathing. Reattach the remote to activate or press the power button for 1.5 seconds to turn off.

Note: In standby mode, pressing the power button beeps three times, indicating an invalid operation, and the device will shut down after 2 minutes of inactivity.

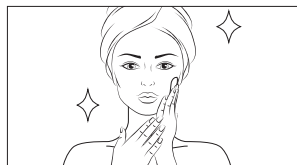

STEP 7

Apply serum and/or day/night cream on your face after treatment.

5.3 Mode Selection

Our light therapy mask features four unique modes, each targeting different skin concerns for a comprehensive solution to healthier, more radiant skin. See below for the functions of each mode.

Mode	LED	Best For	Recommended Treatment Duration
Anti-aging Mode	630nm R + 850nm NIR	Reduces fine lines and wrinkles, stimulates collagen, and helps to improve blood circulation.	10-20 Minutes
Acne Mode	460nm B + 850nm NIR	Reduces acne-related redness and imperfections for clearer skin.	
Brightening Mode	610nm O + 850nm NIR	Suppresses melanin formation, evenly increases the skin brightness, and softens the skin.	
Smoothing Mode	630nm R + 460nm B + 850nm NIR	Combats aging and reduces skin imperfections to enhance skin vitality and radiance.	


5.4 Skin Treatment Regimes

If you are using the product for the first time, follow the schedule for guidance:

- a. Week 1: Use once every 3 days for 10 minutes at a time.
- b. Week 2: Use once every 2 days for 10 minutes at a time.
- c. Week 3: Use once every 2 days for 10-15 minutes at a time.
- d. Week 4: Use once every day for 10-15 minutes at a time.

6. Cleaning

Ensure the face mask is disconnected from the controller. Use a damp cloth or a water-based wipe to wipe off the face mask. Then wipe it dry with a dry cloth.

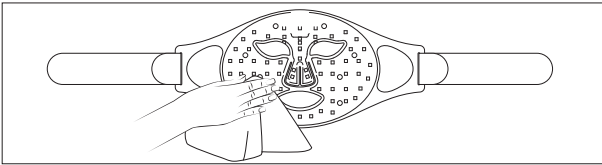

Warning: Do not use any scouring pads, abrasive, or chemical cleaners as these may cause permanent damage to the cosmetic surfaces.



Warning: Do not immerse any of the components in water or any other liquid.

7. Storage

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.

8. Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact INIA team. Do not repair the product yourself.

Problem	Solution
The device can not turn on.	Make sure the wireless remote controller is properly charged.
The product does not charge.	Make sure the USB charging cable is securely connected between the charging port and the power source.
The controller indicator light blinks and emits five quick beeps during use.	The battery is running low. Please charge the controller.

9. Explanation of Symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

Symbol	Explanation of Content
 Contraindications	This signal symbol/word designates a hazard with moderate risk, which may result in severe injury if not avoided.
 Warnings	Dangerous Content Due to Incorrect Operation
 Dangers	Contents Related to Personal Safety Due to Incorrect Operation
	Keep Dry
	Keep Away From Fire
	Disposal instructions for electronic and electrical products. Disposal of the product should be done in accordance with regulations.
	FCC Federal Communications Commission
	Compliant With CE's EMC Standard Requirements
RoHS	Restriction of Hazardous Substances

10. Disposal



The product should be disposed of separately from household waste. Should the product no longer be capable of being used at some point in time, dispose of it in accordance with the regulations in force in your city or county.

11. Warranty

The warranty period for this product may vary depending on the country of sale. We will strive to extend the warranty period to the maximum extent permitted by the laws of that country. If the instrument experiences quality issues within the warranty period, we will provide free repair or replacement. However, damages caused by improper use or non-compliant operations are not covered under warranty. In such cases, we will replace the damaged parts or offer a trade-in at cost price. The ownership of the old instrument belongs to our company.

Please present your order purchase receipt and our official website's after-sales registration information when requesting warranty service. The warranty period for our products starts from the date of receipt of the order at the time of purchase.

Our company may not offer warranty or return services under the following conditions:

-
1. Product damage caused by unforeseeable events such as fire, earthquake, flood, etc.
 2. Instrument damage caused by unauthorized disassembly or improper storage (e.g., rodent infestation, liquid ingress).
 3. Damage caused by the forced use of the product in non-product usage environments, including but not limited to exceeding workloads. Equipment damage caused by abnormal external forces (falling, compression, collision, etc.).
 4. Equipment damage caused by the use of non-manufacturer accessories and expansion components.

Our company shall not be liable for consequential damages or losses arising from your negligence, misuse, operation beyond the instructions, or other conditions beyond INIA's control, including but not limited to loss of business profits, business interruption, other losses caused by the use of the instrument, equipment replacement costs, or customer claims.

Our company will assume the types and extent of liability for unexpected events in accordance with relevant laws and regulations.

*All rights of final interpretation belong to VITALITY NEXUS INC.

FCC Statement

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- *Reorient or relocate the receiving antenna.
- *Increase the separation between the equipment and receiver.
- *Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- *Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.