


# Spectral Mask Instrument

## User Guide

Model 2344

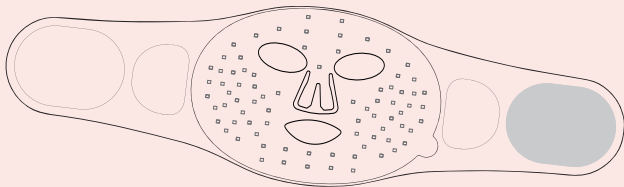

Please be sure to read this manual before using this product.

Thank you for purchasing  
our Infrared Light Beauty Mask.

This mask uses infrared light technology to  
help improve skin health and appearance.

Please read this manual carefully before  
use to ensure safe and effective use.


## Contents

Product Introduction	01
Packing List	02
Product parameters	03
Installation steps	03
Controller operating instructions	05
Introduction Of Each Mode	06
Steps for Use	07
Precautions before Use	08
Precautions for Use	09
After-Sales Service	12

## Product Introduction

This is committed to developing a powerful, portable and easy to use light therapy beauty instrument, beauty equipment does not work most often because the instrument is cumbersome to operate customers can not stick to use, so we developed a medical silicone LED light therapy mask instrument, bendable and portable, 10~15minutes can be a key skin resurfacing. 100\*4 EU CE certification of the hospital grade light source, intensive care of the face, its source of continuous the skin rejuvenation energy can effectively improve the epidermis and dermis of the skin, and restore the ideal skin of women "filter pulling full".

## Packing List

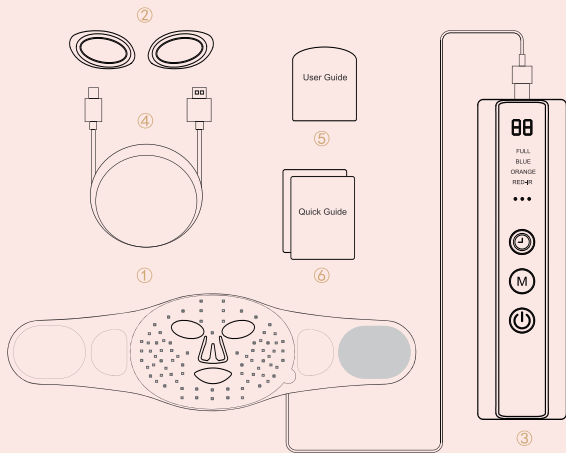

① Spectral Mask Instrument

② Blackout eye mask

③ Controller

④ USB-C Cable

⑤ User Guide

⑥ Quick Guide


## Product Parameters

Product name	Spectral Mask Instrument	Number of Leds	100*4
Material	ABS+PC+silicone	Charging Cable Length	1m
Face Mask Dimension	665X7X193mm	Controller Dimension	321*262*118mm
Charging Voltage	5V	Operating Current	0.85A
Rated power	3.2W	Battery Capacity	2000mAh
Charge time	<2.2H	Perating Time	10 to 30 mins
Lights Wave	BLUE: 455-465nm;ORANGE: 600-610nm; RED:625-635nm;NIR:830-850nm.		

## Installation Steps

### 1.Installation of eye shields


Please install the eye masks into the corresponding positions.  
(divided into left and right)

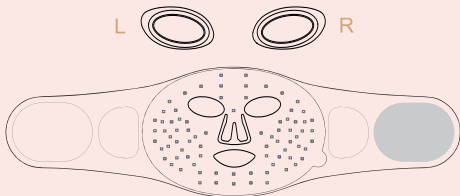


③

## 2. Connect the remote control

Please plug the USB port on the mask into the remote control.

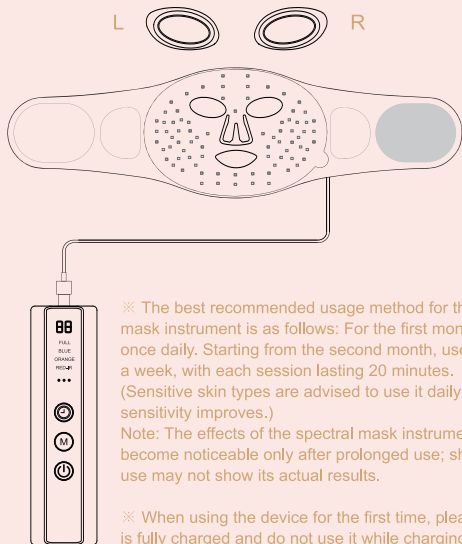

## Controller operating instructions

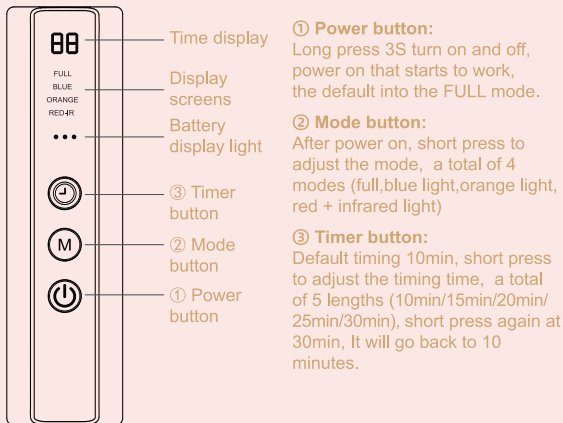

### Battery level indication description:

"•••" represents a battery level of more than 75%.

"••" represents a battery level of more than 45%.

"•" represents a battery level of less than 30%. Flashing indicates a need for charging.

## Introduction Of Each Mode

### • Major Functions

#### •FULL Mode (Blue light,Orange light,Red light+ infrared light)

The three light modes work in turns to meet your daily skincare needs.

#### •BLUE Mode(455-465nm)

Blue light can help improve acne, acne marks, and balance oil secretion.

#### •ORANGE Mode(600-610nm)

Orange light helps improve skin redness, gives you delicate skin.

#### •RED Mode(625-635nm) & infrared (830-850nm) Mode

The combination of red light and infrared light can act on subcutaneous tissue and is the best combination to tighten the skin and increase skin elasticity and luster.


### ● Different Color , Different Function

● Blue light (wavelength 455-465nm)

● Orange light (wavelength 600-610nm)

● Red light (wavelength 625-635nm)

● Infrared light (wavelength 830-850nm)

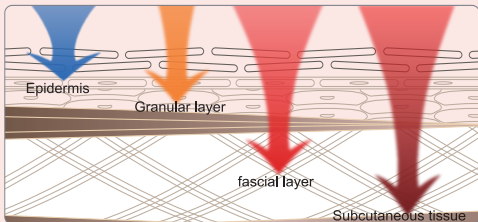

※ NASA research in the 20th century found that different wavelengths of light have different skin care effects on the human body

## Steps for Use

### How to use and steps

**Step 1** Thoroughly cleanse and dry the face;

**Step 2** Install the protective eye shield, Insert the USB interface into the remote control, place the mask on the face, Adjust and fasten the Velcro so that the mask is fixed on the face. Press and hold the power button of the remote control to turn on the device;

**Step 3** The device contains 4 modes according to the mode, including full, blue light, orange light, red light + infrared, which can be adjusted according to your needs. Turn on the default full mode, press the mode button to switch;

**Step 4** The product features a timer function with options for 10 minutes, 15 minutes, 20 minutes, 25 minutes, and 30 minutes. You can adjust the time using the "Timer button". After the timer ends, the device will automatically shut off. Remove the mask to stop the treatment.



1. Clean and dry your face.



2. Install the blindfold .



3. Wearing a face shield and switching on.



4. Freely choose a skin care mode for light therapy.



5. Remove the mask and perform skin care.



6. Wipe the silicone parts with a soft towel or wet wipe.

## Precautions before Use

### 1.Wearing the LED face mask

- Please do not use this product with makeup on.
- This Silicone LED face mask comes with an adjustable velcro strap that allows you to comfortably wear the mask.
- Put the Blackout eye mask into the eye area of the mask to avoid the bright light.

### 2.Sensitization test

Before carrying out your first treatment you must carry out a sensitization test. Follow the step by step guide below:

- 1) Place the Silicone LED face mask over your inner arm.
- 2) Switch on the controller and allow the Silicone LED face mask to run its full 10-minute treatment program.
- 3) At any time if your skin feels HOT or uncomfortable. STOP using the device immediately.
- 4) After the treatment wait at least 6 hours.If the skin on your inner arm appears red and/or itchy Silicone LED face mask is unsuitable for you.

### Prohibited Use

Note: If you fall into any of the following categories, please do not use this product:

1. Individuals who are pregnant, trying to conceive, or breastfeeding.
2. Individuals with photosensitivity or those currently taking any medications that may cause photosensitivity.
3. Individuals with hereditary skin conditions (such as lupus, photosensitive eczema, or albinism).
4. Individuals who have been advised by a doctor not to use this product.

## Precautions for Use

### **Warning**

- Do not use the equipment in a manner other than that described in the manual, as this may cause the machine to malfunction or cause an accident.
- If the user is suffering from physical, sensory or mental illness, or is undergoing treatment for a dermatological condition, the product must be used only under supervision or under the guidance of a responsible person (such as a dermatologist).
- Do not wrap the power cord around the power connector.
- Do not damage the power cord, bend, twist, stretch or place heavy objects on it.
- Do not place the product in a place where the temperature is high, as this may lead to fire or electric shock.
- Do not put the equipment into the fire.
- Do not use it for infants or children and store it out of the reach of children.
- Do not immerse the equipment in water or wash it with water.

## **Attention**

- This product is intended for daily use only; do not disassemble, modify, or repair it without authorization.
- Do not use or store the device and its accessories in high-humidity places such as bathrooms.
- If you experience pain or discomfort during use, stop using the product immediately.
- Do not share this product with others.
- Do not throw the device, as strong impacts may damage the product.
- Do not touch the product with sharp or high-temperature objects.
- When using, wear protective eye mask to ensure a proper fit and prevent light leakage.
- Do not stare directly at the LED lights for an extended period, as it may cause eye damage.
- Strictly adhere to the recommended usage time.
- If the product is damaged, stop using it immediately.

## Warm reminder:

LED light will affect skin cells, and the skin will be slightly heated and dry, which is normal.

Using moisturizing masks and skin care products after light therapy can relieve it.

The USB port of the mask needs to be plugged into the remote control for the mask to function properly. The remote control acts as a power bank. The mode and timing functions of the mask can be adjusted via the remote control. If the mask is connected directly to a charger, it can only maintain a single mode and cannot change light colors or other functions.

When the remote control is out of battery , the battery indicator light will flash.

The remote control can only provide power to the face mask. Please do not use the remote control to charge your phone or other electronic devices

## After-Sales Service

Thank you for choosing our Infrared Light Beauty Mask. We hope it enhances your skincare routine. For any questions or further assistance, please contact our customer service team.

If you have any questions after purchasing the product, please contact the after-sales email and we will handle it properly for you. We sincerely hope you can have an amazing shopping experience.



