



LUNDERG[®]

DESIGNED SOLUTIONS FOR BETTER HEALTH

QUICK SETUP GUIDE

Early-Alert Bed Alarm System

If you have any issue or question, contact us
directly at questions@lunderg.com

**SCAN
QR CODE**



PLEASE
watch the **EASY TO FOLLOW**
instructional video at
www.Lunderg.com/videoguide



Step 1

Preparation (bands)

Secure the bands to the clips. Repeat for all 4 clips



Loop the bands with the clips through the eyelets



Repeat for all 4 eyelets



Step 2 Preparation (stickers)


Get the 4 anti-slip stickers



Stick them on the BACKSIDE of the pad



Pull back the sheet / fitted sheet

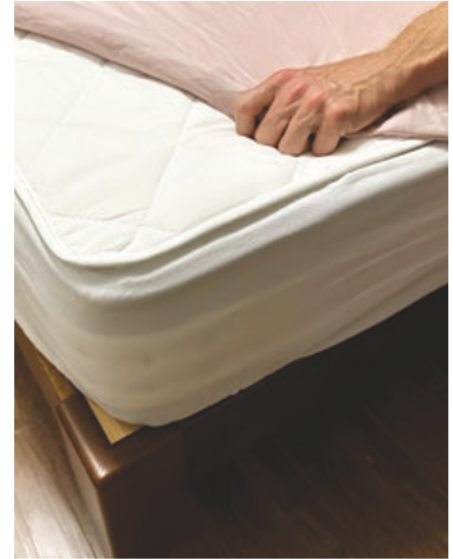

Step 3

Turn ON the monitor

Put batteries following + / -



Close cover with screwdriver



Turn ON




Step 4 Turn ON the Pad

Locate transmitter pocket and open it



Set transmitter to TX-Only (left) and close pocket

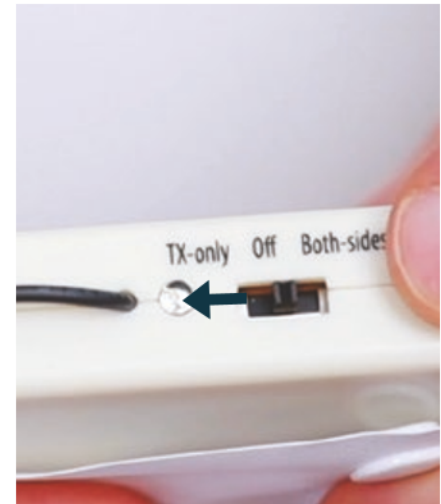

Step 5 Pairing


Hold "PAIR" button for 2-6 sec until status light turns STEADY GREEN



Press & Release the center of the pad while holding the monitor as image below



Sensor 1 light blinks once to indicate pairing was successful

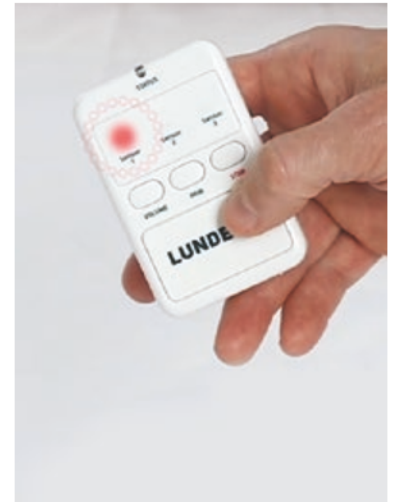

Step 6 Test

Press & Release
again the center of
the pad



⚠️ ALARM MUST SOUND
Cancel alarm by FIRMLY
pressing "STOP" button



If the Alarm
didn't sound
turn OFF the
pad & monitor
and go back
to Step 3.

Step 7 Pad Placement

Place pad across the bed at hips' height



Secure pad to bed using the clips (bands must be tensed)



Cover pad with the fitted sheet





**THE SYSTEM IS
READY TO USE**



Test system by lying down and exiting bed before use with loved one/patient.



This is just a quick guide to get you started. To unlock the full potential of your system, please read the user manual and check the instructional video!

Questions?



Contact us at
questions@lunderg.com

We will reply you within 24 hours,
7 days a week!



LUNDERG®

Lunderg.com



QUICK SETUP GUIDE

Chair Alarm System

If you have any issue or question, contact us
directly at questions@lunderg.com

SCAN
QR CODE



PLEASE
watch the **EASY TO FOLLOW**
instructional video at
www.Lunderg.com/chairvideo



Step 1

Turn ON the Monitor

Put batteries following + / -



Close cover with screwdriver



Turn ON



Step 2

Turn ON the Pad

Locate transmitter pocket
and open it



Turn ON



Step 3 Pairing

Hold "PAIR" button for 2-6 sec until status light turns STEADY GREEN



Press & Release the center of the pad while holding the monitor as image below



Sensor 1* light blinks once to indicate pairing was successful



*Sensor 2 or 3 if other devices are paired

Step 4 Test

Press & Release again
the center of the pad



⚠️ ALARM MUST SOUND
Cancel alarm by **FIRMLY**
pressing "STOP" button



If the Alarm
didn't sound
turn OFF the
pad & monitor
and go back to
Step 1.



**THE SYSTEM IS
READY TO USE**



Test the system by sitting down on the chair, wheelchair, recliner, or sofa, and standing up before use with a loved one/patient.

Questions?



Contact us at
questions@lunderg.com

We will reply you within 24 hours,
7 days a week!



LUNDERG®

Lunderg.com