

Homedics®

ModMassage Shiatsu Foot Massager with Heat



INSTRUCTION MANUAL AND WARRANTY INFORMATION

1-year limited warranty 

FMS-282H-AU

IMPORTANT SAFETY INSTRUCTIONS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK:

- **NEVER** leave the appliance unattended when plugged in. Unplug when not in use and before putting on or removing attachments.
- **DO NOT** operate under a blanket or pillow. Excessive heat may cause fire, electric shock, or injury.
- Close supervision is required when the appliance is used by, on, or near children or individuals with disabilities.
- Use the appliance only as described in this manual. **DO NOT** use attachments not recommended by Homedics.
- **NEVER** operate if the cord or plug is damaged, or if the product has been dropped or submerged in water.
- **DO NOT** attempt to repair the unit. There are no user-serviceable parts.
- **DO NOT** carry the unit by the cord, or use the cord as a handle.
- Keep the cord away from heated surfaces.
- **NEVER** operate the appliance with the air openings blocked. Keep openings free of lint, hair, and debris.
- **DO NOT** drop or insert any object into any openings.
- **DO NOT** use outdoors.
- **DO NOT** operate where aerosol (spray) products are being used.
- Use only on the body; **DO NOT** use near the eyes or other highly sensitive areas.
- **DO NOT** use this product while in bed.
- **NEVER** use while sleeping or drowsy.
- **ALWAYS** unplug after use.
- **DO NOT** use attachments not recommended by Homedics, specifically any attachments not provided with the unit.
- **NEVER** operate this product if it is damaged, or not working properly. Contact the Homedics customer service team for assistance.

SAVE THESE INSTRUCTIONS

CAUTION — PLEASE READ ALL INSTRUCTIONS BEFORE OPERATING

- Not intended for medical use.
- Consult your doctor before using if you:
 - Are pregnant
 - Have a pacemaker
 - Have concerns with your health
- **NEVER** leave the unit unattended, especially around children.
- **DO NOT** use the appliance for longer than instructed; extended use can shorten product life.
- **NEVER** use this product on swollen, inflamed areas or over skin eruptions.
- **DO NOT** use while bathing or in a steamy environment.
- This appliance should not be used by individuals with sensory impairment unless supervised.
- This appliance is intended for household use only.




NOTE: Only use the adaptor that was supplied to power this foot massager.

SPECIFICATIONS


Adapter Input	100-240V AC, 50/60 Hz, 1 A
Massager Input	24 V DC, 1.5 A

Made In China

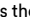INSTRUCTIONS FOR USE

1. Insert the power cord into a power outlet.
2. It is highly recommended that you wear socks while using this foot massager.
3. Place massager on the floor in front of your chair or seat and place feet on the massager.
4. To activate massage, press the  button once. The massage will begin.
5. For soothing heat, press the  button a second time. Allow several minutes for the massager to heat up.
6. To turn off the massage, press the  button a third time.

NOTE: The heat function may turn off after an extended period of use but the massage function will continue. When the temperature cools down again the heating function will turn on again automatically.

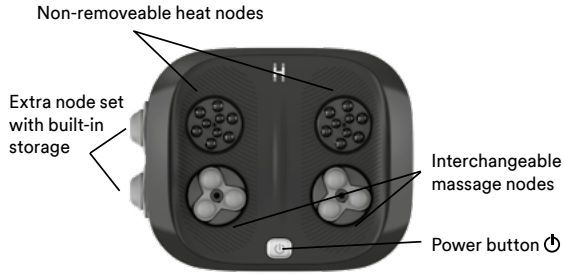

INTERCHANGEABLE MESSAGE NODES

This foot massager includes massage nodes that you can change to create your perfect massage experience.

SET 1.

DUAL-SIDE NODES – these nodes can be flipped over for two different massage sensations



Larger ball tip = gentler massage



Smaller ball tip = deeper massage

SET 2.

LARGE CONCENTRIC NODES – these nodes offer a more intense massage, great for massaging your feet arches.



HOW TO CHANGE THE NODES

1. First turn off the unit.
2. Remove the massage nodes.
3. Grip one of the nodes firmly and pull it straight up.
4. The extra set of nodes is docked on the side of unit. Pull outward to remove nodes.
5. Place the chosen node on the massage base and press firmly down.
Store the unused set of nodes on the docking prongs at the side of your foot massager.



MAINTENANCE

THE STORE

Place massager in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects. To avoid breakage, DO NOT wrap the power cord around the unit. DO NOT hang the unit by the cord.

THE CLEAN

Unplug the unit and allow it to cool before cleaning. Clean only with a soft, slightly damp sponge. NEVER allow water or any other liquids to come into contact with the unit.

- DO NOT immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, or strong chemicals that might or might not be flammable and/or damaging to the product to clean.
- DO NOT attempt to repair the massager. There are no user serviceable parts. For service, contact the Consumer Relations telephone number listed in the warranty section.

We or us means Homedics Australia Pty Ltd ACN 31 103 985 717 and our contact details are set out at the end of this warranty;

You means the purchaser or the original end-user of the Goods. You may be a domestic user or a professional user;

Supplier means the authorised distributor or retailer of the Goods that sold you the Goods in Australia and New Zealand; and

Goods means the product or equipment which was accompanied by this warranty and purchased in Australia and New Zealand.

For Australia:

Our Goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled, subject to the provisions of the Australian Consumer Law, to a replacement or refund for a major failure and for compensation for any other reasonably foreseeable loss or damage. You are also entitled, subject to the provisions of the Australian Consumer Law, to have the Goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure. This is not a complete statement of your legal rights as a consumer.

For New Zealand:

Our Goods come with guarantees that cannot be excluded under the Consumer Guarantees Act 1993. This guarantee applies in addition to the conditions and guarantees implied by that legislation.

The Warranty

Homedics sells its products with the intent that they are free of defects in manufacture and workmanship under normal use and service. In the unlikely event that your Homedics product proves to be faulty within 1 year from the date of purchase due to workmanship or materials only, we will replace it at our own expense, subject to the terms and conditions of this guarantee. The warranty period is limited to 3 months from the date of purchase for products used commercially/professionally.

Terms and Conditions:

In addition to the rights and remedies that you have under the Australian Consumer Law, Consumer Guarantees Act of New Zealand or any other applicable law and without excluding such rights and remedies warranty against defects: The Goods are designed to withstand the rigors of normal household use and are manufactured to the highest standards using the highest quality components. Whilst unlikely, if, during the first 12 months (3 months commercial use) from their date of purchase from the Supplier (Warranty Period), the Goods prove defective by reason of improper workmanship or materials and none of your statutory rights or remedies apply, we will replace the Goods, subject to the terms and conditions of this warranty.

We **DO NOT** have to replace the Goods under this Additional Warranty if the Goods have been damaged due to misuse or abuse, accident, the

attachment of any unauthorised accessory, alteration to the product, improper installation, unauthorised repairs or modifications, improper use of electrical/power supply, loss of power, malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance, transportation damage, theft, neglect, vandalism, environmental conditions or any other conditions whatsoever that are beyond the control of Homedics.

This Warranty does not extend to the purchase of used, repaired or second-hand products or to products not imported or supplied by Homedics Australia Pty Ltd, including but not limited to those sold on offshore internet auction sites.

This Warranty extends only to consumers and does not extend to Suppliers.

Even when we **DO NOT** have to replace the Goods, we may decide to do so anyway. In some cases, we may decide to substitute the Goods with a similar alternative product of our choosing. All such decisions are at our absolute discretion. All such replaced or substituted Goods continue to receive the benefit of this Additional Warranty for the time remaining on the original Warranty Period (or three months, whichever is the longest).

This Additional Warranty does not cover items damaged by normal wear and tear including but not limited to chips, scratches, abrasions,

discolouration and other minor defects, where the damage has negligible effect on the operation or performance of the Goods.

This Additional Warranty is limited to replacement or substitution only. As far as the law permits, we will not be liable for any loss or damage caused to property or persons arising from any cause whatsoever and shall have no liability for any incidental, consequential, or special damages.

This warranty is only valid and enforceable in Australia and New Zealand.

Making a Claim:

In order to claim under this Warranty, you must return the Goods to the Supplier (place of purchase) for replacement. If this is not possible, please contact our Customer Service department by email: cservice@homedics.com.au or at the address below.

All returned Goods must be accompanied by satisfactory proof of purchase which clearly indicates the name and address of the Supplier, the date and place of purchase and identifies product. It is best to provide an original, legible, and unmodified receipt or sales invoice.

You must bear any expense for return of the Goods or otherwise associated with making your claim under this Additional Warranty.

DISTRIBUTED BY

AUSTRALIA: Homedics Australia Pty Ltd,
14 Kingsley Close, Rowville, VIC 3178

NEW ZEALAND: CDB Media Ltd, 4 Lovell Court, Albany,
Auckland, New Zealand 0800 232 633

Homedics Customer Service:
Email: cservice@homedics.com.au



DISTRIBUTED BY

**AUSTRALIA: Homedics Australia Pty Ltd,
14 Kingsley Close, Rowville, VIC 3178**

**NEW ZEALAND: CDB Media Ltd, 4 Lovell Court, Albany,
Auckland, New Zealand 0800 232 633**