Sonic Electric Toothbrush

Operation Instructions

Enjoy Your Life Every Day!

Important safety information

Before using the product, please read this important information carefully and keep it properly for future reference.


Danger

Keep the charger away from the water source. Do not place or store it above or near any container with water, such as bathtub, washbasin, pool, etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before connecting it to the power supply.

Warning

Do not replace the power cord. If the power cord is damaged, please discard the charger. You should always use the original model charger for charging to avoid danger. If the product is damaged (brush head, toothbrush handle or charger), please stop using it.
The charger of this product can not charge other brands of electric toothbrush or any product with inductive/wireless charging function, either.
If possible, avoid charging for more than 12 hours. Do not let children play with this product.
Do not wash the brush head, toothbrush handle or charger in the dishwasher.
If you have had oral cavity or gum surgery in the past 2 months, please consult your dentist before using this product.
If the gum bleeds a lot after using this product or it will still bleed after 1 week, please consult your dentist.
If you feel discomfort or pain when using this product, please stop using it.
If any questions related to power/maker or other implantable device, please consult your doctor or the device manufacturer of the implantable device before use.
If you have medical problems, please consult your doctor before using this product.
This product is designed only for cleaning teeth and gum only. Do not use it for any other purpose.
If you feel any discomfort or pain, please stop using this product.
This product is for personal use, so it can not be used again by other patients in the dentist office or institutions.
Please stop using the brush head that splintered or with bent bristles. For bristles of traditional and tufting type, the bristles in the tuft will become lighter with use when a brush head should be replaced. It is recommended to replace the brush head every 3 months. If there are signs of wear, it should be replaced immediately. Do not use other brush heads that are not recommended by the manufacturer.
If your toothpaste contains peroxide, sodium bicarbonate or other bicarbonates (a common ingredient in whitening toothpastes), be sure to brush thoroughly with the brush head with soap and water after each use, which can prevent possible paste breakage.

Parts

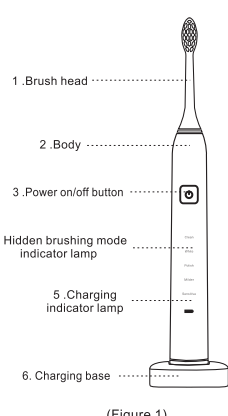1. Brush head
2. Body
3. Power on/off button
4. USB charging base
5. Charging indicator lamp
6. USB charger base

Introduction

Intelligent sonic-wave electric toothbrush (Figure 1)
1. Brush head
2. Body
3. Power on/off button
4. Hidden brushing mode indicator lamp
5. Charging indicator lamp
6. USB charger base
Note: Items in the packaging box may vary depending on the model purchased.


Product Information

Model: S1000-DC IV	Rated power: 2W
Product material: PC/ABS	Brush head: Nylon

Preparation for use

Brush head installation steps:
1. Align the brush head so that the direction of the bristles is the same as the front of the toothbrush handle (Figure 2).
2. Press the brush head down hard to the metal axis until it will not go down anymore.


Instructions

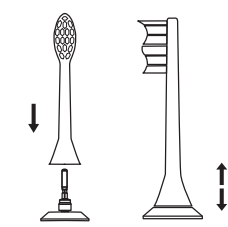Charging the product

Plug the charger into the power socket (Figure 3).
Please make sure that the machine body and the charging base are perpendicular to the horizontal plane during charging, otherwise the machine body will fall down and cause safety hazards.
Place the body on the charger base, and the charging indicator lamp will flicker until it is fully charged when the indicator is on.
The body is pre-charged, so you don't need to fully charge it before you use it for the first time. This charger can only charge one toothbrush handle at a time normally.

Charging base

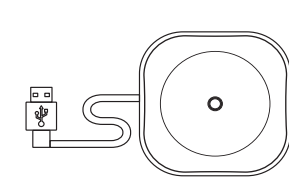Use this product

Use steps:
1. Insert the brush head carefully into the oscillating axle, and leave a certain gap between the brush head and the main body to reserve a working space for the normal work of the brush head.
2. To prevent toothpaste splashing, take a small amount of toothpaste on the brush head before use (generally speaking, sonic-wave electric toothbrush saves more toothpaste than manual toothbrush).
3. First hold the brush head and put it into the mouth, then start the toothbrush.
4. Move the toothbrush slowly with the bristles at 45° to the gum line.
5. Hold the toothbrush handle gently and move slowly in a circular motion. Do not press it hard!
6. Please stop using water with water after brushing.

The product can be used safely with:

Restoration of teeth (filling, crown, porcelain teeth).

Brushing mode

Cleaning mode
Effectively remove periodontal plaque (default mode). The duration of cleaning mode is 2 minutes, including four 30s intervals.

Whitening mode

2 minutes to deep whitening teeth, effectively improving tooth surface color.

Polishing mode

Brushing and polish the front teeth, and users can press and hold the power button to stop polishing in advance according to their own needs.

Milder mode (beginner mode)

It is suggested that mild mode can be used for the beginner which can effectively clean teeth and gently massage gum to improve gum health.

Sensitive mode

It is a super gentle mode for sensitive teeth and gums which lasts 2 minutes.

Toothbrush handle

1. Remove the brush head and wash the metal axle with warm water. Make sure to remove all remaining toothpaste.
Caution: Do not press the handle surface with rubber sealing ring on the metal axle, otherwise it may cause damage.
2. Wipe the entire surface of the toothbrush handle with a wet cloth.
Caution: Do not use rubbing alcohol of isopropyl, vinegar or bleach to clean the handle, otherwise it may cause discoloration of the handle.

Cleaning, service and storage

Clean the brush head: wash the brush head and bristles after each use, and dry them with natural air. Do not put the bristles, so this will deform them and cause them to fall off.
Clean the body: clean the handle surface with clean water, and scour it with soft cloth, and dry it with natural air. Do not press the objects on the rubber sealing ring on the metal axle, otherwise it may cause damage.

Replace brush head

Professional dentists recommend replacing the brush head every 2-3 months for the best results.

Clean the bottom of the brush handle: before charging, be sure to clean the bottom of the brush handle, wipe the bottom surface with a wet cloth, and then charge after the cleaning.

Storage

If you do not use this product for a long time, please fully charge the toothbrush handle, using the charger from the power socket, clean the two parts and store them in a cool and dry place away from direct sunlight.

Recovery

When discarding this product, please do not put it together with the general domestic waste. Instead, hand them to the recycling center designated by the government. This will be beneficial to the environment.

This product contains a built-in rechargeable battery, do not discard the rechargeable battery together with general domestic waste.

Follow your country's regulations on classifying and recycling electronic and electrical products and rechargeable batteries. Proper discarding of products will help to avoid adverse consequences on the environment and human health. Warranty period of intelligent sonic-wave electric toothbrush.

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How to Contact Us

Log in your amazon account > click "Your orders" > find the order ID > click "Contact Seller"

Description of charging time: it can last for 30-40 day (2 times) / day. 2 minutes/ times after continuous charging for 10 hours, it takes about 20-25 hours to fully charge, which can be used for 10-15 days (2 times/ day) / minutes/ times.

Description of battery decharge time: Charge in contrary pendant 8-10 hours, 30-40 pure 12 times/ day. 2 minutes/ times. / Full power 20-25 hours pure 10 times/ day. 2 minutes/ times. / Full power 12 times/ day. 2 minutes/ times.

How to Contact Us
Log in your amazon account > click "Your orders" > find the order ID > click "Contact Seller"

CE FC FDA