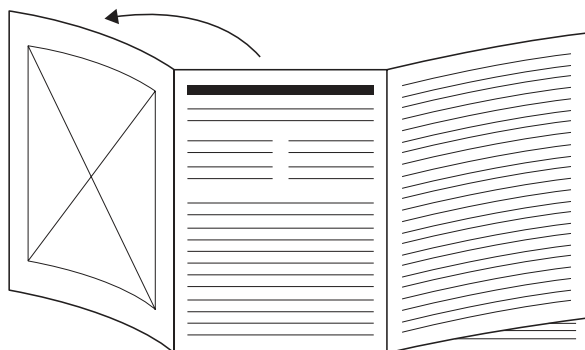
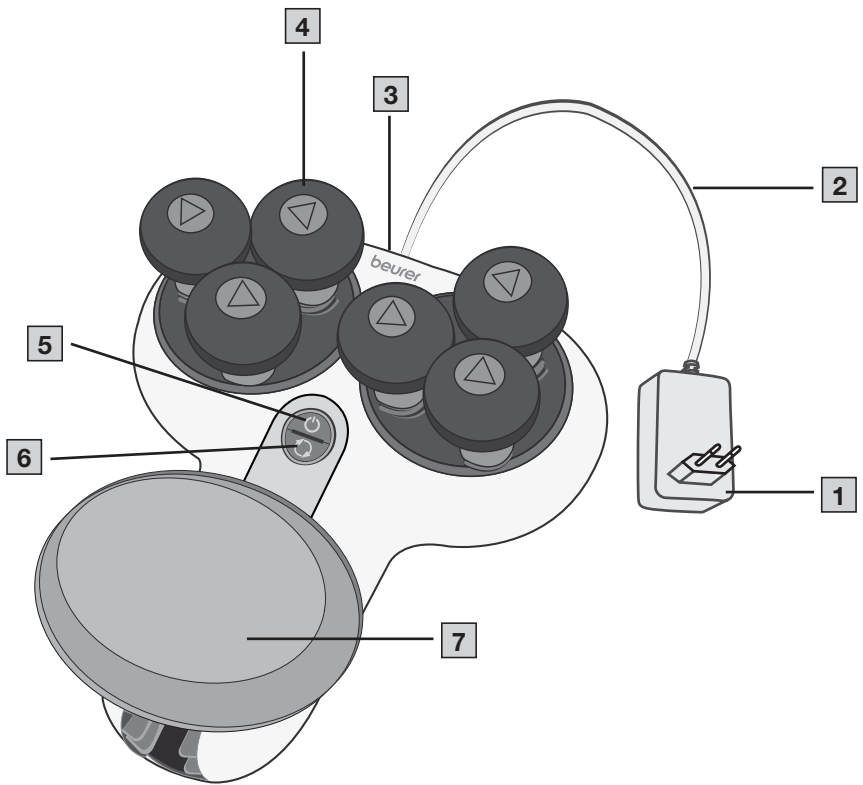
<b>DE</b>	Achillessehnen-Massagegerät Gebrauchsanweisung .....	4	<b>PL</b>	Urządzenie do masażu ścięgna Achillesa Instrukcja obsługi .....	45
<b>EN</b>	Achilles tendon massager Instruction for Use .....	10	<b>NL</b>	Massageapparaat voor de achilles- pezen Gebruiksaanwijzing .....	51
<b>FR</b>	Appareil de massage du tendon d'Achille Mode d'emploi .....	15	<b>DA</b>	Achillessene-massageapparat Betjeningsvejledning .....	57
<b>ES</b>	Aparato de masaje del tendón de Aquiles Instrucciones de uso .....	21	<b>SV</b>	Massageapparat för hälsenan Bruksanvisning .....	63
<b>IT</b>	Apparecchio per massaggi per il tendine d'Achille Istruzioni per l'uso .....	27	<b>NO</b>	Massasjeapparat for akillessenen Bruksanvisning .....	69
<b>TR</b>	Aşil tendonu masaj cihazı Kullanım kılavuzu .....	33	<b>FI</b>	Akillesjanteiden hierontalaite Käyttöohje .....	75
<b>RU</b>	Прибор для массажа ахиллова сухожилия Инструкция по применению .....	39			

- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Unfold page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- ES** Despliegue la página 3 antes de leer las instrucciones de uso.
- IT** Prima di leggere le istruzioni per l'uso aprire la pagina 3.
- TR** Kullanım kılavuzunu okumadan önce 3. sayfayı açın.
- RU** Перед чтением инструкции по применению разложите страницу 3.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit om de gebruiksaanwijzing te kunnen lezen.
- DA** Fold side 3 ud, før du læser betjeningsvejledningen.
- SV** Vik ut sid. 3 innan du läser bruksanvisningen.
- NO** Åpne side 3 før du leser bruksanvisningen.
- FI** Käännä sivu 3 auki ennen käyttöohjeen lukemista.

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 Read these instructions for use carefully. Follow the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

 **WARNING**

- The device is intended only for domestic/private use, not for commercial use.
- This device can be used by children over the age of eight and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

**1. INCLUDED IN DELIVERY**

- Achillomed®
- Mains part
- These instructions for use

**Contents**

1. Included in delivery .....	10	7. Use .....	13
2. Signs and symbols .....	11	8. Cleaning and storage .....	14
3. Intended use .....	11	9. What if there are problems? .....	14
4. Safety notes .....	11	10. Technical data .....	14
5. Device description .....	12	11. Warranty .....	14
6. Initial use .....	12		

## 2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	<b>Warning</b> Warning instruction indicating a risk of injury or damage to health
	<b>Important</b> Safety note indicating possible damage to the device/accessory
	<b>Product information</b> Note on important information
	For indoor use only
	Read the instructions
 	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard
	Protection class II device
	The electronic device must not be disposed of with household waste
	Polarity of d.c. power connector
	Manufacturer
	<b>CE labelling</b> This product satisfies the requirements of the applicable European and national directives.
	United Kingdom Conformity Assessed Mark
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.

	Safety isolating transformer, short-circuit proof
	Energy efficiency level 6
	Switch mode power supply unit
<b>IP 30</b>	Device protected against foreign objects $\geq 2.5$ mm
	Do not use the power adapter if the plug is damaged.
	Importer

## 3. INTENDED USE

This device is exclusively intended to massage the human Achilles tendon. The device is only suitable for **private use**, and not for medical or commercial purposes. It cannot replace medical treatment. Do not use the device if one or several of the following warnings applies to you. If you are unsure whether the device is suitable for you, please consult your doctor.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

## 4. SAFETY NOTES

Read these instructions for use carefully! Non-observance of the following information may result in personal injury or material damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the device to third parties.

### WARNING

Do not use the device

- On animals
- If you suffer from a medical abnormality or injury of the Achilles tendon (e.g. open wounds)
- If you have a pacemaker
- During pregnancy
- On the face (eyes), on the larynx or any other particularly sensitive body parts
- For foot massages
- While sleeping
- After the consumption of medicines or alcohol (reduced cognition!)
- In vehicles

You should consult your doctor before using the device, especially

- If you suffer from a severe illness or have recently had surgery on your Achilles tendon
- If you have a pacemaker, implants or other aids
- If you suffer from thrombosis, diabetes or unexplained pain

## **WARNING**

Keep packaging material away from children. There is the risk of suffocation.

## **Electric shock**

### **WARNING**

As with every electrical device, use the Achillomed® carefully and cautiously in order to avoid hazards due to electric shock. Therefore, observe the following instructions for use:

- Only use the device with the mains part provided and at the mains voltage specified on the mains part
- Never use the device if it or its accessories show visible signs of damage
- Do not use the device during a thunderstorm

In the event of defects or operational faults, switch the device off immediately and disconnect it from the mains socket. Do not pull on the mains cable or on the device to disconnect the mains part from the socket. Never hold or carry the device by the mains cable. Keep the mains cable away from hot surfaces. Do not crush, bend or twist the mains cable. Do not insert needles or other sharp objects into the device.

Ensure that the device, the mains part and the mains cable do not come into contact with water, steam or other liquids. Only use the device in dry rooms (e.g. never in the bath, sauna).

Never reach for a device that has been submerged in water. Unplug the mains part from the socket immediately. Do not use the device if it or its accessories are visibly damaged.

Do not drop the device or subject it to any impacts.

## **Repairs**

### **WARNING**

- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact Customer Services or an authorised retailer.

## **Risk of fire**

### **WARNING**

There is a risk of fire if the device is used improperly or these instructions for use are ignored.

Therefore, observe the following instructions for use:

- Never use the massager underneath a cover, such as a blanket, pillow, etc.
- Never use the device near petrol or other highly flammable substances.

## **Handling**

### **IMPORTANT**

The device must be switched off and disconnected from the mains after every use and before cleaning.

- Do not place any objects between the rotating massage heads. Ensure that the moving parts can move freely at all times.
- Do not sit, lie or stand with your full weight on the moving parts of the device.
- Do not place any objects on top of the device.
- Do not expose the device to high temperatures.

## **Disposal**

### **IMPORTANT**

For environmental reasons, do not dispose of the device in the household waste at the end of its service life. Dispose of the device at a suitable local collection

or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

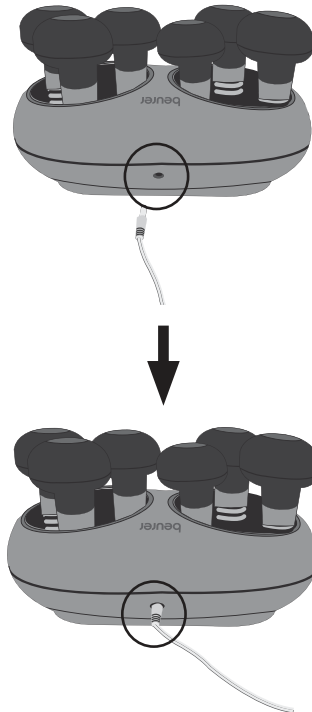
## **5. DEVICE DESCRIPTION**

The associated illustrations are on page 3.

- 1** Mains part
- 2** Mains cable
- 3** Mains part socket
- 4** 2x3 rotating massage heads
- 5**  -button  
Press once = low massage setting  
Press twice = high massage setting  
Press 3 times = switch off device
- 6**  -button  
(2 massage directions)
- 7** Footrest with height and width adjustment

## **6. INITIAL USE**

- Remove the packaging.
- Check the device, mains part and mains cable for damage.
- Insert the mains cable connector plug into the mains part socket on the device.
- Ensure that the mains cable does not pose a trip hazard.


## 7. USE

### **i** Note


- Do not use the Achillomed® for longer than 15 minutes a day. Longer massages can lead to over-stimulation of the Achilles tendon, creating tension instead of aiding relaxation.
- The Achillomed® features an automatic switch-off function, which is set to the maximum operating duration of 15 minutes. If you change the massage speed or direction during the application, the automatic switch-off timer starts again from the beginning (15 minutes).
- You can use the Achillomed® with bare feet or with socks on. To ensure an optimal massage effect, we recommend that you use the Achillomed® with bare feet.

To start the application with the Achillomed®, proceed as follows:

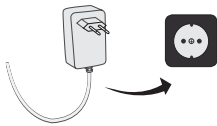1. Plug the mains part into a suitable socket. Ensure that the mains cable does not pose a trip hazard.

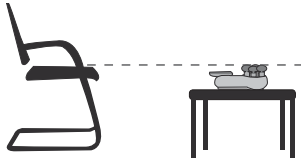
2. Place the Achillomed® on a raised surface. Ensure that the Achillomed® is positioned roughly at the height of your seat (e.g. chair or sofa).


3. Sit down on the chair (or any kind of seat) in front of the positioned Achillomed®.

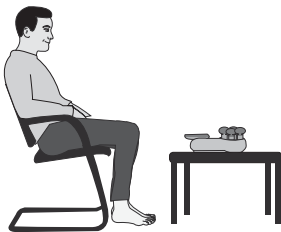**⚠ WARNING** Trip hazard! Never use the Achillomed® when standing up.

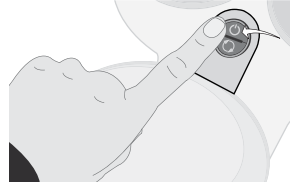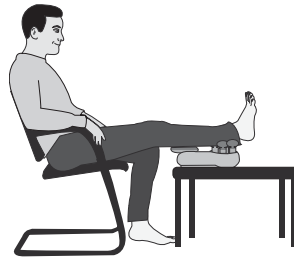
4. Press the  button on the Achillomed®. The massage heads start to rotate.

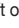
**⚠ WARNING** Danger of crushing! Do not force or jam your hands, fingers or hair between the rotating massage heads.

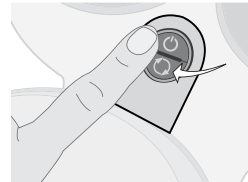
5. Position your leg carefully with the Achilles tendon to be treated on the Achillomed®. Ensure that your Achilles tendon is in the centre between the rotating massage heads. The massage should feel pleasant and relaxing at all times. If the massage is painful or uncomfortable, stop the massage or change the position of your foot.

You can individually adjust the position of your foot using the footrest with height and width adjustment. You can place your other leg on top of the leg being treated. This way, you can also regulate the massage intensity.

6. To change the massage direction, press the  button. Depending on the massage direction, the left or right LED of the  button lights up. If you press and hold the  button for half a second, the Achillomed® switches to a massage programme. The left and right LEDs of the  button light up. The massage direction is now controlled automatically. Press the  button to switch the automatic massage programme off again.

7. To increase the massage speed, press the button  again. Pressing the button  again switches the Achillomed® off.

**i** Note We recommend that you start with the lower massage speed so that you get used to the Achillomed® gradually.

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Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany
Beurer UK Ltd., Suite 16, Stonecross Place Stonecross Lane North WA3 2SH Lowton