MERACH

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
PANTONE 2394 C

USER MANUAL

Massage Gun

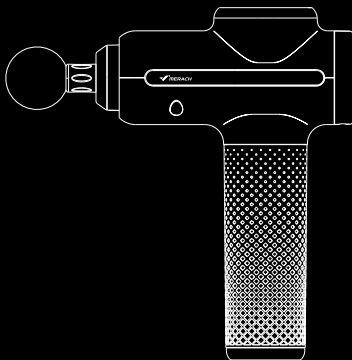Questions or Concerns?

IMPORTANT!

Please read all details before use, and keep this user manual for future reference.
PLEASE CONTACT US BEFORE THE RETURN: Mon-Fri, 9:00 am-5:00 pm PST/PDT
Phone: 1- (877) 3563730 Mon-Fri 8am-5pm (PST)
support.eu@merach.com

WHO WE ARE

Our customers aren't just important; their health, fitness, and pure enjoyment for life is at the core of everything we do at MERACH. From innovative fitness research and cutting-edge product development, to our exceptional purchasing & after-sales service, we are our customers' life-long partner in crafting a better fitness future.

Through quality fat-burning & muscle-building fitness equipment such as our quality rowing machines and exercise bikes, plus our innovative MERACH APP, and digital courses, we give the opportunity and support to all to pursue their health and fitness goals and discover a new way of healthy living.

Welcome to MERACH

WARRANTY REGISTRATION

Please register your warranty to keep yourself protected !

Register your warranty via
merachfit.com/pages/warranty-europe

Email your item model to
support.eu@merach.com

If purchased outside US/Canada: To register your product warranty, contact your local distributor.

For details regarding product warranty or if you have questions or problems with your product, please contact your local distributor.

CONTENTS • ENGLISH

Warnings and Important Safety Instructions.....	01
Product Overview.....	05
Product Usage Instructions.....	06
Full-body Massage.....	10
Warranty Information.....	13

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WARNINGS AND IMPORTANT SAFETY INSTRUCTIONS

It is important to know the following before use, otherwise, it may cause damage or become dangerous!

- **WARNING:** To reduce risk of burns, fire, electric shock, or injury to persons, read the following warnings carefully.
- Please check the main unit and cable before using, if damaged, it must be replaced by the manufacturer or its appointed professional to avoid danger.
- Please use the 5V=2A charger.
- Consult your doctor before using this product. Do not use if injured.

DO NOT use under any of the following conditions:

- Aneurysms, bleeding, use of blood thinners.
- Heart disease, pacemaker, or defibrillator.
- Pregnancy
- Cancer
- For adult use only. Keep away from children and pets. This device is not intended for people (including children) who are physically weak, slow or have mental disorders, except under the guidance or help of people responsible for their safety. Children shall be supervised to ensure that they will not play with this product.
- Use the MERACH massage gun only as instructed in this manual.
- Keep away from liquids and heat sources.
- Keep ventilation ports free from dust and debris.

- Never leave the MERACH massage gun unattended while it is running. Unplug the unit after charging.
- **DO NOT** remove screws or attempt to disassemble.
- Close supervision is necessary when this unit is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never drop or insert any object into any opening.
- **DO NOT** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **DO NOT** put fingers, hairs or any body parts to touch between the massage head and machine while the machine is running to avoid an accident.
- This product is non-waterproof.
- **DO NOT** use the product when charging.
- Use a towel or soft cloth to wipe the surface to moisture.
- **CAUTION:** Risk of Fire and Burns. Do Not Open, Crush, Heat Above 140° F or Incinerate. Follow Manufacturer's Instructions.

Do not touch this part while the machine is running without massage head.

Do not touch this position when the device is running with the massage head.

Maintenance:

- In storage or traveling, keep it away from hard objects and avoid bumps and scratches.
- The external cleaning of the equipment can be wiped with a soft cloth dipped in alcohol.
- If not used for a long time, it is recommended to charge once a month. Maintain power at least 20% when in use.

Technical Details:

Name : MERACH Massage Gun

Model : MR-M019

Peak Power: 40W

Battery Capacity : 2000mAh

Size : 2.4*6.5*8.9inch/60*165*225mm

Charging Time : 3 Hours

Input : 5V $\overline{\text{---}}$ 2A

This symbol stands for "Conformité Européene", which means "Conformity with EU directives". With the CE-marking the manufacturer confirms that this product complies with applicable European directives and regulations.

This symbol stands for "UK Conformity Assessed", which means "Conformity with British Standards 11". With the UKCA-marking the manufacturer confirms that this product complies with applicable British Standards.

SIMPLIFIED EU DECLARATION OF CONFORMITY

- Hereby, Zhejiang Yulu Electronic Technology Co.,Ltd. declares that the fitness equipment MERACH MR-M019 is in compliance with Directive 2014/53/EU.
- The full text of the EU declaration of conformity is available at the following internet address: <https://merachfit.com/pages/declarations-of-conformity>

NOTES ON DISPOSAL

Packaging materials are raw materials and can be recycled. Separate the packaging materials correctly and dispose of them properly in the interest of the environment.

Your local administration can provide additional information.

Used devices do not belong in household waste!

Dispose of used devices properly! Doing so helps fulfill our responsibility to protect the environment. Your local administration can provide information on collection points and opening hours.

Batteries do not belong in household waste!

As the end user, you are legally required to return discharged batteries to a collection point. You can take your old batteries to the public collection points in your community, or you can take them to any place where batteries of that type are being sold.


Manufactured by:

Zhejiang Yulu Electronic Technology Co.,Ltd.
Room 805, 8th Floor, Xianfeng
Technology Building, 298 Weiye Road,
Binjiang District, Hangzhou,
Zhejiang Province, China
yulu_mrk@merach.com

Imported by:

HANGZHOU JINGGE E-COMMERCE LTD
Room 1205, Xianfeng Technology Building,
298 Weiye Road, Puyan Subdistrict,
Binjiang District, Hangzhou,
Zhejiang Province, China
jingge_mrk@163.com

PRODUCT OVERVIEW

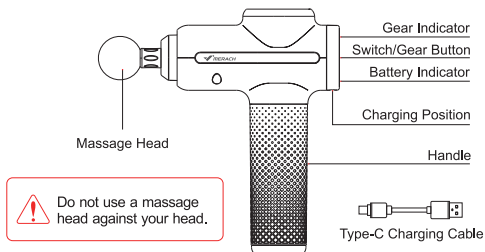The massage gun is a portable rechargeable health relaxation device with strong power, long range, multi-stop adjustment and other characteristics. The device vibrates through different frequencies to relax the muscles and relieve the discomfort of the muscles in various parts of the sedentary and recumbent population.

Large Ball Head
(Foam)

D-type Head
(PP)

Bullet Head
(PP)

Small Ball Head
(Foam)

U-type Head
(PP)


Flat Head
(PP)

PRODUCT USAGE INSTRUCTIONS

Control Panel Description


Power ON/OFF

- When the power is off, press and hold the switch for 2 seconds to enter standby mode, and the power indicator light will light up to show the current power.
- In standby mode, short press the power button once to enter working mode, the motor starts, and the gear 1 light is on; press the power button in sequence to adjust the gear, and the corresponding gear light will light up.

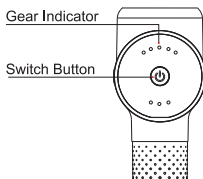Battery Indicator

Displays the current battery level after power on; when the battery is fully charged, the battery indicator light is always on.

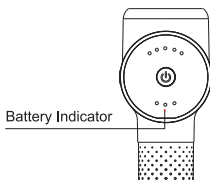
Gear Adjustment

In standby mode, short press the switch button to enter gear 1, and the corresponding gear light 1 will light up; continue to short press the switch to enter gear 2: short press the switch again to enter gear 3; and so on, the highest gear is gear 5; short press the switch button in gear 5, and the gear will enter standby mode. Adjust the gear, and the corresponding gear light will light up. (Standby → Gear 1 → Gear 2 → Gear 3 → Gear 4 → Gear 5 → Standby, cycle in sequence.)


NOTE:

In order to protect bones and fascia and avoid long-term massage and damage, this product is set to automatically stop the intelligent protection mode after 10 minutes of opening.

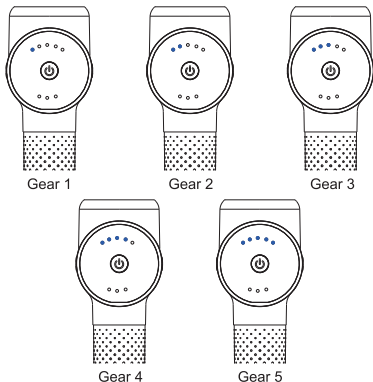Use/Replace Massage Head

When the Massage Gun is powered off, select the corresponding massage head according to the massage needs of different parts, and insert it into the round hole in front of the Massage Gun. When replacing the massage head, please turn off the Massage Gun first, take out the original massage head and replace it.

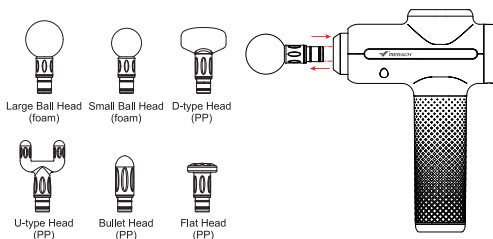Ball Head: suitable for massaging large muscle groups such as thighs, waist and back, arms, etc.



D-type Head: suitable for linear muscles such as arms and legs.

U-type Head: suitable for massaging shoulders, necks, both sides of the spine, and Achilles tendons.

Bullet Head: suitable for impacting deep tissues such as meridians, palms, soles of feet, etc.

Flat Head: suitable for relaxation and shaping of various muscle parts.

About Charging

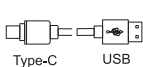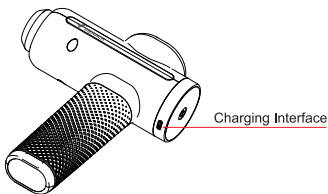(User Provided)

Note: 5V $\rightarrow$ 2A and above cell phone charging head is recommended.

- Be sure to fully charge before first use.
- Charging time is about 4h.
- Insert one end of the Type-C charging cable into the Massage Gun and the other end into the charger to enter the charging state.
- When charging: the power indicator light flashes white; when fully charged: the power indicator light is always white.
- Please unplug the charger in time after fully charged.
- Please charge in time after the power is used up. If you do not charge for a long time, the Massage Gun may not be used.

To ensure safety, do not use the Massage Gun while it is charging

FULL-BODY MESSAGE

BACK

We recommend the Ball Head, or Flat Head here. Run the level of massage gun stalls according to your body force. Avoid keeping the massage head in one area for more than 3 minutes.


ARM

The Ball Head and Flat Head attachments are recommended for arms. Always avoid direct contact with your bones. Generally, you should spend about 2 minutes massaging this area.


HIP

The Ball Head and Flat Head attachments are recommended for hips. Avoid sensitive area. Stay in an area for about 30 seconds before moving along

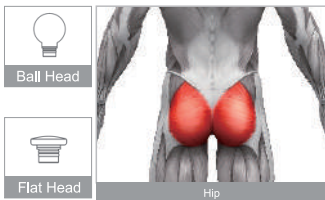LEG

We recommend the Flat Head, with legs you want to target the quads and the calves. Especially the quads may need something a little more intense stay in an area for about 30 seconds before moving along. You can also massage your calves after a leg workout session or after running. Massage one leg for 1-2 minutes and then move on to the other leg.

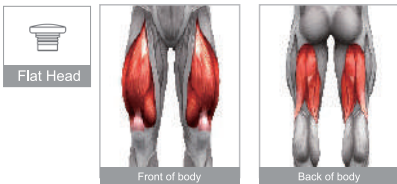CALF

We recommend massaging your calves with a Bullet Head, especially if you need to relieve calf pain during a run. Start with low intensity.

You can cross one foot over the other knee to get good posture.

Bullet Head

Calf

WARRANTY INFORMATION

Product Name	MERACH Massage Gun
Model	MR-M019
Default Warranty Period	12 Months
For your own reference, we strongly recommend that you record your order number and date of purchase.	
Date of Purchase	
Serial Number	

Terms & Policy

The MERACH products are premium in material, craftsmanship, and service. MERACH 1-year limited warranty is effective from the date of purchase. MERACH will replace defective products based on eligibility. Refunds are available to the original purchasers of our products within the first 30 days of purchase. This warranty extends only to personal use, doesn't apply to the commercial, rental, or any other uses for which the product is not intended for. There are no warranties other than the warranties expressly set forth with each product.

This warranty is non-transferable. MERACH is not responsible for any damages, losses or inconveniences caused by equipment's failure, user's negligence, user's abuse, or improper operation that does not follow the included user manual.

This warranty does not apply to the following situations:

- Damage due to abuse, accident, alteration or vandalism.
- Improper or inadequate maintenance.
- Damage in return transit.
- Unsupervised use by children under 14 years old.

MERACH and its subsidiaries are only responsible for damages caused by the intended use or as instructed in the user manual. Some states do not allow this exclusion or limitation of incidental or consequential losses, so the above disclaimer may not apply to you. This warranty gives you specific legal rights and you may also have other rights, which may vary from state to state.

All expressed and implied warranties, including the warranty of merchantability, are limited to the period of the limited warranty.

Additional 12-month Warranty:

You can enjoy an additional 12-month warranty. Just active your warranty by sending the product's model name via **register@merach.com**, or log onto **merachfit.com/pages/warranty** and enter your order number (i.e. Amazon) within the first 14 days of purchase to register your new product for the extended warranty. **If** you are unable to provide the order number for your product, please type a short note in the order number field along with the date you received your product.

Defective Products & Returns:

Should your product prove defective within the specified warranty period, please contact the Customer Support via **support.eu@merach.com** with your invoice and order number. **DO NOT** dispose of your product before contacting us. Once our Customer Support team has approved your request, please return the product with a copy of your invoice and order number.

Customer Support

Your satisfaction is our goal! Should you encounter any issues or have any questions about your new product, feel free to contact our Customer Support Team via **support.eu@merach.com** or Phone Support:

US/CA ☎ **1- (877) 3563730 Mon-Fri 8am-5pm(PST)**

TikTok Shop Exclusive after-sales: **service@merach.com**

***Please have your Order ID and serial number (usually on the equipment) ready before contacting Customer Support.**

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