


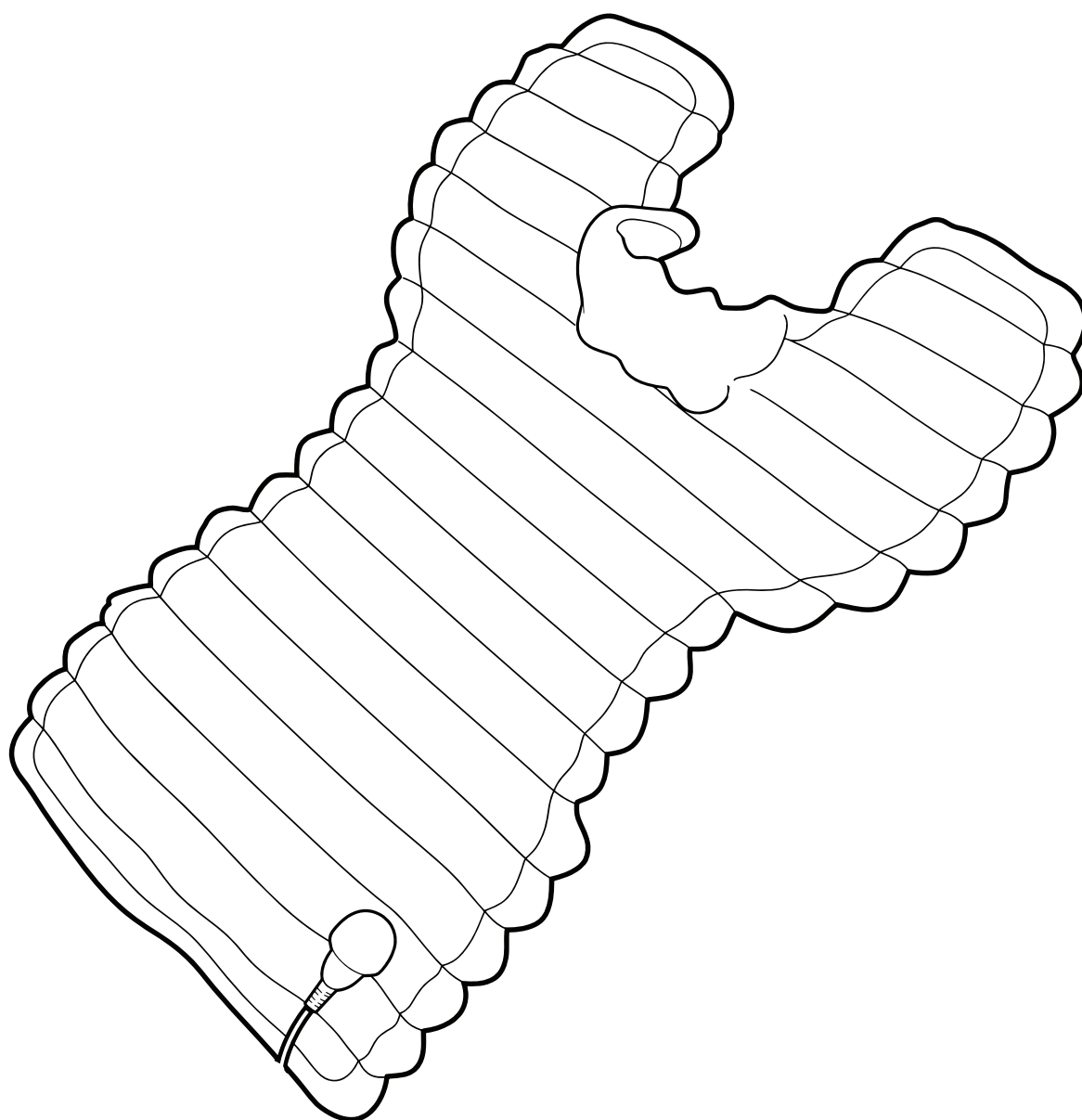

HEATED BACK PAD

Product Instructions

Product No. 2000346

vonhaus


Let's get started!

Visit our product hub for useful tips and tricks on getting the most out of your new product.

SCAN ME



Product Diagram

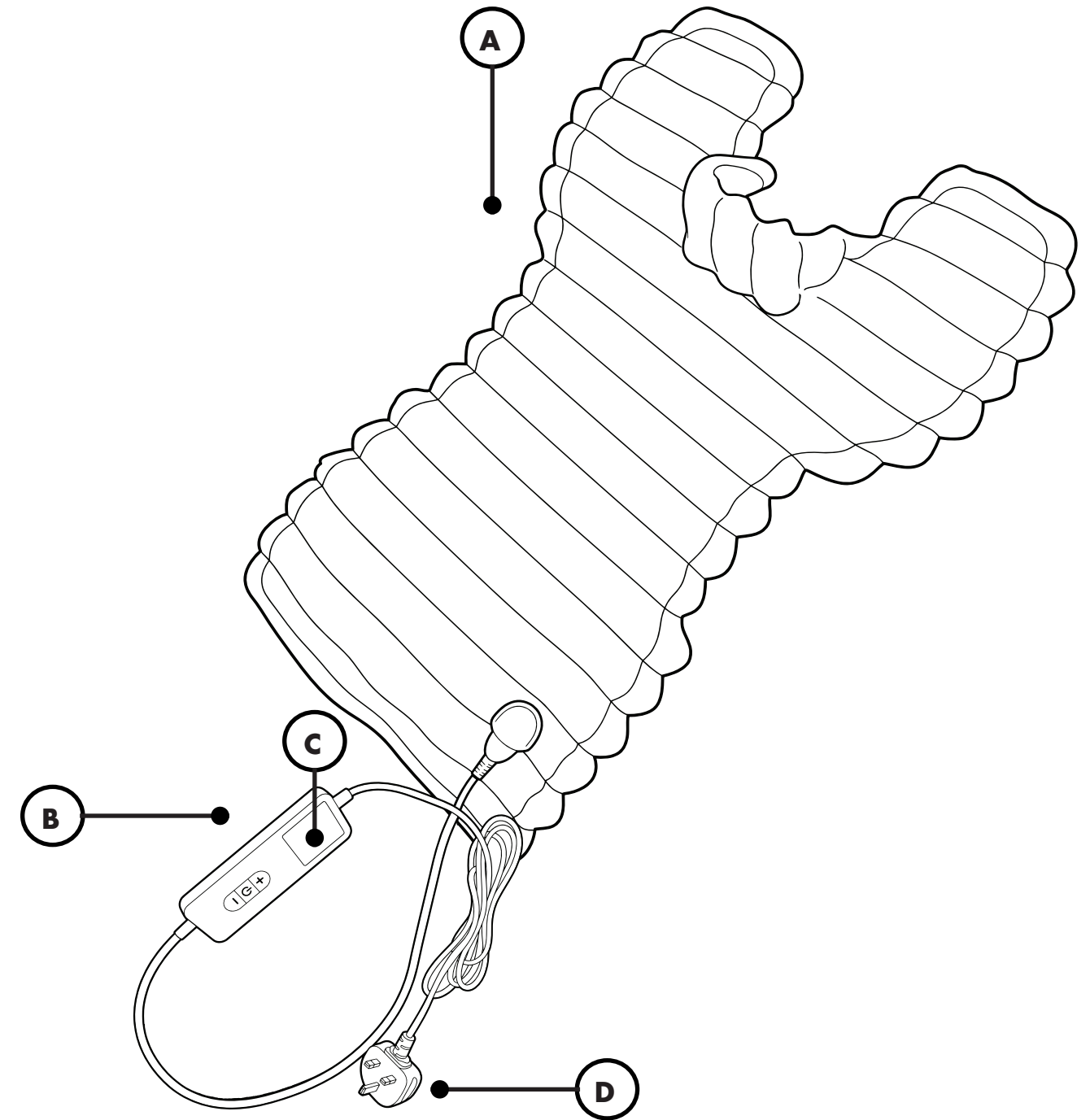


A. Heated Back Pad
B. Controller

C. Digital Display
D. Plug

Need further help?
support@vonhaus.com

See our full range at
www.vonhaus.com

Operation

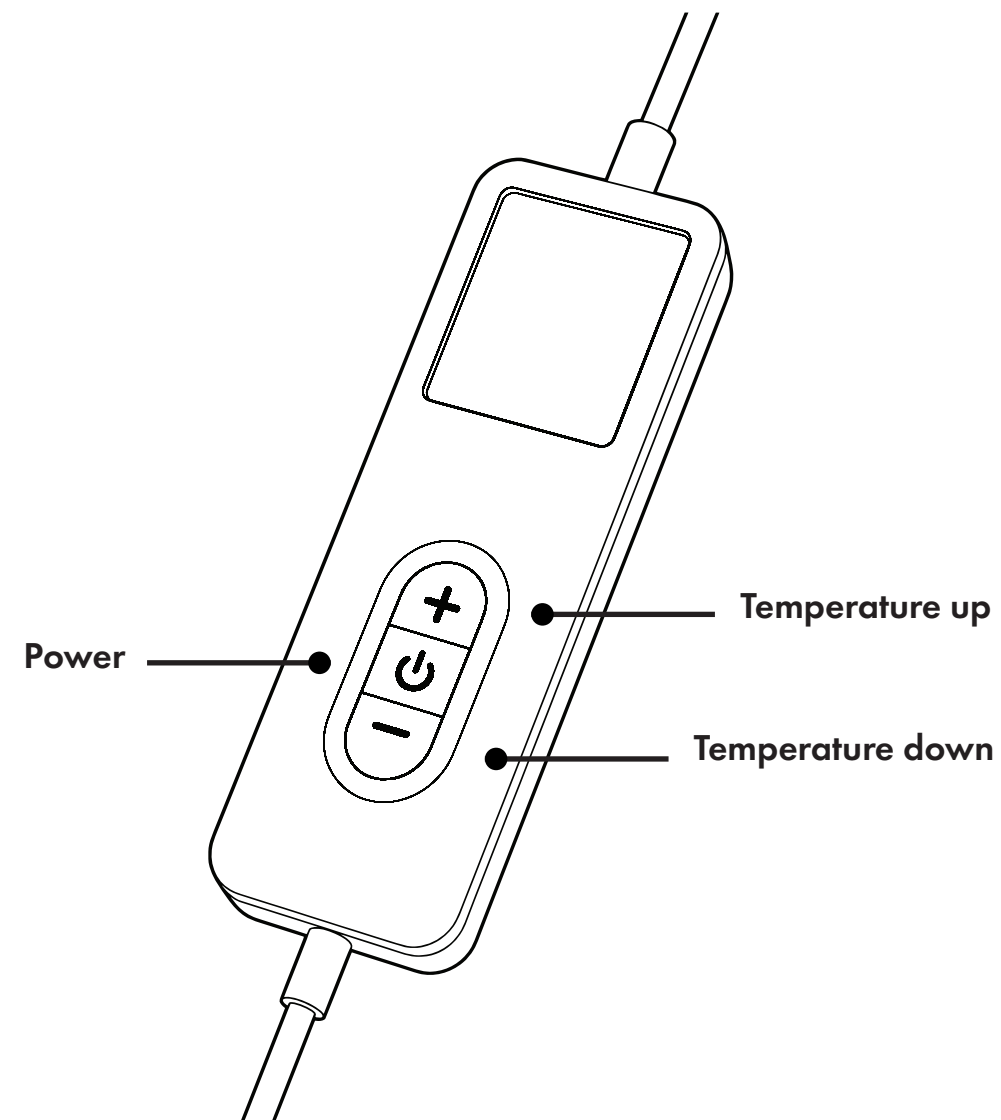

Getting Started

1. Plug the back pad in to the mains unit and switch on
2. Press the **power** button and the display screen will illuminate. The back pad will default to the lowest heat setting (level 1), and it will begin to gradually heat up
3. Press **+** to increase the heat setting (between 1-9)
4. Press **-** to decrease the heat setting (between 1-9)

Note: For continuous use, only use heat setting 1.

5. Press the **power** button to switch off. Unplug from the mains and allow to cool down before storing away

Usage & Safety




Please read all instructions carefully before use and retain for future reference.



Do not use outdoors.



Unsuitable for children under 3.

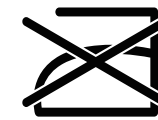


Do not iron.



Do not wring dry.



Do not use folded or tucked.

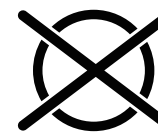


Do not dry clean.



Do not bleach.



Do not pierce or use sharp objects near the blanket.

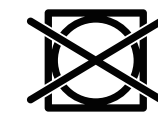

Do not tumble dry.



Do not wash above
30°C / 86°F

Information

Cleaning and Maintenance

1. Unplug the back pad from the mains units and allow to cool before cleaning or storing away.
2. Disconnect the cable from the back pad and put in the washing machine. Wash at 30°C on a delicate cycle.
3. Do not tumble dry. Lay the back pad flat and allow to thoroughly dry before connecting the cable again.

Technical Specifications

- Rated Voltage: 220-240V~
- Rated Power: 100W
- Rated Frequency: 50/60Hz
- Off Mode Consumption: 0.27W

Safety Warnings

- Do not fold or crease the back pad during use
- Do not pierce with sharp objects
- Must not be used by children under the age of 3, those insensitive to heat or if the user is asleep
- Inspect regularly for signs of wear or damage. If there are signs of damage or it does not work, please contact support@vonhaus.com
- Prolonged use at a high heat could lead to burns
- Arrange the cables correctly to avoid risk of entanglement, strangulation, or tripping
- Do not fall asleep while the heat pad is switched on. If extended use is expected, set the temperature to position 1 before sleeping.
- Prolonged use at a high setting may cause skin burns.

Information

Disposal



Do not dispose of this product as unsorted household waste. Separate collection for special treatment is required. Return your used device via a collection and return scheme.

Customer Service

If you are having difficulty using this product and require support, please contact support@vonhaus.com

Warranty

To register your product and find out if you qualify for a free extended warranty please go to www.vonhaus.com/warranty

Please retain a proof of purchase receipt or statement as proof of the purchase date.

The warranty only applies if the product is used solely in the manner indicated in the warnings page of this manual, and all other instructions have been followed accurately.

Any abuse of the product or the manner in which it is used will invalidate the warranty.

Returns

Returned goods will not be accepted unless re-packaged in its original packaging and accompanied by a relevant and completed returns form. This does not affect your statutory rights.

No rights are given under this warranty to a person acquiring the appliance second-hand or for commercial or communal use.

Copyright

All material in this instruction manual are copyrighted by DOMU Brands.

Any unauthorised use may violate worldwide copyright, trademark, and other laws.

Thank You

Thank you for purchasing your product/appliance.

Should you require assistance with your purchase, please contact us at support@vonhaus.com

Vonhaus is a registered trademark of DOMU Brands Ltd.

Made in China for DOMU Brands Ltd.
Unit 30, Stakehill Industrial Estate
Middleton, Manchester, M24 2RW

EU Authorised Representative

Authorised Representative Service
77 Camden Street Lower
Dublin, D02 XE80
Ireland

Let's get started!

Visit our product hub for useful tips and tricks on getting the most out of your new product.

SCAN ME

