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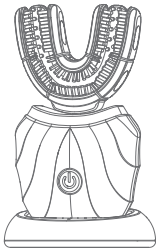
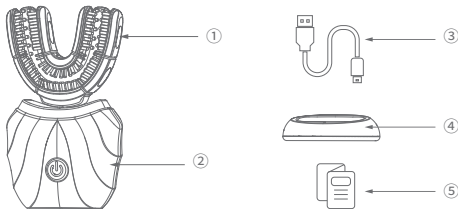
③ USB Type-C 充电线 × 1

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

360° U-Shaped Automatic Sonic Toothbrush


User Manual

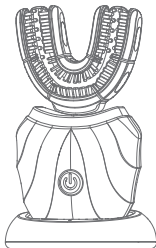Package Contents:

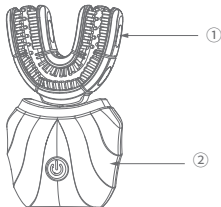- ① Brush head × 1
- ② Handle × 1
- ③ USB Type-C charging cable × 1

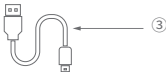③

④

⑤

- ④ Charging base × 1
- ⑤ User manual × 1

Note:

① Brush head type may vary by model. Available in nylon bristle or silicone. Please refer to the actual product received.

Product Features:

Wireless Inductive Charging

3-hour charging, 15-day battery life for convenient use

Five Cleaning Modes

Gentle, Clean, Whitening, Polish and Massage

45-Second Full-Mouth Cleaning

360° coverage. Completes a full brushing cycle in just 45 seconds

Powerful Sonic Motor

Wide swing vibration with high-quality nylon or silicone bristle brush heads for a deep cleaning experience

IPX7 Waterproof

Safe for rinsing under water. Easy to clean without fear of water damage

How To Use:

- ① Align the brush head hole with the motor shaft and insert until securely in place.

- ② Use special foaming toothpaste; if using regular toothpaste, dampen the brush head first, then apply toothpaste to both sides of the brush head

- ③ Place the brush in your mouth and gently bite down so that the upper and lower teeth fit into the mouthpiece.

- ④ Press the power button to start cleaning.

- ⑤ Gently move the handle left and right to allow deeper cleaning.

Note:

Silicone brush heads can be trimmed to suit individual mouth size.

- ⑥ Rinse the brush head with clean water after use and allow it to air dry. Wipe the handle with a dry cloth and keep it dry.

How To Charge:

- ① Insert the USB Type-C plug into the charging base, and connect the Type-A end to a 5V USB power source.
- ② The indicator light on the charging base will turn on, showing the power is connected.
- ③ Place the handle vertically on the base slot. The RGB light on the handle will pulse while charging and remain solid when fully charged.

Specifications:

U-Shaped Automatic Sonic Toothbrush

- Model: JS508
- Motor: Wide-swing sonic motor
- Vibration: Adjustable, 4,800–36,000 VPM
- Modes: 5 modes
- Additional features: 480nm blue light for antibacterial care
+ 640nm red light for restorative care
- Ambient Light: RGB color-changing ambient lights
- Battery: 3.7V lithium battery
- Battery Life: 15–20 days
- Waterproof rating: IPX7
- Size: 6.2 × 6.8 × 2.9 cm, 6.2 × 11 × 2.9 cm (with brush head)
- Timer: 45 seconds
- Surface Material: Silicone + Spray Coating
- Brush Head: Various silicone brush heads / nylon bristle brush heads

Wireless Charging Dock

- Charging: Inductive charging dock (fully charged in 4 hours)
- Input power: 5V / 2.5W
- Waterproof rating: IPX1 (not waterproof)

Safety Instructions:

- ① Operating temperature range: 0 – 50 .
- ② Charge at room temperature. Keep both the handle and charging base dry.
- ③ Store in a cool, dry place away from direct sunlight.
- ④ After each use, wipe the handle clean. Do not soak in water.
- ⑤ This toothbrush has a fully sealed waterproof design. Continuous use may cause slight heating and minor expansion of the silicone sleeve, which is normal. For best performance, let the product rest for 1–2 minutes after three consecutive uses to prevent the motor from overheating and damaging the product.
- ⑥ If not used for an extended period, charge and discharge the product once every 3 months.

- ⑦ Keep the charging base away from water. Do not immerse it or place it in dusty environments.
- ⑧ Use only a DC 5V 1A adapter (e.g., USB Type-A power supply). Avoid incompatible fast chargers, which may damage the product or pose safety risks.
- ⑨ Do not place metal objects on the charging base while charging, to avoid potential hazards.
- ⑩ If the brush head is difficult to remove, soak the connection between the brush head and the handle in hot water at 50–70°C and heat for 30 seconds, then try to remove it. Do not twist the brush head forcefully, as this may cause damage to the product.

- ⑪ This is not a toy. Children must use under adult supervision. If you have oral conditions or medical restrictions, please follow your doctor's advice.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.