


USER MANUAL

Foot Massager

MODEL NO.: MD-FM12

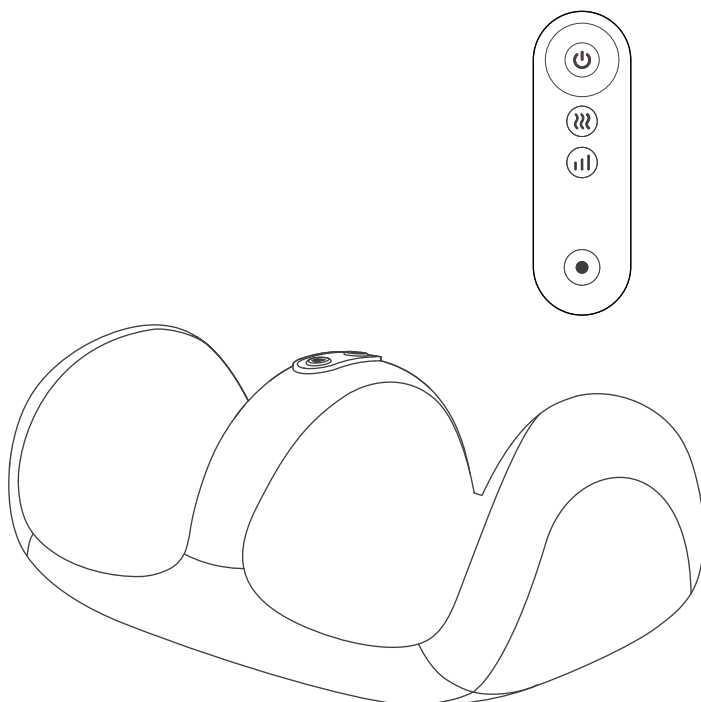

Thank You

Thank you for purchasing Foot Massager. Medcursor is committed to providing our customers with high-quality products and great services. If you have any questions or concerns about using your new product, feel free to contact our helpful customer support team at support@medcursor.com.

We hope you enjoy your new product!

Explore more surprises at
 www.medcursor.com

The following pages of this manual have the instructions for all the functions and features of your Foot Massager. Please read the user's manual carefully.

TABLE OF CONTENTS

IMPORTANT SAFEGUARDS 1

GETTING TO KNOW YOUR FOOT MASSAGER..... 3

BEFORE FIRST USE 4

OPERATION INSTRUCTIONS 5

CARE & MAINTENANCE 5

TROUBLESHOOTING 6

FCC STATEMENT 7

WARRANTY INFORMATION 8

CUSTOMER SUPPORT 8

IMPORTANT SAFEGUARDS

IMPORTANT:

When using the electrical appliance, basic safety precautions should always be followed including the following:

1. Read all instructions.
2. Close supervision is necessary by or near children.
3. To protect against electric shock, do not immerse cord, plug, or any parts of the Foot Massager in water or other liquids.
4. Do not let the cord hang over the edge of the table or counter.
5. Do not operate the appliance with a damaged cord or plug, or after the appliance malfunctions or has been damaged in any way.
6. The use of accessory attachments not recommended by the appliance manufacturer may cause hazards or injury.
7. Unplug from outlet when not in use.

MATTERS NEEDING ATTENTION:

1. This machine cannot be used in humid environments or bathrooms. Please keep it in a dry place away from strong light. The kneading massage area should not be used during sleep to prevent long-term stimulation to a certain part from affecting health.
2. If excessive stimulation occurs during use, causing discomfort, you should stop using the machine immediately.
3. When a child or a person with limited mobility uses this product, a suitable supervisor must be nearby.
4. If the power cord of the machine is damaged, please contact customer support.
5. You must not replace the damaged cord yourself.

WARNINGS:

1. This device is not intended for use by individuals (including children) with physical disabilities, sensory or neurological impairments, or a lack of relevant experience and knowledge, unless they are supervised and guided by a guardian.
2. Children must be supervised to ensure that they do not use this appliance as a toy.
3. If the surface of this machine heats up during use, extra caution must be taken when using the it for people who are sensitive to heat.

WARNING – Button Cell Battery Safety


- **INGESTION HAZARD.** This product contains a button cell or coin battery.
- **DEATH** or serious injury can occur if swallowed.
- **A** swallowed button cell or coin battery can cause internal chemical burns in as little as 2 hours.
- **KEEP** new and used batteries **OUT OF REACH OF CHILDREN.**
- Seek **immediate medical attention** if a battery is suspected to be swallowed or inserted inside any part of the body.
- **Battery Replacement**
Replace the battery with the same type and rating only.
Insert the battery with the correct polarity (+/-).
Keep the battery compartment securely closed after replacement.
Dispose of used batteries according to local regulations.



**SAVE THESE INSTRUCTIONS
HOUSEHOLD USE ONLY**

GETTING TO KNOW YOUR FOOT MASSAGER

1.PRODUCT DIAGRAM

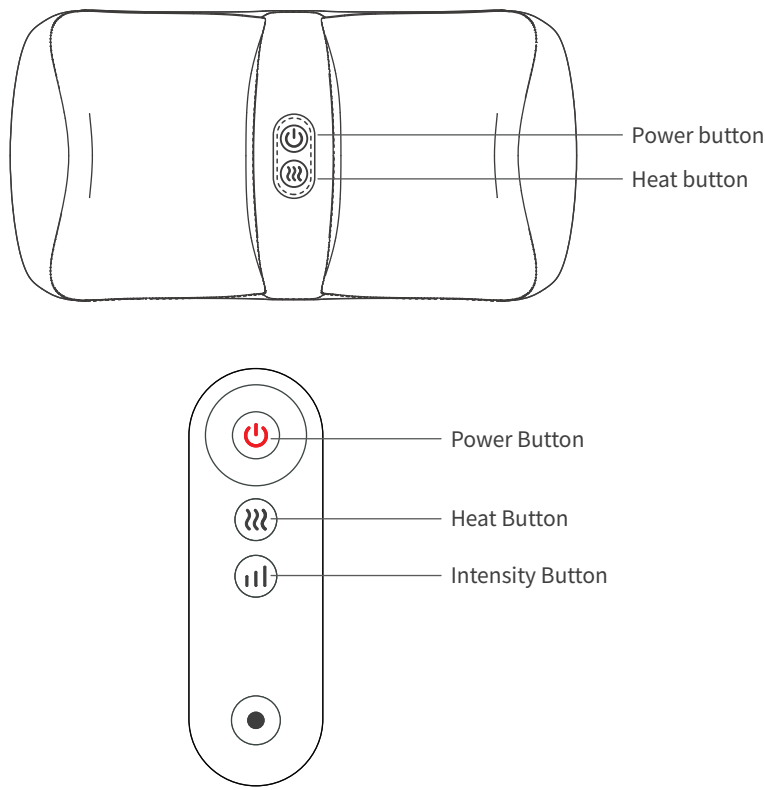


2.SPECIFICATIONS

Product Name:	FOOT MASSAGER
Product Model:	MD-FM12
Rated Voltage:	12V2A
Rated power:	24W
Website:	www.medcursor.com
Email:	support@medcursor.com

BEFORE FIRST USE

USING THE MASSAGER

1. Remove the Shiatsu Foot Massager and its accessories from the packaging box. Kindly check if all package contents are included.
2. Connect the DC power adapter to the massager and plug it into an outlet.

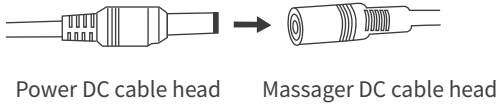

3. If you feel mild pain during the initial use of the massager, we recommend that you wear thick socks or put a towel under your feet to reduce pressure and improve your massage experience. Once you get used to the massage intensity, you can remove the towel or socks for a more comfortable and effective massage experience.

OPERATING INSTRUCTIONS

- Press the "⏻" button to turn on the machine. The default timer is set to 15 minutes, the heating function is turned on by default, and the massage intensity is set to low intensity. After turning on the machine, press the button to switch the massage intensity from low intensity to medium intensity, press again to switch to high intensity, press once more to stop the massage, and press again to return to low intensity, and so on. Long press the button to turn off the machine, and all functions will stop.
- After turning on the machine, press the "🌀" button to control the heating function. The heating function is turned on by default. Press once to turn off the heating, press again to turn the heating function back on, and so on.
- ★ The program is set to automatically change the massage head direction every 60 seconds.

REMOTE CONTROL

- Remote allows full function control of the device within operating range
- Wireless control for convenience during use

BUTTON FUNCTIONS

BUTTON	FUNCTION	OPERATION DESCRIPTION
Power	Power On / Off	Short press to turn on or off the device. When powered on, the device starts with the default settings: 15-minute timer, massage head rotating forward, heating function ON, and massage intensity set to low intensity.
Intensity	Intensity Adjustment	Short press to cycle massage intensity levels: Low Intensity → Medium Intensity → High Intensity → Stop → Low Intensity.
Heat	Heating Control	Short press to turn the heating function ON/OFF independently.

CARE & MAINTENANCE

MAINTENANCE

1. When not in use, disconnect the power supply and store the product in a safe place.
2. Avoid contact with sharp objects to prevent piercing or damaging the product.
3. Do not store this product near high temperatures, open flames, or direct sunlight, as it may cause discoloration or damage.
4. Avoid contact with benzene, diluents, or other corrosive liquids on the product and its accessories.
5. When storing, do not bind the product with wires or hang it using wires.

CLEANING

1. Before cleaning, please turn off the power supply and wipe the product gently with a soft, slightly damp sponge or cloth.
2. The product should not be soaked in any liquid and should not be washed directly with water or any other liquid.
3. Do not use sharp brushes, gasoline, kerosene, furniture cleaning fluid, etc., to clean this product.

TROUBLESHOOTING

COMMON PROBLEMS	POSSIBLE SOLUTIONS
There is a clattering noise during operation.	This is perfectly normal and is caused by the operation of the product. It is a normal condition.
The massage device switches off suddenly.	The device may have been used for 15 minutes, activating the auto-shutoff function. To continue using it, press the "power" button after resting the unit for 15 minutes.
The massager is turned on, but the massage functions stopped working.	1. Make sure the power outlet is working. 2. Unplug the power adapter from the massager and plug it back in. 3. If you have been using the massager for a long time, the overheat protection may have activated. Turn the massager off and allow it to cool down for at least 40 minutes, then turn on again.
Remote does not work	Check battery / ensure within range

FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC RADIATION EXPOSURE STATEMENT

The equipment complies with FCC Radiation exposure limit set forth for uncontrolled environment. The device shall be operated and installed without restriction.

WARRANTY INFORMATION

We offer a 1-year warranty on the product, starting from the date of purchase. Within the warranty period, we will address any defects in the appliance resulting from faults in materials or workmanship. Please contact us at support@medcursor.com to initiate a warranty claim.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original product parts are not used.

CUSTOMER SUPPORT

1. If you have any product questions or concerns, please don't hesitate to contact us directly at support@medcursor.com or visit our website at: www.medcursor.com.
2. We have a hassle-free warranty and customer service team for satisfactory solutions to any issue you might have.



Scan to Follow Us



@Medcursor
www.medcursor.com