PHILIPS

One

for Kids

HY1130


English

Important safety information


WARNINGS

MEDICAL WARNINGS


1

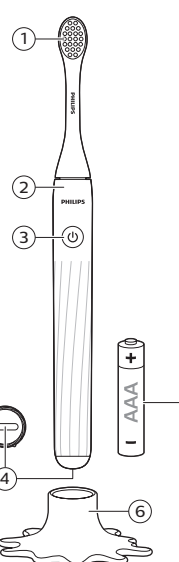
2

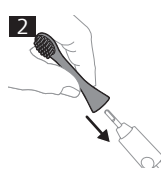
3

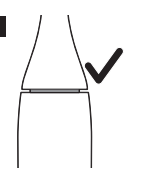
4

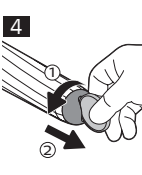5

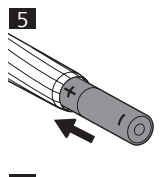
6

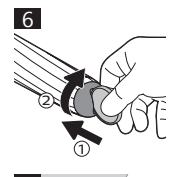
7


8

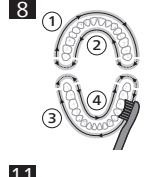
9

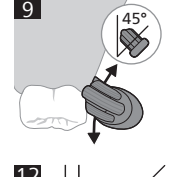
10

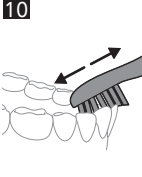11

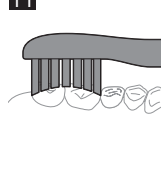12

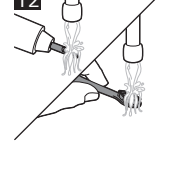Battery replacement

Your power toothbrush is powered by a replaceable AAA-battery. A single AAA battery should provide enough power for two 2-minute brushing sessions every day for 3 months. If used more often, the battery will drain before 3 months.

Note: Make sure your hands and the product are dry when you insert the battery.

- 1 Insert a coin into the coin slot at the bottom of the handle and turn counterclockwise (Fig. 4). Gently shake the battery out of the handle.
- 2 Insert a fresh AAA battery into the battery compartment (Fig. 5). Make sure the + and - poles of the battery point in the right direction.
- 3 Place the removable end cap back on the handle and use a coin to turn clockwise until it is secure and the lock indicator markings on the handle and end cap line up (Fig. 6).

Note: If the indication on the handle and end cap do not line up, it may result in your toothbrush not working.

Caution: The battery replacing operation should only be conducted by an adult. Keep the end cap, battery, and coin far away from children.

Using your Philips One for Kids toothbrush

If you are using this toothbrush for the first time, it is normal to feel more vibration than when using a manual toothbrush. The Philips One for Kids uses microvibrations to help gently clean your teeth. Please follow the brushing instructions below for the best experience.

Philips One for Kids comes with a SmartTimer and QuadPacer to help you brush evenly throughout your mouth.

Note: This product, including the toothbrush stand, is not a toy.

SmartTimer

The SmartTimer indicates that your 2-minute brushing cycle is complete by automatically turning off the toothbrush at the end of 2 minutes.

Note: Dental professionals recommend brushing no less than 2 minutes twice a day.

QuadPacer

The QuadPacer is an interval timer that beeps every 30 seconds to remind you to brush all teeth evenly by dividing your mouth into 4 sections.

Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

Note: If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the SmartTimer and QuadPacer resets.

Brushing instructions

Please follow the brushing instructions for the best experience.

Note: For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.
- 3 Press the power button to turn brush on.
- Note:** Keep the center of the brush in contact with the teeth all the times.
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with back and forth strokes so that the bristles reach all surfaces of the teeth. Change to another section when the QuadPacer beeps. Continue this motion throughout the 2-minute brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

- 5 To clean the inside surfaces of the front teeth, tilt the brush vertically and make several up-and-down strokes.
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your power toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Battery status

This product does not come with a battery status indicator. As the battery gets weaker, you may experience a drop in vibration power. This is normal. While the microvibrations are still effective, you can increase the power by replacing your battery.

Cleaning

You should clean your toothbrush and accessories regularly to remove toothpaste and other residue. Failure to clean your product may result in an unhygienic product and damage to it may occur.

Warning: Do not clean product or accessories in dishwasher, microwave, with chemicals or with boiling hot water.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.
- 2 Wipe the entire surface of the handle with a damp cloth. Dry after each use.

Note: Do not tap the handle on the sink to remove excess water.

Cleaning brush head and toothbrush stand

- 1 Rinse the brush head and bristles and dry after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the toothbrush stand weekly with warm water and wipe clean with a damp cloth.

Storage

If you are not going to use the product for a long time, clean it (see 'Cleaning') and store it in a cool and dry place away from direct sunlight.

Warning: If product will not be used on a regular basis, remove the battery to prevent possible damage due to chemical leakage. If battery does leak, remove it carefully. Do not allow bare skin to touch leaking fluid.

Replacement

Brush head

Replace the brush heads every 3 months to achieve optimal results. Only Philips One for Kids brush heads are compatible with the Philips One for Kids toothbrush.

Locating the model number

Remove brush head from handle. Look at the top of the handle near the metal shaft for the model number (HY111X).

Support and Limited Warranty

Philips provides a limited warranty for its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided.

Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: www.philips.com/support.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Disposal

Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit www.recycle.philips.com.

This product contains batteries:

- Dispose of batteries properly. Do not incinerate. Batteries may explode if overheated.
- Do not wrap in metal or aluminum foil. Tape the waste battery terminals before discarding.
- It is suggested that you contact your local town or city to determine proper battery redemption site(s) in your area.
- To remove the battery, follow the instructions in the battery replacement chapter.

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical and electronic products. Correct disposal helps prevent negative consequences for the environment and human health.

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