


SMART EYE MASSAGER

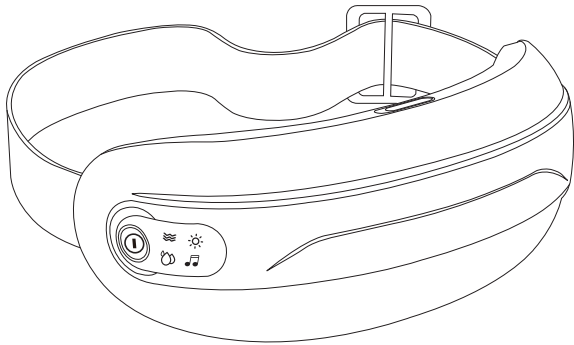

INSTRUCTION MANUAL

—Please read this manual carefully before using this product and keep it for future reference—

Instructions Before Use

- When turning on the spray function, please make sure that there is enough water in the tank to avoid dry burning and consequently causing instrument failure.
- If there is any discomfort in the skin or eyes during use, please stop using immediately.
- This product is neither a substitute for medical equipment nor for therapeutic use.
- Please adjust the appropriate length and frequency of use according to your own state.
- Do not wear contact lenses during use.
- Please do not use this product with makeup on, or your makeup may run.
- Please do not use the product in the following situations:
 - When feeling sick or having a fever;
 - There are wounds, eczema, or swelling on the skin;
 - Skin abnormalities such as dermatitis or excessive sunburns
- The following people should not use this product:
 - Those with allergies or facial skin sensitivities;
 - Those with eye diseases, eye inflammation, swelling, injury, eczema, cataract, glaucoma, retinal detachment and other abnormalities;
 - Those with facial nerve disorders or skin sensory disorders;
 - Those with high myopia (greater than -6 D);
 - Pregnant women;
 - Those who have or have had eye diseases, or had eye surgery resulting in temperature sensitivity or retardation, or are undergoing eye treatments.

Our company is not be liable for any problems due to failure to follow the Instruction Manual.

Cautions for Use

* Frequency of use

Use as needed with 15 minutes per use. It is advised to rest for a while after one use before using it again.

* Cautions Before Use

When using this product, please adjust to your personal comfortable position to allow your body to completely relax. Stop using if your skin becomes sore or you feel uncomfortable.


* Cleaning of the Product

- When cleaning the product, please make sure that it is turned off without charging;
- If its appearance needs to be cleaned, use a dry cloth with water to wipe the polished surface and then let it dry naturally;
- The main unit cannot be washed with water. Do not use water to rinse its body. The spraying port and water tank inlet can be wiped gently with a cotton swab.
- When the product is not to be used for a long time, clean and dry its body and components and place them in a cool and dry place. Avoid direct sunlight and keep away from heating devices and open flames.

* Status after Use

After using, empty the liquid in the tank and keep the tank clean. Avoid leaving the tank for a long time when there is liquid in it.

Product Introduction

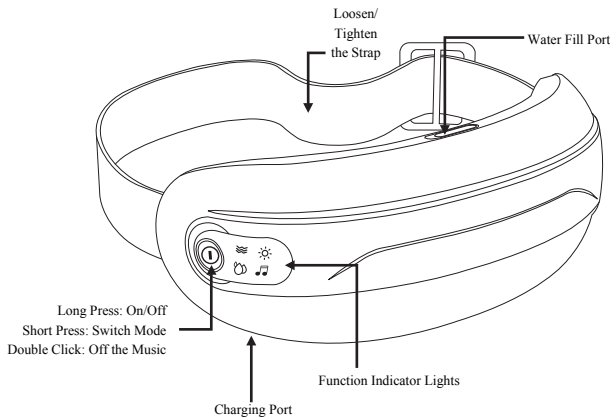

Product Parameters

- [Product Name] Smart Eye Massager
- [Product Model] W6
- [Rated Voltage] 5V
- [Rated Power] 5W
- [Smart Timer] 15min
- [Massage Mode] Vibration massage, eye moisturizing spray, and eye warming thermal mist

Operating Instructions

- [Water Filling]
- Add mineral water or tap water into the complimentary small mouth bottle;
 - Open the soft rubber plug at the filling port of the water tank upwards, fill with the right amount of mineral water or tap water with the small mouth bottle, and then fasten the plug.

Note: The water tank should not be overfilled each time.

Mode Description

- Energizing Mode: Thermal mist + Eye moisturizing spray (large amount) + Vibration massage + Soothing music
- Eye Exercises Mode: Vibration massage + Soothing music
- Soothing Mode: Thermal mist + Eye moisturizing spray (continuous spray) + Soothing music
- Sleep Mode: Vibration massage + Soothing music

Steps for Use

- Adjust the strap to the appropriate position, and wear the massager correctly to start using;
- Note: If there is too little water in the tank, to avoid dry burning caused by instrument failure, do not lie down to use;
- Long press ① for 2 seconds, and the instrument will speak “Welcome to use this eye massager, and please select the mode”, and then press the button ① again to enter the energizing mode. Short press ① to enter the next mode. There are 4 modes that can be switched (for details, please refer to the Mode Description);
 - The massager will be automatically turned off after 15 minutes of use. Please remove the residual water in the tank and dry the instrument, and then place it in a cool and dry place.

*** Press the button ① twice to turn off/on the music.**

*** If you need to stop during using, long press the button ① for 2 seconds to turn off the instrument.**

Note: Bluetooth will be automatically turned on after the host is started, and you can just open your cell phone's Bluetooth and search for the Bluetooth name “YBY” to connect.

Charging Description

When the instrument speaks in a low battery tone of “Low battery, please charge”, please charge in time to avoid affecting the normal use.

This product uses rechargeable lithium batteries, and the indicator light will show “Red” when charging, and “green” when the charging is completed. Please use the power adapter meeting the rated voltage (5V=1A) of this product to charge, and do not turn on the power when charging and do not put it near open flame or heat it.

Precautions

- Store the product out of the reach of infants and children.
- Do not use the product in the bathroom.
- Do not modify or repair the product privately.
- If any of the following malfunctions occurs, please stop using immediately and disconnect the power supply:
 - The instrument, adapter or charging cable is deformed or abnormally hot;
 - The instrument, adapter or charging cable emits a burning odor;
 - The instrument, adapter or charging cable makes a strange noise during using or charging

Failures	Possible Causes	Solutions
Unable to power on	Low battery level	Please charge promptly
Unable to charge	Incorrect adapter	Check if the adapter parameters match the product requirements
Unable to spray	Calcium substances in the water clog the mesh holes of the misting tablets	1. Place 4-5 drops of household white vinegar in the tank 2. Turn the instrument on for about 2 minutes 3. Turn off the instrument, pour out the water in the tank and replace it with mineral water
		Wash it with clean water and turn it on again
	There are bubbles in the water	Shake it up or restart it

FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.