beurer FT 09

ENGLISH Digital thermometer Instruction manual

Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.

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ENGLISH

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer, Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested

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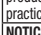
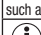
Make them accessible to other users and observe the information they contain.

With kind regards, Your Beurer team

Signs and symbols Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.
READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	- FOR HOUSEHOLD USE ONLY -
	- This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
	WARNING – Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
	CAUTION – Indicates a potentially hazardous situation, which, if not avoided, may result in minor/moderate injury or product/property damage. It also alerts against unsafe practices.
	NOTICE – Addresses practices not related to personal injury, such as product and/or property damage.
	Note Note on important information.
	Follow instructions for use
	Type BF applied part
IP 27	Protected against foreign objects equal to or greater than 12.5 mm and against temporary immersion in water.
	Permissible storage and transport temperature and humidity

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	Permissible operating temperature and humidity
---	--

WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- Keep the product out of the reach of children and pets. This product is not intended for use by children.

CAUTION – To reduce the risk of personal injury or product/ property damage:

- This product is intended solely for household, non-commercial, non-industrial use in measuring body temperature.
- Do not leave the product unattended while it is in use, especially when children are present. The product should only be used by an adult.
- Before each use, check the thermometer for signs of damage or wear. Damaged or worn thermometers should not be used.
- Always wait for the acoustic signal which indicates the minimum time required for measurement. Discuss the temperature with your family physician.
- If the temperature is very abnormal, please contact your family doctor promptly.
- When the thermometer is switched on, it carries out a self-test. Therefore the accuracy of measurement need not be tested at regular intervals.
- Do not walk, run, or talk while measuring a temperature.
- Do not drop the product or expose it to strong shocks.
- Avoid exposing the product to dust, high humidity, or extreme temperatures.
- Replace the battery when the Low Battery symbol shows on the display, or inaccurate readings could result.

NOTICE – To reduce the risk of product/property damage:

- Do not attempt to disassemble the product, except to replace the battery.
- Do not allow the product to be bitten or bent.
- Using the thermometer in strong electromagnetic fields such as next to a mobile phone, can cause a malfunction.
- Clean the product after each use following the instructions provided under Care and Maintenance.

- Do not use or store the thermometer outside the manufacturer's stated temperature and humidity range.
- Portable and mobile HF communication systems may interfere with this device.

- Notes on electromagnetic compatibility**
- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the other devices in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Remove discharged batteries from the product and dispose/ recycle in compliance with all applicable laws.
- If your skin or eyes come into contact with battery fluid, flush out the affected areas with water and seek medical assistance.
- Do not charge or short-circuit batteries.

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- A waiting period of 20 seconds between testing is recommended to ensure complete accuracy.

Rectal Temperature Method

- Apply a water-soluble lubricant to the Probe Tip. Do not use petroleum jelly.
- Gently insert the Probe Tip into the rectum approximately ½ inch (1 cm). Normal rectal temperature is 36.6°C (97.8 °F). **NOTE:** Do not use a thermometer that has been used for rectal temperature measuring for armpit or oral temperature measuring.

Oral Temperature Method

- CAUTION:** Do not talk, eat, drink, bathe, exercise, or smoke before a temperature measurement.
- Place the Probe Tip under the tongue as near as possible to a heat pocket, and keep the mouth closed (see Fig. 1).

Armpit Temperature Method

CAUTION: Wipe the armpit with a dry towel before placing the Probe Tip under the arm.

- Place the Probe Tip under the arm, positioning the thermometer vertically.
- Position the arm across the chest, causing the Probe Tip to be completely covered by the arm.
- Measuring time is approximately three minutes, depending on the placement of the Probe Tip and the constant closure of the arm against the body.
- The measured temperature will be lower than an oral measurement.

Mode of measurement	Recommended time required for measurement	Mean divergence to the rectal method
in the anus	until the acoustic signal	–
in the oral cavity	until the acoustic signal	-0.4 °C to -1.5 °C (31.2 °F to 29.3 °F)
in the armpit	minimum 5 minutes	-0.7 °C to -2.0 °C (30.7 °F to 28.4 °F)

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Care and Maintenance

- Clean the thermometer after each use with a soft cloth and isopropyl alcohol diluted with water. Wipe again with a clean, water-dampened cloth, and let dry.

NOTICE: Do not soak the thermometer or sterilize by boiling or steaming it.

- Store the thermometer where it will not be exposed to direct sunlight, dust, humidity, or extreme temperatures.

Battery Replacement

CAUTION: Do not replace the battery in the presence of small children or pets, or during conditions of extreme heat.

- The LCD Display will flash when the battery needs replacing.
- Pull the Battery Cover away from the base.
- Use a non-metal pointed tool to gently pry loose the old battery.
- Insert a new 1.5V LR41 alkaline battery with the positive (+) side up.
- Replace the Battery Cover.

WARNING: This thermometer contains a button cell battery containing Mercury. Do not dispose of the used battery in municipal waste or in fire. Dispose or recycle discharged batteries in compliance with applicable laws.

Troubleshooting Guide

Error Message	Problem	Solution
ERR	The system is not functioning properly.	Unload the battery, wait for 1 minute and repower it. If the message reappears, contact the retailer for service.
Hi	Temperature is above 42.9° C (109.9 °F).	Turn off, wait one minute and take a new temperature via close contact and sufficient rest.

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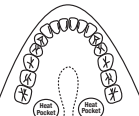
Fig. 1

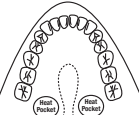
Fig. 1

Error Message	Problem	Solution
⚠	Temperature below 32° C (90° F).	Turn off, wait one minute and take a new temperature via close contact and sufficient rest.
🔋	Battery low or depleted.	Replace battery.

Type	FT 09
Range of measurement	32 °C to 42.9 °C (89.6 °F to 109.2 °F)
Accuracy of measurement	± 0.1 °C (32.1 °F) in a water bath with a temperature of 35.5 °C to 42.0 °C (95.9 °F to 107.6 °F) <p>± 0.2 °C (32.3 °F) in a water bath with a temperature of 32.0 °C to 35.4 °C (89.6 °F to 95.7 °F)</p> ± 0.2 °C (32.3 °F) in a water bath with a temperature of 42.1 °C to 42.9 °C (107.7 °F to 109.2 °F)
Ambient temperature	+10 °C – +40 °C (+50 °F - +104 °F), with 30 % to 85 % relative humidity
Storage temperature	-10 °C – +60 °C (+14 °F - +140 °F), with 25 % to 90 % relative humidity

This device meets the requirements of the EU Directive 93/42/EEC concerning medical devices, as well as those of the Medizinproduktegesetz (German Medical Devices Act), the European Standard EN 12470-3: Clinical thermometers – part 3: Performance of compact electrical thermometers (non-predictive and predictive) with maximum device and the European Standard EN60601-1-2 (In accordance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3 and IEC 61000-4-8), and is subject to particular precautions with regard to electromagnetic compatibility. Please note that portable and mobile HF communication systems may interfere with this unit.

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This appliance conforms to the following standards: ASTM E1112 Standard Specification for Electronic Thermometer for Intermittent Determination of Patient Temperature, ISO 80601-2-56 Medical electrical equipment – Part 2-56: Particular requirements for basic safety and essential performance of clinical thermometers for body temperature measurement, IEC 60601-1-11 Medical electrical equipment – Part 1-11: General requirements for basic safety and medical electrical systems used in the home healthcare environment and complies with the requirements of IEC 60601-1-2 (EMC), AAMI/ANSI ES60601-1 (Safety) standards. And the manufacturer is ISO 13485 certified.

Warranty
Limited Lifetime Warranty For Original Purchaser
Your Beurer Digital Thermometer, model FT 09, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Digital Thermometer, model FT 09, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory;

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alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer’s recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or

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limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

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Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

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This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: www.beurer.com

For Canada only:
Manufacturer:
Beurer GmbH
Söflinger Str. 218
89077 Ulm Germany

Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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