
UREVO Wireless Recovery Boots

User Manual

Mode.UCRM003

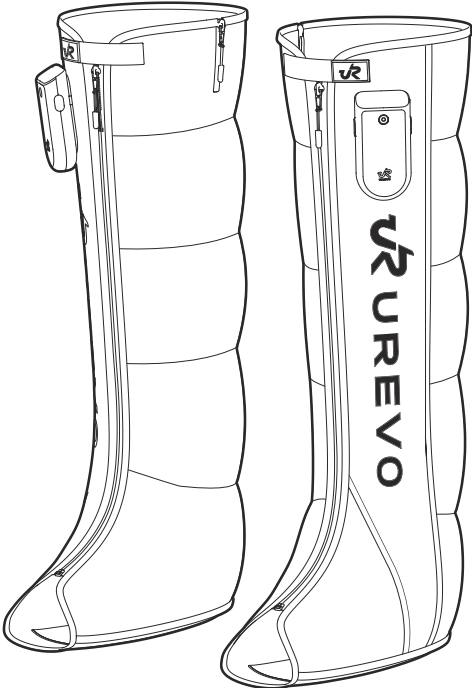IMPORTANT: Read all instructions carefully before using this product. Retain this owner's manual for future reference. The specifications of this product may vary from this photo, subject to change without notice.

UREVO WELLNESS US CO., LTD
22635 NE MARKETPLACE DR STE 100, REDMOND, WA,
98053-5898, UNITED STATES

Version: 1.0 FCC ID: 2BQ4M-UCRM003

MADE IN P.R.C.

Activate UREVO App Now
For Longer Warranty &
Smart Features Access

Table of Contents

01	EN	IMPORTANT SAFETY INSTRUCTIONS
03		CUSTOMER SERVICE
03		WARRANTY
04		UREVO APP
06		SPECIFICATIONS
06		PACKING LIST
07		PRODUCT INTRODUCTION
10		HOW TO USE
13		ERROR CODE
13		CLEANING AND MAINTENANCE
13		ACCESSING PRODUCT INFORMATION

Important Safety Instructions

WARNING!

Read all instructions before using (this appliance).

DANGER – To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
2. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
6. Never drop or insert any object into any opening.
7. Do not use outdoors.
8. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
9. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
10. Do not stand on or in appliance. Use only while seated.
11. This appliance contains batteries that are non-replaceable.
12. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
13. Do not immerse this appliance in water or any liquid.
14. Do not use this appliance in humid environments such as bathrooms and saunas.

SAVE THESE INSTRUCTIONS

	Note: There is a possibility of personal injury and damage when used incorrectly		Handle fragile items with care
	Please read the instructions carefully before using this product		Moisture-proof and rainproof
	Do not throw rubbish carelessly. Please dispose of it according to local regulations		

1. Charging with adapters that do not meet the requirements of this product is prohibited.
2. Non-professional repair technicians are prohibited from decomposing and repairing this product.
3. If the power cord is damaged, in order to avoid danger, it must be replaced by the manufacturer or its maintenance department and similar full-time personnel.
4. Do not use in wet areas or when the product has penetrated into water or other liquids.
5. After use, please place the massager in a dry place, away from humidity and direct sunlight.
6. Do not place in fusible, explosive, chemical solution, damp place.
7. Do not touch the product with sharp objects to prevent damage to the product parts.
8. If you feel uncomfortable or abnormal during use, please stop using immediately.
9. Do not use in high temperature environment.
10. Please keep the product clean and use a dry towel during maintenance.
11. Before using this product, it is recommended to consult a doctor or use under the guidance of a professional:
 - High-risk pregnant women and children should not use.
 - Mental patient
 - Suffer from skin infection, skin swelling, scalding burns, etc
 - High-risk patients who have a serious illness and are not suitable for massage
 - The injured have weak or damaging bones
 - People with cardiac pacemakers
 - The body is fitted with other metal devices, such as prosthetics
 - People with serious heart disease or abnormal blood pressure
 - People with spinal abnormalities or spinal curvature disorders

SAVE THESE INSTRUCTIONS

FCC STATEMENT:

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
 - (2) This device must accept any interference received, including interference that may cause undesired operation.
- (3) Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with RSS247 of Industry Canada. Cet appareil se conforme à RSS247 de Canada d'Industrie. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device. appareils radio exempts de licence. Son fonctionnement est sujet aux

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Radiation Exposure Statement: The product complies with the Canada portable RF exposure limit set forth for an uncontrolled environment and are safe for intended operation as described in this manual.

The further RF exposure reduction can be achieved if the product can be kept as far as possible from the user body or set the device to lower output power if such function is available.

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This Class B digital apparatus complies with Canadian ICES-003.

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Hours of Operation:

Monday - Friday (9 am - 12 pm; 1 pm - 6 pm PST)

For more content, please visit : urevo.com

Warranty

You'll receive a **12-month** warranty and lifetime customer service after activation. This warranty applies only to eligible purchases from UREVO authorized distributors.

For more details, visit urevo.com/warranty.

UREVO APP

Discover exciting new features and massage modes by connecting to the UREVO App.

App Connection Steps:

Turn on Bluetooth on your phone/tablet.

Open the UREVO App and follow the on-screen instructions to add and connect your UREVO Wireless Recovery Boots Via Bluetooth.

Activate Your UREVO App for Extended Warranty & Smart Feature Access!

The UREVO App offers:

- Smart control for an immersive relaxation experience.
- Personalized smart massage intensity, tailored just for you.
- A diverse range of massage modes to suit various life scenarios.
- Intelligent data reports for visible health benefits.



OTA Firmware Update Instructions

To Update Your Firmware:

1. Navigate to the Settings page within the massager app.
2. Tap on the [Firmware Update] option.
3. Select [Update] to begin the firmware update process.

Important Notes During Upgrade:

- The appliance will be temporarily unavailable (estimated 5-10 minutes).
- Do not power off your phone or the massager.
- Ensure your Bluetooth connection remains active throughout the entire process.

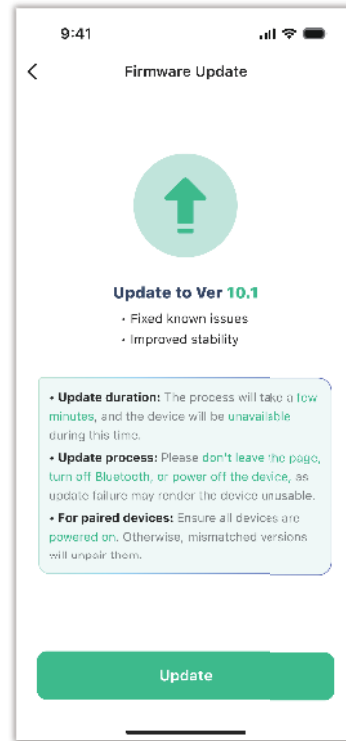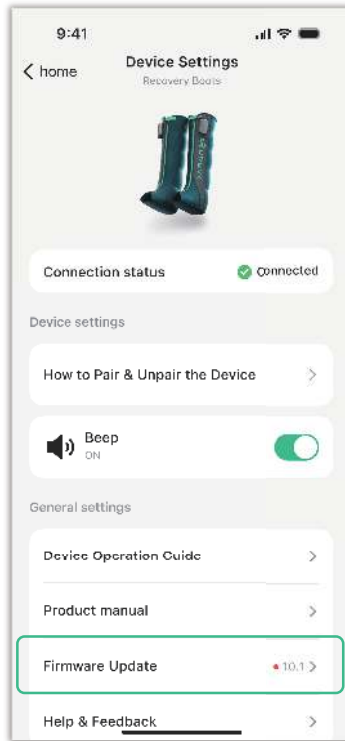
Critical Requirements for OTA Update:

- Ensure battery level is above 25%.
- Keep the appliance powered ON throughout the entire procedure.

Specifications

Product configuration	Detailed parameters
Product name	UREVO Wireless Recovery Boots
Product model	UCRM003
Modes	4 modes
Pressure	Low\Mid\High
Decibel rating	<65dB
Battery	5000mAh
Charging time	<3H
Working time	>4H
Operating Conditions	Temperature 5°C to 35°C/ 33°F to 95°F
Product Weight	4.2 kg / 9.24lb
Input power	20W
Output Power	20W
Charging port	Type C



Packing List

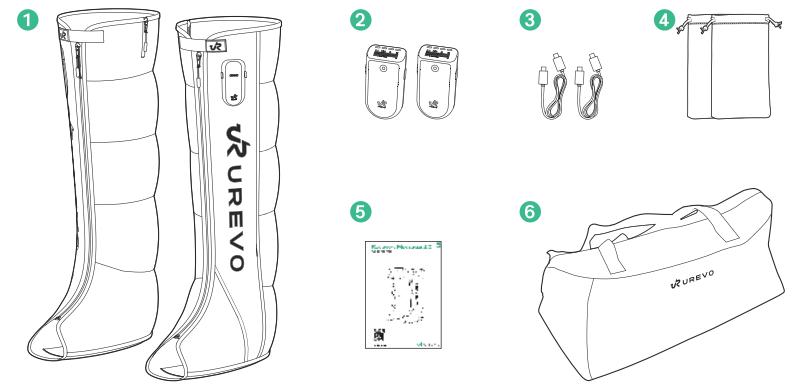- 1 UREVO Wireless Recovery Boots x 1
- 2 Controller x 2
- 3 Type-C to Type-C fast charging cable x 2

- 4 Flannel bag x 2
- 5 User Manual x 1
- 6 Storage Bag x 1

Product Introduction

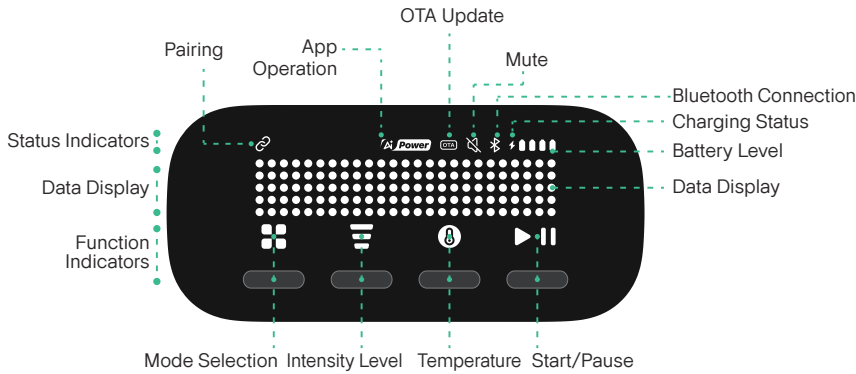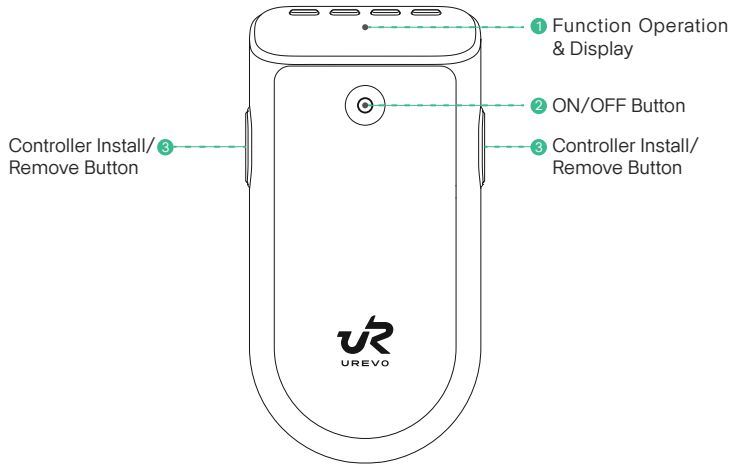a. Icon Description

Pairing Icon

Default: Solid (indicating both leg units are paired).
 Status: Lit when paired, off when unpaired.
 Note: Two left leg units or two right leg units cannot be paired together.

App Icon

Default: Unlit.
 Status: Lights up when connected to the app and the Smart Massager function is active.

OTA Icon

Default: Unlit.
 Status: Flashes during the Over-the-Air (OTA) firmware update process.

Mute Icon

Buzzer Default: On (icon off).
 Sounds: A tone sounds for abnormal codes during power switching, pairing completion, OTA, massage end, or when the controller is properly inserted. There's no sound for function selection or massage start/stop.
 Disabling Buzzer: Can be disabled via the app (icon remains lit).
 Reactivation: If reactivated, the icon will turn off as soon as the "beep" sound is heard (due to the appliance's memory function).

Bluetooth Icon

Unconnected: Flashes (on for 0.4s, off for 0.2s). If connected to the app within 2 minutes, it stays lit; otherwise, it turns off.
 Connected: Remains steadily lit.
 After Disconnection: Flashes (on for 0.4s, off for 0.2s). If reconnected to the app within 2 minutes, it stays lit; otherwise, it turns off.

Charging Icon

Status: Remains steadily lit during charging (both self-charging and reverse charging). Turns off when not charging.
Note: The product cannot be operated during charging or discharging.

Battery Level

0% - 10%: Leftmost LED segment flashes (0.5s on / 0.5s off).
 10% - 25%: One segment steadily lit.
 25% - 50%: Two segments steadily lit.
 50% - 75%: Three segments steadily lit.
 75% - 100%: Four segments steadily lit.


b. Description of buttons and corresponding icons

	ON/OFF button	<p>Power On/Off To Turn On: Press and hold the button for 2 seconds. Accidental Power Off Prevention: If you continue pressing the button for 3 seconds or more without releasing, the appliance will turn off. To Turn Off (when on): Press and hold the button for 2 seconds. Indicator Light Colors: White, Red, Orange, Green. Normal Operation: Remains White. Charging Status: Switches from Red (low charge) to Orange (mid-charge) then to Green (fully charged).</p>
	Mode Selection	<p>Number of Modes: 4 in total. Default Mode: Mode 1. Indicator: Four blocks light up sequentially based on the selected mode: bottom-left, bottom-right, top-left, top-right.</p>
	Intensity Selection	<p>Number of Levels: 3 in total: LOW, MID, HIGH. Default Level: Level A (LOW intensity). Pressure Values: LOW (80 mmHg), MID (140 mmHg), HIGH (180 mmHg). Ready Status: Remains at LOW level, and its corresponding indicator light is always on.</p>
	Temperature Selection	<p>Number of Levels: 4 in total: Off, 38°C (100°F), 43°C (109°F), 48°C (118°F). Unit Switching: The app supports switching between °F and °C. Indicator Light: Semi-lit under "Ready" status (default OFF), and fully lit at active temperatures.</p>
	Start Massage	<p>Ready Status: The Start icon is always on. The Pause icon is off.</p>
	Pause Massage	<p>In Operation Status: The Start icon is off. The Pause icon is always on. To End Program Quickly: Press and hold the Start/Pause button for 2 seconds.</p>


How to Use

1. Normal use


Step 1: Wear the appliance according to the instructions in the picture.

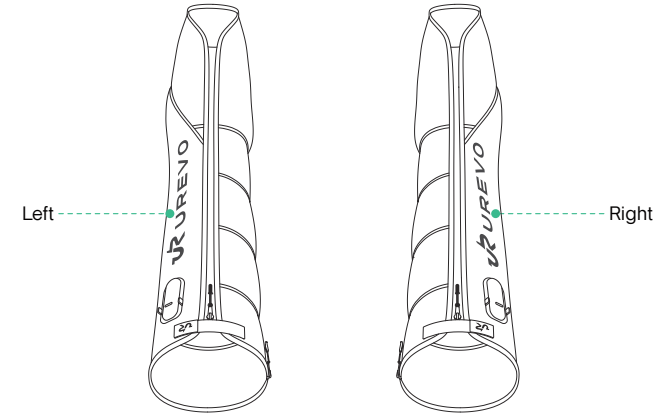
Step 2: Adjust the width and length adjustment zippers accordingly and fasten the straps.

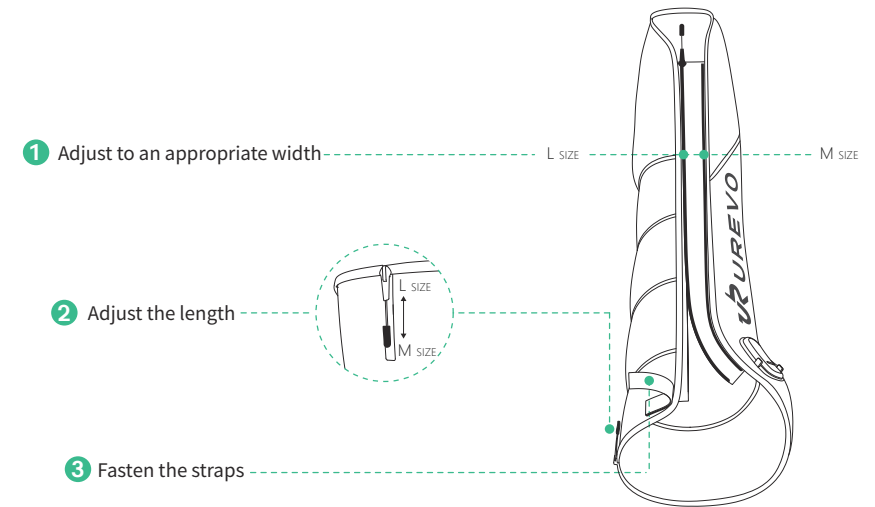Step 3: Install the controller, then press and hold the power button for 2 seconds to turn on the appliance.

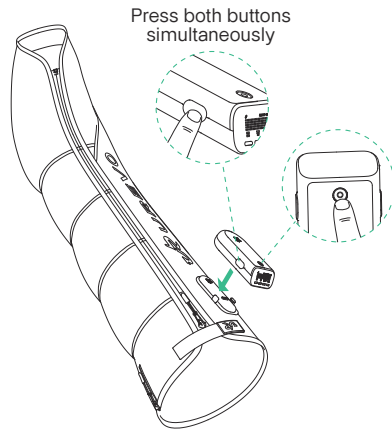Step 4: Select the program, massage pressure, and bottom heating temperature (OFF by default) from left to right, and then press the Start button.

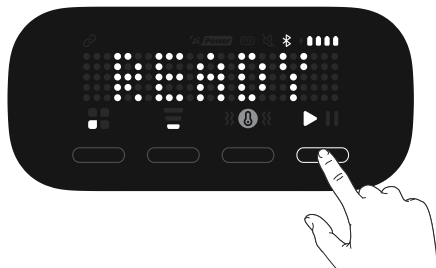
2. Targeted Relief: Localized Intensive Massage (The massage heads and straps are sold separately.)

How to Apply:


1. Select a suitable massage head featuring protrusions.
2. Apply the massage head directly to your target area or secure it with the provided strap.
3. Begin the massage at the LOWEST intensity setting.

Important Safety Warning: If you experience any discomfort during the massage, discontinue use immediately.

a. The massage head is adhered to the base unit

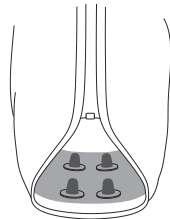
b. The massage head is attached to the securing strap

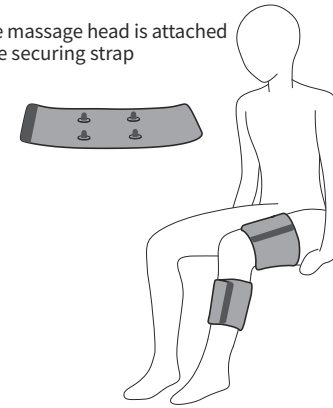3. Additional product functions


Pairing / Unpairing

In the "Ready" state: Press the Intensity Level button on both units at the same time for 3 seconds to pair or unpair.

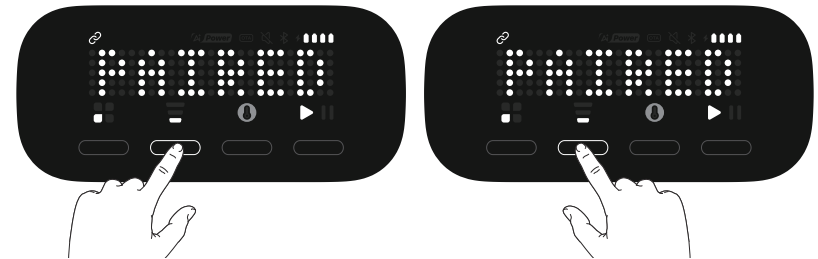Power Bank

Provides up to 20W output power for charging external electrical devices.

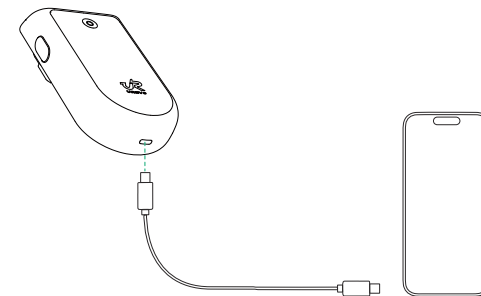Error Code

Error Code	Resolutions	Recommended Action
E02	To ensure proper function, please reinstall the controller, then turn off and restart the appliance.	Contact customer service.
E03	Disable heating function, , then turn off and restart the appliance.	
E04	Turn off and restart the appliance.	
E05	To ensure proper function, please reinstall the controller, then turn off and restart the appliance.	
E06	To ensure proper function, please reinstall the controller, then turn off and restart the appliance.	
E07	/	
E08	To ensure proper function, please reinstall the controller, then turn off and restart the appliance.	
W01	1. Inconsistent firmware versions of paired products. The hardware cannot operate until you unpair them. 2. Please upgrade the firmware in the APP.	

Cleaning and Maintenance

Before Cleaning: Ensure the product is turned off and disconnected from the power supply.

Cleaning Steps: Simply wipe the surface of the appliance with a clean cloth and water.

After Cleaning: Allow the surface to completely dry for at least 10 minutes before using the appliance.

Important Notes: Do not use irritating liquids such as gasoline or solvents to clean the appliance, as this may cause malfunction, damage to parts, or discoloration.

Ensure no liquid enters the appliance.

Accessing Product Information

To view your product's basic information, follow these steps:

1.While the device is charging, insert the charging cable.

2.Press and hold the Mode button for 5 seconds.

The screen will then display the product version and controller code alternately.

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