

Product Details

- Two Shiatsu Zones – Shiatsu zones have three finger-like rollers on each side allowing for a total of 6 massage nodes; Get a massage from various heights as well as dual rotation from all angles
- Finger-Like Kneading – Rollers with dual rotation at various heights combined with finger like kneading allows for a unique massage experience, comparable to the feeling of getting a hand massage from someone else
- Whole Body Massage – This versatile heated whole body massaging wrap is your solution to easing pain in your neck, shoulders, back, legs, or calves; Add a wrap to your existing pain management to promote faster healing
- Find Your Comfort – We recommend leaning back into a chair or sofa and letting the massager do the work; Using the handled straps allows you to massage at just the right spot with customizable pressure
- Massage for Health – Massaging promotes health, adding massage to help with pain management for injury recovery or chronic pain can help decrease inflammation, help relieve pain and discomfort, relieve muscle tension, tightness, and stress
- **Power Source:** Corded Electric
- **Item Weigh:** 1.96 Pounds (0.88 kg)
- **Product Dimensions:** 14.5"L x 11.25"W x 3.4"H (36.8 x 28.6 x 8.6 cm)

Wahl Deep-Knead Heated Massage Wrap Quickstart Guide

Wahl's Shiatsu Wrap is the Most Versatile Solution to Ease your pain. All Body massage for all areas such as neck, shoulder, back, legs and calves. Two-Shiatsu Zones with 3 Rollers each allowing 6 Massage Nodes of Various Height and Dual Rotation to Massage from All Angles. The Unique Rollers provide a True Finger-like Kneading. You can customize pressure by using the Handled Straps to massage at just the right spot. You can also simply lean back in a chair to adjust intensity as well. Made from super soft flexible fabric for comfort while using. Comes with a wall power adapter that can be used corded for continues use. Relive your pain & fatigue with Soothing Heat to Penetrate muscles for faster relief. Product in packaging may not be exactly the same as photo due to continuous improvements to the product.



Insert adaptor jack into the plug on the left side of the unit.



Plug into a AC 120-volt outlet.



To turn the massage rollers on, press the power button on the controller once.



This unit has a dual rotation feature. Press the power button again to make the massage rollers rotate in the opposite direction.




To turn the heat on, press the heat button once. To turn the heat off, press the heat button again.



To turn the massage rollers off, press the power button a third time. If the heat is on, it will remain on if you turn the rollers off. This unit automatically turns off after 20 minutes.

Wahl Deep-Knead Heated Massage Wrap Key Features

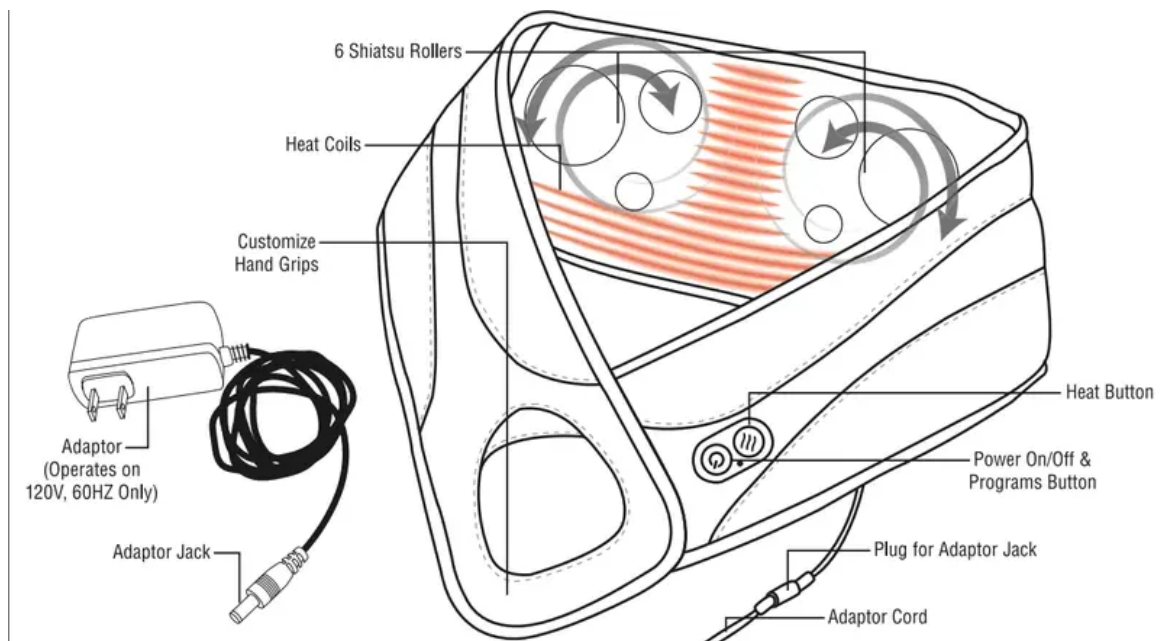

FAQ - Frequently Asked Questions

Question:Has anyone unzipped this for cleaning? How do you get it to unzip?

Answer:Have not unzipped

Question:Is it noisy? Trying to find a quiet one

Answer:It's a little bit noisy..and I don't like it very much

Question:How do you know when it is fully charged

Answer:I think it has a light. But it's kind of crappy, I'd send it back and find something better. I'm not a fan, wish I hadn't been too lazy to return mine.

Question:Does it have an adapter to be use in a car, van, truck?

Answer:You might be able to rig something, but it would be uncomfortable to use in a vehicle seat or supine in bed.

Question:How is this wrap heated? Electric or battery? Are all your products electric

Answer:This one is electric

Question:How do you clean this item?

Answer:Good but the massager could be a little more stronger.

Question:Length of cord?

Answer:By memory about 3ft

Question: Has anyone unzipped this for cleaning? How do you get it to unzip?

Answer: Have not unzipped

Question: Is it noisy? Trying to find a quiet one

Answer: It's a little bit noisy..and I don't like it very much

Question: How do you know when it is fully charged

Answer: I think it has a light. But it's kind of crappy, I'd send it back and find something better. I'm not a fan, wish I hadn't been too lazy to return mine.

Question: How is this wrap heated? Electric or battery? Are all your products electric

Answer: This one is electric

Question: Does it have an adapter to be use in a car, van, truck?

Answer: You might be able to rig something, but it would be uncomfortable to use in a vehicle seat or supine in bed.

IMPORTANT SAFEGUARDS

When using your electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS AND SAFEGUARDS BEFORE USING THIS MASSAGER.

DANGER

To reduce the risk of death or injury by electric shock:

1. Do not reach for an appliance that has fallen into water.
Unplug immediately.
2. Do not use while bathing or in the shower.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
4. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
5. Never use pins or other metallic fasteners with this appliance.
6. Carefully examine the covering before each use.
Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
7. Keep Dry – Do not operate in a wet or moist condition.

WARNING

To reduce the risk of burns, electrical shock or injury to persons:

1. This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Do not operate this appliance under a blanket or pillow. Excessive heating may occur and cause fire, electric shock or injury to persons.
3. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
4. This massager contains a motor. To avoid risk of injury, it must not be used within six inches of an implantable cardiac device.
5. Use this appliance only for its intended use, as described in this manual. Do not use attachments not recommended by the manufacturer.
6. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return to Wahl for examination and repair.
7. Do not carry this appliance by supply cord or use cord as a handle.
8. Keep the cord away from heated surfaces.
9. Never drop or insert any object into any opening on this appliance.
10. Do not use outdoors.
11. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
12. To disconnect, turn all controls to the off position, then remove plug from outlet.
13. Do not use massager in close proximity to loose clothing or jewelry.
14. Keep long hair away from massager while in use.
15. In case of pregnancy, diabetes or illness, consult your physician before using this massager.
16. Do not use massager on:
 1. Sensitive skin, swollen or inflamed areas, in the presence of poor circulation, on insensitive skin, where skin eruptions are present.
 2. The front of the neck.
 3. Genital areas.
 4. Sleeping or unconscious persons.
17. Do not use this massager for more than 20 minutes in any one area.
18. Before you unplug from outlet, hold the power button for seconds to turn OFF and press the heat button to turn OFF, then remove plug from outlet.

19. Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or disabled person or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
20. When storing this product, do not crush and avoid sharp folds.
21. Never operate the appliance with the air openings blocked. Keep the air openings (if any) free from lint, hair or the like.
22. This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Keep children under the age of 13 away from this machine.

Some of our massagers have heat attachments and gel packs.

For a unit that has heat:

23. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

For a unit that has a gel pack:

24. Never apply the frozen gel pack directly to skin as direct application may cause frostbite. Never apply the heated gel pack directly to skin as direct application may cause burns.

Warning

This content is compiled from multiple sources and is provided for reference purposes only. It may not be complete or fully applicable to all situations. If you are unable to resolve your issue, please contact the product manufacturer or an authorized service provider for official support.