

About the Renpho Mini Massage Gun

1. What's in the box

- Renpho Mini Massage Gun x 1
- Massager Head x 5 Type C
- Charging Cable x 1 Rubber
- Sleeve x 2 User Manual x 1
- Warranty Card x 1

2. Product Introduction



TIP: None of the massage heads are to be used for your head. DO NOT use on your head

<p>Ball head</p> <p>Ball head: For large muscle groups like quads, glutes, extreme muscles, and joint soreness.</p>	<p>U-shaped head</p> <p>U-shaped head: for spine, neck, and shoulder for sensitive muscles and soft tissues.</p>
<p>Air cushion head</p> <p>Air cushion head: for sensitive body meat relaxation.</p>	<p>Flat head</p> <p>Flat head: for all parts of the body</p>
<p>Bullet head</p> <p>Bullet head: for joints, deep tissues, trigger points and small muscle areas like feet, and wrist.</p>	<p>Rubber Sleeve</p> <p>Each massage head is pre-installed with a rubber sleeve, if it is damaged, replace it with a spare one.</p>

Using Mini Massage Gun

1. Important tips before use

1.1 Charging

- Please charge the massage gun before use.
- To charge the massage gun, plug the type C cable into an adapter and into the unit's charging interface.
- The indicator light on the massage gun will BLINK BLUE during the charging process. The indicator light will be SOLID BLUE when charging is complete.
- For safety considerations, the massage gun will NOT work while charging. Additionally, the massage gun will automatically shut off after 10 minutes of continuous use. If you want longer massage, please turn it on then.



1.2 Battery Life

- This Massage gun has a built-in rechargeable battery. A full charge takes about 3 hours. Total battery life is up to 150 minutes after a full charge. Please fully charge it before first use.

- Battery life can vary by application and/or operating conditions. If the product is not in use for a long time, please keep it powered off and charge it once a month so that the battery will last longer.
- Please keep the massage gun away from heat and avoid using it in places where the surrounding temperature is relatively hot.

2. Preparation


- Please fully charge the massage gun for 3 hours before use.
- Select the appropriate massage head to connect with the host.



3. Operation Instruction


1. Connect the device with the massage head

Insert the massage head completely into the host when using, and plug out when stop using.

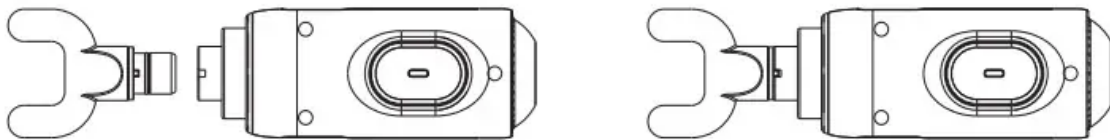Note: The installation of the U-shaped head should be positioned. When the buckle is inserted into the slot, it means the installation is in place. If the massage gun's head is not installed in place, there will be a risk of falling off while use.



2. Long press the  button on the panel for 2 seconds to power on and start the first speed, short press the  button on the panel, the speed will switch between 1-5 level.

Long press the  button for 2 seconds to power it off.

* After 10 minutes of operation, the device automatically powers off.


3. Use the massage gun on any area of your body to relieve stiffness, soreness, and pain.



Shoulder Massage




Back Massage

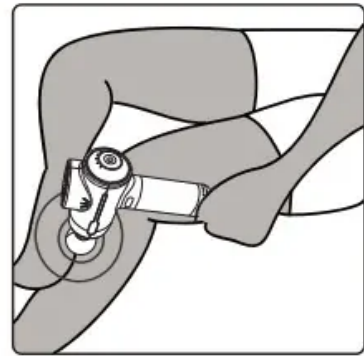

Foot Massage



Arm Massage

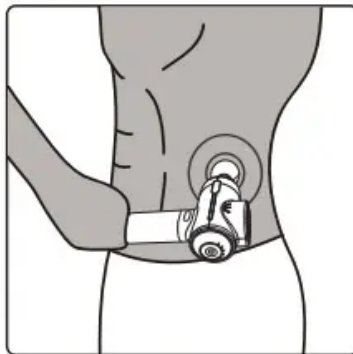

Waist Massage



Leg Massage

Product Care and Maintenance

- Place the massage gun in a safe, cool, and dry place when not in use. •Always unplug the massage gun prior to cleaning.
- After use, wipe massage gun with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
- DO NOT submerge the massage gun in water or any other liquid.
- Keep away from all solvents and harsh detergents.
- DO NOT attempt to repair this massage gun. There are no user-service- able parts.
- DO NOT dispose the massage gun in a land fill. It contains a Li-ion battery and other electrical components, which is considered hazardous E-waste.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	Power Is off.	Long press the Power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically	The massage gun's overheat protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before powering it back on.
	Massage gun shuts off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is to inserted too tight.	Plug out the massage head a little, then power it on.

Specification

Product Name: Renpho Mini Massage Gun Model: RP-GM171 Rated Voltage: 7.4V=

Rated Power: 25W

Input: 5V~2A

Battery Capacity: 2500mAh

Amplitude: 10mm

RPM: Max 3200r/min

Auto Shut-off Timer: 10 minutes

Product Dimensions: 7.3*5.7*2.3 inches

Product Weight: 1.5lbs

Color: Black / Blue / Silver

Important Safety Instructions

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Any individual who may be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor's care should consult with a medical professional before using a massaging device designed for home use.
- DO NOT use the massage gun on:
 - Sensitive skin, swollen or inflamed areas, in the area of poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain.
 - Your head.
 - Varicose veins.
 - The throat area.
 - Genital area (DO NOT use this for sexual purposes. This is a powerful massage gun and may cause severe bodily harm).
 - An unconscious or sleeping person.
- DO NOT use this product for more than 30 minutes at a time.
- If you feel unwell or in any unusual discomfort, please stop using the product immediately.
- Extensive use could lead to the product's excessive heating and shorter product life.
- DO NOT use this product before bed. The massage gun has a stimulating effect and can delay sleep.
- Please keep the massage gun away from heat and avoid using it in places where the surrounding temperature is relatively hot.
- DO NOT use the massage gun with wet hands.
- Never cover the product when it is in operation.
- This unit should not be used by children or invalids without adult supervision.
- It is prohibited for non-professional personnel to repair the massage gun or replace parts.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if

they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- The use of a third-party massage head may cause damage to the massage gun.
- Use heated surfaces carefully to avoid serious burns.
- DO NOT drop, apply excessive force or heavy weight over the product to avoid malfunction and/or damages.
- DO NOT keep product near a high magnetic field or where it can be affected by a magnetic field, for it can cause batteries to discharge.
- WARNING Chemical Exposure: If the battery pack is leaking, DO NOT allow the liquid to come in contact with skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice. Place the battery in a sealed plastic bag and dispose of safely according to local environmental regulation.
- DO NOT disassemble or drill into the product and/or batteries. DO NOT take apart, modify or change the parts. This may cause product malfunction, electric shock, or fire. And the product will no longer be covered under manufacturer warranty.
- DO NOT throw the battery into fire or heat it.
- The used battery should be placed in a sealed plastic bag and disposed of safely according to local environmental regulations.
- Battery charging time and usage time may decrease the more the battery is used.

Warning

This content is compiled from multiple sources and is provided for reference purposes only. It may not be complete or fully applicable to all situations. If you are unable to resolve your issue, please contact the product manufacturer or an authorized service provider for official support.