



coro user guide


What's inside?

- Coro nipple shield
- Coro charging and storage pod
- USB cable
- User guide
- Cleaning syringe

Size: Coro nipple shield is designed with the same dimensions of a standard commercial nipple shield of medium size.

Please review our size guide at www.coroflo.com.

Main compartment

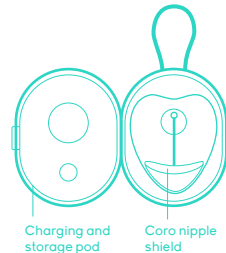

Charging and storage pod


Coro nipple shield



Cleaning syringe

What's inside?


Front compartment

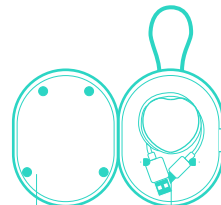

Charging and storage pod

USB cable

Charging Coro

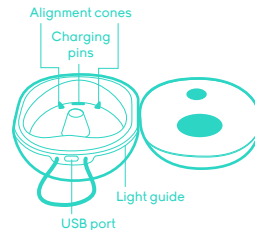

Alignment cones

Charging pins

Light guide

USB port

Charging Coro

Step one

Insert Coro into the Coro charging and storage pod. Use the alignment cones to guide the shield onto the charging pins.

Step two

Plug the USB cable into the USB port and connect the other end of the USB cable to a USB port or to a wall charger. If your phone is compatible, you can charge Coro from your phone (USB cable not provided). The average charging time is approx 30 min.

Step three

The light guide will indicate the charging status of the shield with a dim light and it will switch off when fully charged.

You can store the USB cable in the front compartment of the Coro charging and storage pod.

The battery life for the fully charged Coro shield is approximately 3+ hours of continuous use. We recommend fully charging Coro before each feed.

In order to protect it and preserve its materials, the Coro shield should be stored out of direct sunlight. It is recommended to store the Coro shield clean and dry in its charging and storage pod.




The Coro explainer video is available at www.coroflo.com.

Download the App



The Coro app is compatible with iOS and Android. Go to www.coroflo.com and download it for free from the App or Google Play stores.

It is recommended that you register in order to access all your breastfeeding data. Coroflo's statement about data protection is available at www.coroflo.com.

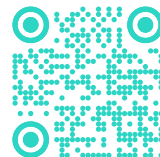

Important safety information



Read all the instructions before using this product.

CAUTION! Coro is not a medical device and should not be used to diagnose, treat or prevent any medical condition. Always remember to contact your lactation consultant or GP may you experience any issues or pain while breastfeeding. Guidance on the use of nipple-shields can be sought from a qualified lactation consultant or GP. Coroflo can also provide sources of information on request. Additional information is available at www.coroflo.com.

Cleaning Coro

IMPORTANT: Ensure Coro has been cleaned fully before using.

SYRINGE

Before use, separate all syringe parts and wash in warm soapy water. Leave to dry on a clean surface.

Cleaning Coro

CORO



01

Prepare two small bowls: one with clean warm water mixed with a few drops of dishwashing liquid, one with clean drinking quality water.



02


Fill the syringe fully with water and dishwashing liquid.



03

Place the syringe fully into the Coro teat and carefully inject all the water solution. Be sure to push down the syringe slowly for a few seconds to ensure the Coro teat is cleaned fully. Rinse the soapy syringe.

Cleaning Coro

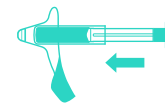

04

Repeat the same process with the clean water until no soapy water remains and then rinse it.



05

Get rid of any remaining water by pushing air in with the syringe. Repeat as necessary.



06

Leave all parts on a clean surface to air dry.

Cleaning Coro

Coro must be rinsed thoroughly to avoid dry up of breast milk and prevent bacterial growth. **Under no circumstances should any other implement or tool be inserted into the teat as this will damage the sensor; this also excludes the use of any type of cleaning brush (soft-bristled or otherwise).** If so desired, sterilisation liquid or tablets such as Milton may be used on the Coro shield, but only in cold, drinking-quality water.



CAUTION! In order to avoid infections, always inspect the Coro shield and the inside of the Coro charging and storage pod before use to check for dirt. Do not use the Coro shield if there is visible dirt or contamination on it.



The Coro Shield and its charging and storage pod are intended for use by a single user only and should not be shared between babies for hygiene reasons.



CAUTION! The expected lifetime of the silicone parts is 6 months of continuous use, this can be longer with intermittent use but if there are any signs of degradation or wear and tear, Coro needs to be replaced.

Start a feeding session



- 01 If not already on, enable the Bluetooth on your phone.



- 02 Remove the Coro shield from the charging and storage pod (this activates Bluetooth on the device). If any connection issue arises, return the Coro shield to the charging and storage pod and retry.



- 03 Apply the nipple-shield section of the device to the breast in accordance with professional advice.

Start a feeding session



- 04 Open the Coro app, select either left or right breast.



- 05 Commence feeding as normal. If you are switching breasts, simply press the other breast button.



- 06 When your feeding session is finished, press the Stop button.


Start a feeding session



- 07 Remove the Coro shield from your feeding breast.



- 08 After each feed, wash the Coro shield according to our cleaning Coro guidelines.

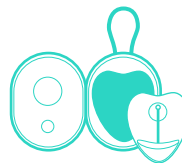

- 09 When dry, return the Coro nipple shield into its charging and storage pod.

Avoiding injury



WARNING! The battery is not user-replaceable and no attempt should be made to open or modify the device in case of a suspected battery issue. Please return the unit to Coroflo in this case.



WARNING! In the event the battery is removed **KEEP OUT OF REACH OF CHILDREN**. Swallowing may lead to serious injury or death in as little as 2 hours due to chemical burns and potential perforation of the oesophagus.



WARNING! The Coro shield and its charging and storage pod must be inspected prior to each use. Do not use the device if there is any visible sign of damage since this may cause injuries. If you notice any cracking or clouding of the silicone area discontinue use and return the unit to Coroflo.



WARNING! Do not attempt to open or repair Coro yourself, return the unit to Coroflo instead. No modifications to the device are permitted.



WARNING! Charge Coro in accordance with the instructions provided. Only use the USB cable provided in the box or a USB cable, smartphone, computer, powered hub or power supply that is certified by a recognized testing laboratory.



WARNING! If you have any trouble using the device or if you experience any malfunction immediately discontinue use and return the unit to Coroflo.

Avoiding injury



WARNING! If the light guide shows a flashing light when the USB cable is connected to the charging pod, this may be an indication of device malfunction or excessive temperature. If the situation persists in a cool environment, disconnect the USB cable immediately, discontinue use, and return the unit to Coroflo.



WARNING! Do not return the Coro shield to the pod if wet or damp.



CAUTION! This product is not a toy. Do not allow children to play with the Coro shield or the Coro charging and storage pod. Avoiding damage to the device

Avoiding damage to the device



WARNING! Do not use microwaves, heat sterilizers, dishwashers, washing machines, boiling water, bleach/cleaning agents or use any cleaning instruments on the shield.



WARNING! Never use water or damp cloths to clean the charging pins on the Coro charging and storage pod as water can damage the electronics inside.



Store the device away from temperature extremes and out of direct sunlight.

If the device stops responding and will not connect, it is possible that the internal battery has discharged. Return the device to the charging pod and allow it to charge via USB before attempting reconnection.

Regulatory information

Coro has passed the consumer product safety requirements and attained the CE & UKCA Marks, declaring conformity to the requirements of applicable Directives.

It has passed independent food-contact testing.

It has passed electromagnetic safety testing and it is equipped with a Bluetooth Low Energy radio transceiver.

All the materials used for the production of the Coro shield and its charging and storage pod are food-contact safe, latex free, and BPA free.

Certification and test details are available from Coroflo.



Disposal and recycling information

The Coro shield and its charging and storage pod are made of different materials and contain electronics which cannot be removed prior to disposal.

Do not discard the device into the household waste; instead follow the disposal guidelines of your local collection system for electronic equipment (including batteries).



Warranty

Coro and its charging and storage pod are backed by a Limited Manufacturer's Warranty against any defects in materials and/or workmanship for two (2) years from the date of original purchase, provided that the product has been used in accordance with the information provided in the product documentation and that all operating and maintenance instructions are strictly respected.

No person or entity is authorised to make any other warranty on behalf of Coroflo.

Exclusions and limitations

This Limited Warranty only applies to this product and does not cover other hardware not supplied in the Coro box. This Limited Warranty also does not apply in the case of:

1. Cosmetic deterioration or damage caused by normal wear and tear included aging of materials [such as discolouration or stretching].
2. Damage caused by failure to follow the instructions for use provided by Coroflo;
3. Use of improper voltage or power supply for charging;
4. Use of the device with, or connecting it to, any product, accessory, software or service not manufactured or supplied by Coroflo;

Warranty

5. Accidental, physical or natural destruction, misuse, neglect or other external causes;

6. Damage caused by service by anyone except a Coroflo representative;

7. Any attempt being made to open or modify Coro or its charging and storage pod [such as to alter its functionality or capability];

8. The serial number on the device being removed, erased, defaced, altered, or rendered illegible in any way;

9. Refusal to return the unit to Coroflo for repair and investigation.

Please note that, for reasons of hygiene, Coro is meant to be used with one baby only; use by more than one person will invalidate this warranty.



Coroflo Ltd.,
DCU Alpha, Old Finglas Road,
Glasnevin, Dublin 11,
D11 KXN4, Ireland
<https://www.coroflo.com>
info@coroflo.com



Warranty

In case a defect is found

Coroflo must be notified of any defective device. If a defect is found in your Coro you must contact Coroflo within the product Limited Warranty period. When a valid claim is received, Coroflo will reserve the right to decide whether to replace or to repair the defective product at no charge [shipping and handling charges may apply]. The duration of this Limited Warranty will not be extended if a repaired/replaced product is received. This Limited Warranty is valid only when Coro is sold by Coroflo or through its authorised resellers.

To obtain warranty service, you must ship the device in its original packaging, or in other packaging providing adequate protection, to the address specified by Coroflo. All claims must be supported by evidence of the defect found and with proof of date of purchase. More information on Coro warranty and on how to receive warranty service is available at <https://www.coroflo.com>. Coroflo will not be responsible for any damage which may arise after replacing your device.

Disclaimer

Coroflo does not assume any responsibility or liability associated with or for the use of Coro. Coroflo reserves the right to revise and/or update this document without notification to the user.

Printing notes

1. Document is 97mm x 387mm
2. Artwork is Pantone 3255C (Turquoise)
3. The background colour of the document is white
4. Bleed is 5mm all round
5. Pink border is not to be printed, these are the keylines
6. The document is to be printed as a double sided print