
OMRON

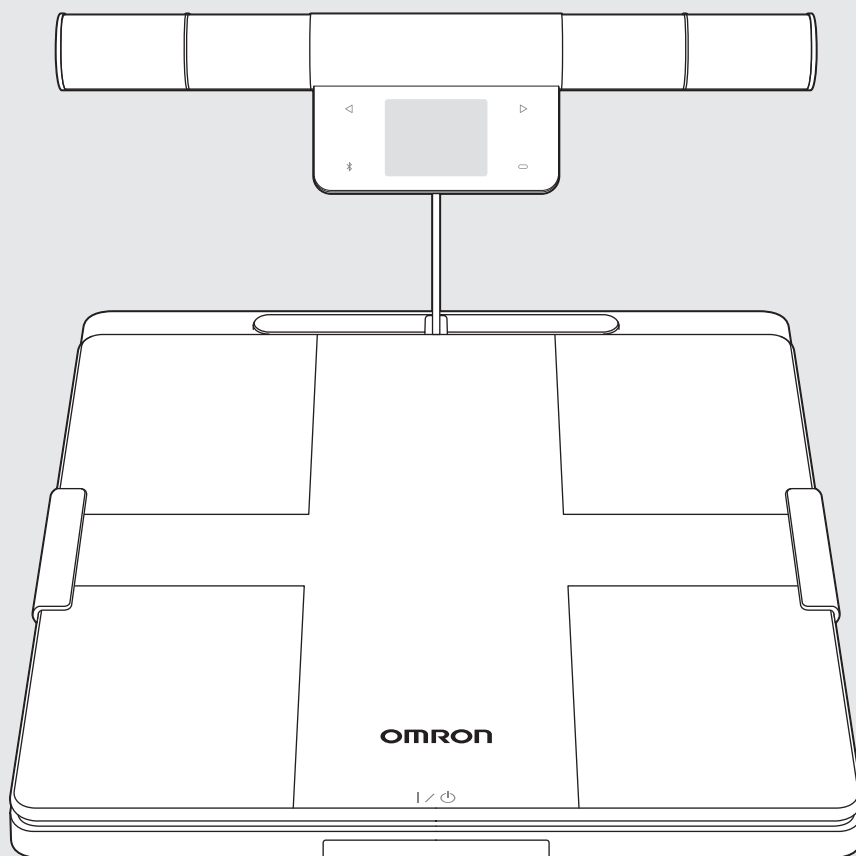VIVA 360+ Smart Scale
Model BCM-750

INSTRUCTION MANUAL

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MANUEL D'INSTRUCTIONS

Table of Contents

Introduction	EN2
Safety Instructions	EN2
Intended Use/Indications for Use	EN2
Receiving and Inspection.....	EN2
Symbols Glossary.....	EN2
Important Safety Information	EN2
1. How To Use	EN4
1.1 Package Contents.....	EN4
1.2 Know Your Device	EN4
1.3 Installing Batteries.....	EN5
1.4 Initial Settings.....	EN6
1.5 For Accurate Measurement	EN8
1.6 Taking a Measurement	EN10
1.7 Automatic Transfer of Readings.....	EN15
1.8 Checking Your Readings Using the Compatible OMRON App	EN16
1.9 Measuring Weight Only.....	EN16
1.10 Measuring without Recording (Guest Mode)	EN17
1.11 Selecting the Personal Number Manually Before Measurement	EN18
2. Settings	EN19
2.1 User Settings/Display Settings/Deleting Data	EN19
2.2 Replacing the Batteries.....	EN21
2.3 Restoring the Unit to Factory Settings	EN21
2.4 Performing Initial Settings on the Unit.....	EN22
2.5 Registering Personal Data on the Unit	EN23
3. Error Messages and Troubleshooting	EN24
3.1 Error Messages	EN24
3.2 Troubleshooting.....	EN25
4. Maintenance	EN27
4.1 Maintenance.....	EN27
4.2 Storage.....	EN27
5. Specifications	EN28
5.1 Specifications.....	EN28
5.2 Important Information regarding Electromagnetic Compatibility (EMC).....	EN29
6. Limited Warranty	EN30
7. FCC/ISED Statement	EN30
8. Information on Body Composition	EN31
8.1 Information on Body Composition.....	EN31
8.2 Measurement Results Interpretation Diagrams	EN33

Introduction

Thank you for purchasing the OMRON VIVA 360+ Smart Scale.

This device is intended for measuring and displaying the following body composition parameters.

Safety Instructions

This instruction manual provides you with important information about the OMRON VIVA 360+ Smart Scale. To ensure the safe and proper use of this device, READ and UNDERSTAND all of these instructions.

If you do not understand these instructions or have any questions, contact 1-800-634-4350 before attempting to use this device.

Intended Use/Indications for Use

The OMRON VIVA 360+ Smart Scale is intended to be used to calculate and display the estimated value of body fat percentage, subcutaneous fat percentage, skeletal muscle percentage, resting metabolism, body water percentage, bone mass, fat free mass and visceral fat levels using the BI (Bioelectrical Impedance) method and indicates the BMI (Body Mass Index).

The OMRON VIVA 360+ Smart Scale is intended to be used by healthy individuals in the age range of 18 to 80 years for the subcutaneous fat percentage and visceral fat level. The age range for body fat percentage, skeletal muscle percentage, resting metabolism, body water percentage, bone mass and fat free mass are intended to be used by healthy individuals between the age of 6 to 80 years. The device is intended for multiple-patient, multiple-use.

Receiving and Inspection

Remove this device and other components from the packaging and inspect for damage. If this device or any other components are damaged, DO NOT USE and contact 1-800-634-4350.

Symbols Glossary

For symbol information, visit: OmronHealthcare.com/symbols-glossary

Important Safety Information

Read the Important Safety Information in this instruction manual before using this device.

Follow this instruction manual thoroughly for your safety.

Warning

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

- Never use this unit in combination with medical electronic devices such as:

- ① Medical electronic implants such as pacemakers.
- ② Electronic life support systems such as an artificial heart/lung.
- ③ Portable electronic medical devices such as an electrocardiograph.

This unit could cause these devices to malfunction, posing a considerable health risk to users of these devices.

- DO NOT use the unit on slippery surfaces, such as a wet floor.
- DO NOT tip, drop, or subject the unit to intense impact.
- DO NOT jump on the unit. Doing so may cause you to fall or slip, resulting in serious injury.
- DO NOT use the unit when your body and/or feet are wet, such as after taking a bath or shower.
- DO NOT step on the edges, display area, or display unit holder of this unit because you may lose your balance and fall, resulting in serious injury.
- Persons with disabilities or persons who are physically frail should be assisted by another person when using this unit to prevent falling when stepping on and off the unit.

Data Transmission

- This product emits radio frequencies (RF) in the 2.4 GHz band. DO NOT use this product in locations where RF is restricted, such as on an aircraft or in hospitals.

**Caution**

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient, or cause damage to the equipment or other property.

- Consult with your physician or healthcare provider before beginning a weight reduction or exercise program.
- Keep the unit out of the reach of infants, toddlers and children.
- DO NOT attempt to disassemble, modify or repair the unit.
- DO NOT place the unit on a cushioned floor surface such as a carpet or mat because the measurement may not be accurate.
- Step onto the unit with bare feet when measuring.
- This unit is not for professional use. DO NOT use this unit in hospitals or other medical facilities.
- Remove the display unit from main unit before stepping on the main unit.
- DO NOT use this unit with other medical electrical (ME) equipment simultaneously. This may result in incorrect operation of the devices and/or cause an inaccurate reading.
- During measurement, make sure that no mobile device or any other electrical device that emits electromagnetic fields is within 12 inches (30 cm) of this device. This may result in incorrect operation of the device.
- OMRON recommends waiting for approximately 2 hours for the unit to warm up or cool down when the unit is used in an environment within the temperature specified as operating conditions after it is stored either at the maximum or at the minimum storage temperature. For additional information of operating and storage/transport temperature, refer to section 5.

Battery Usage

- Replace worn batteries with new ones immediately.
- ONLY use 4 "AA" alkaline batteries with the unit. DO NOT use other types of batteries. DO NOT use new and used batteries together. Replace all 4 batteries at the same time.
- DO NOT use different types of batteries together.
- Remove the batteries if the unit will not be used for 3 months or more.
- DO NOT use batteries after their expiration date.

Data Transmission

- DO NOT replace the battery while your measurement result is being transferred to your smart device. This may result in the incorrect operation of the unit and failure to transfer your measurement result.
- DO NOT place integrated circuit cards, magnets, metal objects, or other devices that emit electromagnetic fields near the unit while your measurement result is being transferred to your smart device. This may result in the incorrect operation of the unit and failure to transfer your measurement result.

**About a wireless communication interference**

The **Bluetooth®** option in the product is used to connect to dedicated apps on mobile devices to synchronize date/time data from mobile device to the product, and to synchronize measurement data from the product to mobile device. Further handling of the data on the mobile device is up to the user's discretion. This product operates in an unlicensed ISM band at 2.4 GHz where any third party can intercept the radio waves, willfully or accidentally, for any unknown purpose. In the event this product is used near other wireless devices such as microwave and wireless LAN, which operate on the same frequency band as this product, there is a possibility that interference may occur. If interference occurs, stop the operation of the other devices or relocate this product away from other wireless devices before attempting to use it.


1. How To Use

1.1 Package Contents

Main Unit, 4 AA Alkaline Batteries, Instruction Manual

1.2 Know Your Device


Main Unit

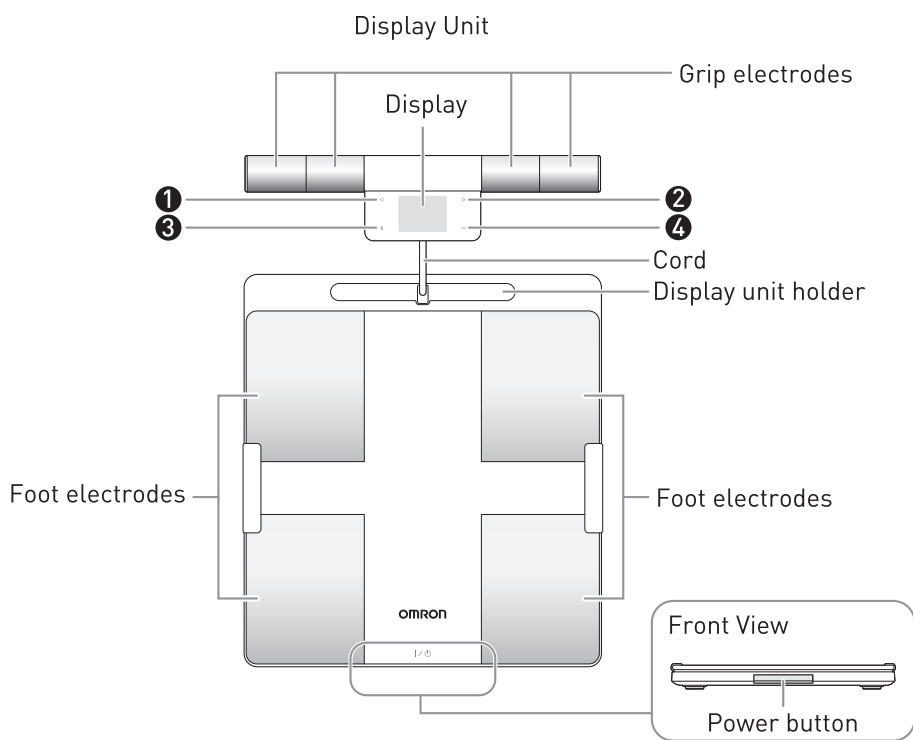Icon	Touch Button Name	Description
① 	Back touch button	Touch this button to confirm the readings, select a personal number, and set your personal data.
② 	Advance touch button	
③ 	Bluetooth touch button	Touch this button to communicate with your smartphone via the compatible OMRON app. When pairing manually this product with your smartphone, touch and hold this button for more than 2 seconds.
④ 	OK touch button	Touch this button to confirm the value.

When you operate a touch button, the unit makes a “beeping” sound.

Display

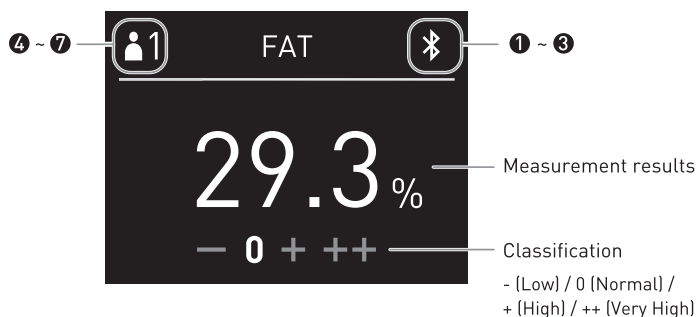
Icon	Name	Description
①		Bluetooth symbol Appears when pairing the device or transferring readings.
②		SYNC symbol Flashes when there are many readings that have not been transferred.
③		Battery symbol When the battery is low, the symbol flashes. When the battery is depleted and needs to be replaced, and "A1" appears. When replacing batteries, replace all four batteries at the same time with new AA alkaline batteries of the same brand. This will not be displayed if there is sufficient battery.
④		Personal number symbol Appears when selecting the personal number 1~4.
⑤		Guest mode symbol Appears when measuring without recording the readings. Select "Guest" for the personal number.
⑥		Weight symbol Appears when measuring only weight.
⑦		Setting symbol Appears when performing initial settings on the unit or display settings.

1.3 Installing Batteries

This product uses AA alkaline batteries. Do not use manganese batteries or rechargeable batteries as they may significantly reduce the number of measurements or cause malfunction.


1. Open the battery cover on the back of the unit.

- A hard object like a pen could be used.

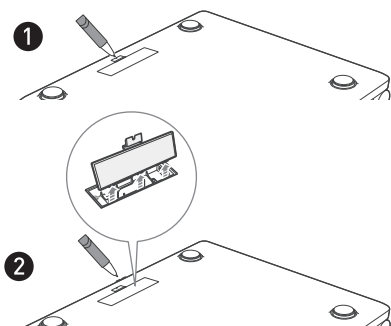2. Install the 4 AA alkaline batteries in the correct polarity as marked inside the battery compartment.

- Press the battery down to make sure it is installed properly.

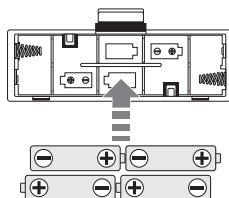3. Close the battery cover.

- Close it until you hear a click.

1.4 Initial Settings

If you prefer not to use the app, settings can be configured directly on the unit. Refer to sub-sections 2.4 Performing Initial Settings on the Unit and 2.5 Registering Personal Data on the Unit for details.

Installing the compatible OMRON app

This unit can transfer your readings to your smartphone, and you can confirm your weight and body composition changes simply by installing the compatible OMRON app on your smartphone. The compatible OMRON app allows you to check various metrics, including weight, body fat percentage, skeletal muscle, BMI, and many more. After installation, register this product and your personal data (date of birth, gender, height and personal number) in the compatible OMRON app.

Note

- Review the list of compatible smartphones at: <https://omronhealthcare.com/compatibility>
Any devices and operating systems not listed on our website are not supported.

1. Go to “Settings” on your smartphone and enable the Bluetooth function.
2. Scan the code below with your smart phone camera or visit OmronHealthcare.com/app to find the compatible OMRON app for this product.

3. Download and install the compatible OMRON app onto your smartphone.

If you already have the compatible OMRON app and have created your account, open the app and add your new unit.

4. Open the compatible OMRON app and register this product and your personal data.
 - Follow the instructions in the compatible OMRON app to pair your smartphone and register.
 - If you are already using the compatible OMRON app, register this product using the following steps.
 - (1) Tap the “Devices” icon at the bottom of the compatible OMRON app screen
 - (2) Tap “Add Devices” and follow the on-screen instructions

5. Confirm that your unit is connected successfully.

When your unit is connected successfully to your smartphone, the pairing complete display will appear.


Register your personal data (current body composition)


The registered personal data will act as a standard for estimating the measured individual (personal number). The readings will be automatically transferred to your smartphone. Keep the compatible OMRON app open until you have completed registering your personal data.

1. Press the power button to turn on the power.

- The weight measurement display will appear.

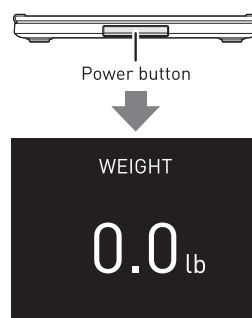
2. Touch  and select the personal number registered in the compatible OMRON app.

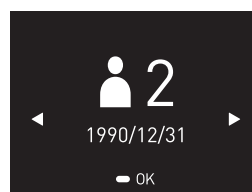3. Touch .


4. After the “0.0 lb” is displayed, take out the display unit and hold the grip electrodes with the upper part of your palms as shown in the image.

- Do not take out the display unit until “0.0 lb” appears on the display.
Otherwise, the weight of the display is added to your body weight which causes an incorrect result.

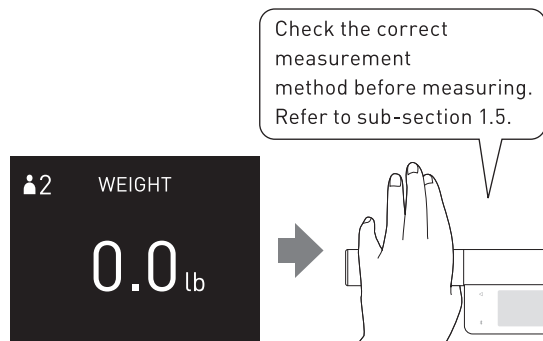
5. Step on the unit to start the measurement.

- Place your bare feet on the foot electrodes with your weight evenly distributed.

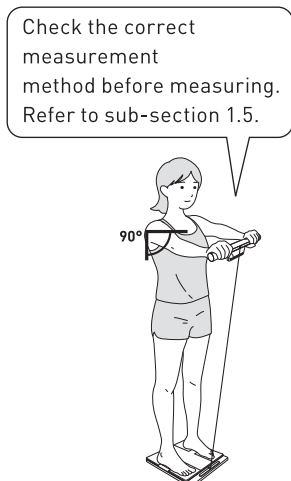- Once your weight is shown on the display, the unit will then start to measure your body composition.

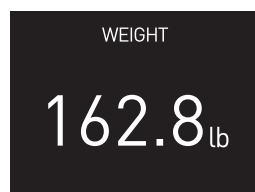6. Extend your arms straight at a 90° angle to your body.

7. After the measurement is complete, the estimated personal number and body composition are displayed. At this point you may step off the unit.

- The unit makes a “beeping” sound.

8. Press the power button to turn off the power.


- Even if you forget to turn off the power, it will automatically turn off after approx. 16 seconds.
- To change or delete your personal data, refer to sub-section 2.1.
- To cancel the setting, press the power button to turn off the power.
- One smartphone can manage the readings of one person. If you are managing the readings of multiple people, prepare one smartphone for each personal number.
- To use automatic transfer, location settings on your smartphone must be enabled.
Even if the location is enabled, your personal information will not be leaked.

1.5 For Accurate Measurement


To take accurate measurements, follow the items below.

Check the Correct Measurement Method


Press your palms firmly on the grip electrodes.

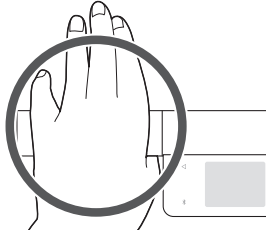- Hold the grip electrodes with the upper part of your palms.

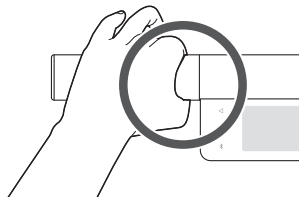
- Hold the inner grip electrodes firmly with your thumb and index finger.

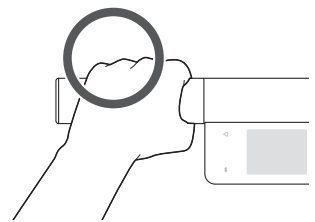
- Hold the outer grip electrodes with your ring finger and small finger.

Correct Posture for Measurement

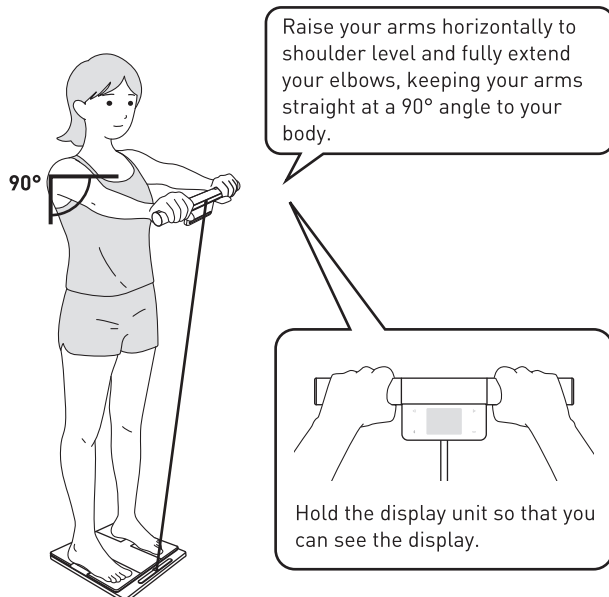
Step on the main unit with bare feet.

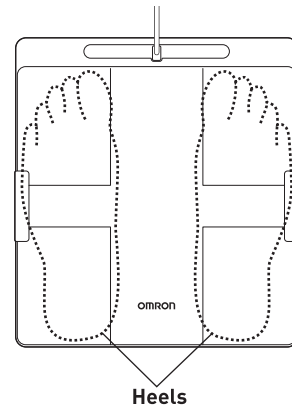- Make sure each of your heels is positioned as shown. Stand with your weight evenly distributed on the measurement platform.


Postures to Avoid During Measurement

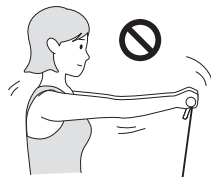- Movement during measurement


- Arms bent

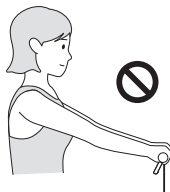- Arms too low or high

- Knees bent

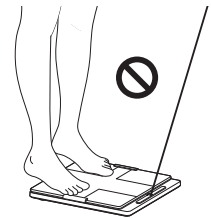- Standing on edge of unit

Recommended Measurement Times

The following times are suitable for measurement:

- After waking up and before breakfast
- Before lunch and about 2 or more hours after breakfast
- Before dinner and about 2 or more hours after lunch
- Before going to bed and about 2 or more hours after dinner or bath




Install the Unit

The unit is regularly calibrated so that the weight is “0 lb” when not in use. When setting up the unit, place it on a hard, flat, and even surface.

When “0 lb Calibration” Is Required


If the unit is placed in conditions as the figures below, “0 lb calibration” is required because the weight displayed may differ from the actual weight.

- The unit is leaning against a wall or other objects
- The unit is placed upside down
- The unit is placed on an object
- An object is placed on the unit

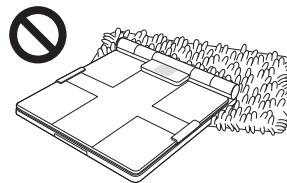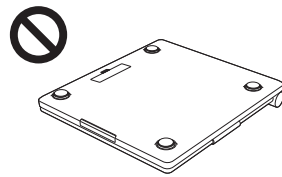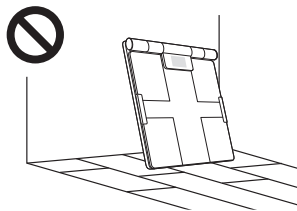Implement the “0 lb Calibration”

1. Place the unit on a hard, flat, and even surface.

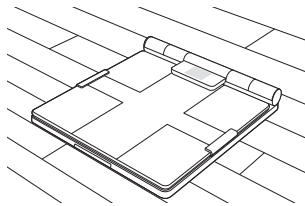2. Press the power button to turn on the power.

- The weight measurement display will appear.

3. With the unit placed on the floor, press the power button.


- Power is turned off.
Even if you forget to turn off the power, it will automatically turn off after approx. 1 minute.
- Turn off the power, wait for at least 5 seconds, then take a measurement.

1.6 Taking a Measurement

Measurements should be taken on a hard, level surface. The readings are automatically transferred to your smartphone. You can check indicators on the compatible OMRON app. Refer to sub-section 1.8.

If you do not transfer the data to your smartphone, it will be saved in the unit's memory. The unit's memory can store up to 30 readings for each personal number. Readings in the unit's memory can be transferred to your smartphone later.

1. Press the power button to turn on the power.

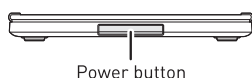
2. After the "0.0 lb" is displayed, take out the display unit and hold the grip electrodes with the upper part of your palms as shown in the image.

- Do not take out the display unit until "0.0 lb" appears on the display.

Otherwise, the weight of the display is added to your body weight which causes an incorrect result.


Check the correct measurement method before measuring. Refer to sub-section 1.5.

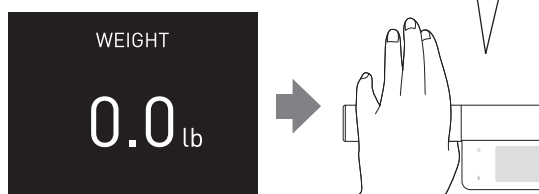
3. Step on the unit to start the measurement.

- Place your bare feet on the foot electrodes with your weight evenly distributed.

Check the correct measurement method before measuring. Refer to sub-section 1.5.

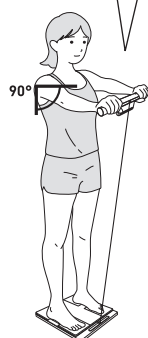- Once your weight is shown on the display, the unit will then start to measure your body composition.

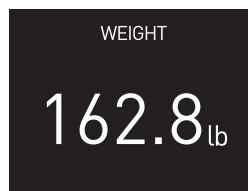
4. Raise your arms horizontally to shoulder level and keep your arms straight at a 90° angle to your body.

5. After the measurement is complete, the estimated personal number is displayed. At this point you may step off the unit.


- The unit makes a "beeping" sound.


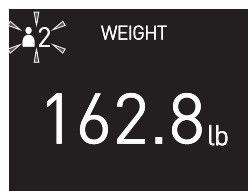
6. Check whether the estimated personal number is correct.

- Check whether your personal number, which you set when registering your personal data, is displayed. This product automatically estimates your personal number based on your weight and body composition using your personal data. However, if someone with a similar body shape or weight takes a measurement, the correct personal number may not be displayed.
- If your personal number is displayed, the power will automatically turn off after approx. 16 seconds without any operation, and the readings will be transferred to your smartphone. Refer to sub-section 1.7.
During transfer, the Bluetooth symbol () on the display will light up.
- If your personal number is different, touch  to select your personal number, then touch . The readings are transferred to your smartphone.

Modification example:

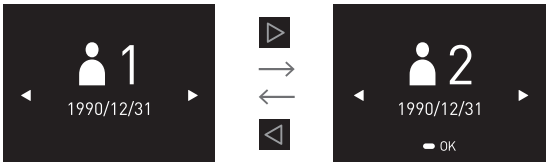- If your personal number is displayed as shown below, it means that your personal number cannot be estimated from the readings.

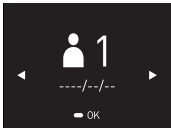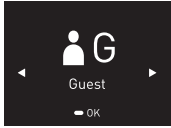
Display	Cause	Action
	Personal data is not registered	Register your personal data using the compatible OMRON app. Refer to sub-section 1.4. Or, touch  and register your personal data. Refer to sub-section 2.5.
	The readings do not match any registered personal data	There may be significant differences between the registered personal data and the readings. Touch  to select the registered personal number, then touch  .

7. Check the measurement results.

- The display on the unit will repeatedly show the measurement results in sequence. (The display will change automatically.)

Display example:

- 4 indicators

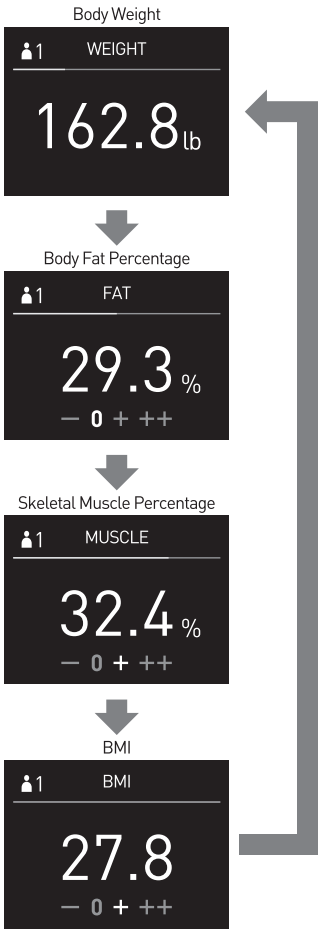The unit will show the above metrics. Other metrics are shown on the app.

1.6 Taking a Measurement

• 18 indicators

- You can set all indicators to be displayed in sub-section 2.1, "Display settings".
- The unit will automatically turn off if no buttons are touched: after 16 seconds for the 4 indicator display, and after 72 seconds for the 18 indicator display.
- Once the  or  button is touched, the display will change, and the unit will automatically turn off after approx. 3 minutes.
- Pregnant women may not be able to measure their body composition properly.
- When the number of readings stored in the unit's memory exceeds 30, the oldest readings will be deleted.
- In the following cases, the readings are saved in the unit's memory.
 - The smartphone on which the compatible OMRON app is installed is turned off
 - Bluetooth communication is disabled
 - The distance between the unit and the smartphone is 16.4 ft or more
- If you enter a different personal number, the correct body composition value will not be displayed.
If you have entered a different personal number, touch , then touch  to select your personal number again.

Calculate Your FIT Index in Athlete Mode

Athlete Mode is used to measure body composition in people who engage in high-frequency, high-intensity, and long-duration training as part of their daily routine, and who have greater muscle mass relative to their BMI.

The estimation formula of body composition values differs between Athlete Mode and Normal Mode.

If you are unsure whether to select Athlete Mode, refer to your FIT Index, calculated based on the following three factors: frequency, intensity, and duration of exercise. It is recommended to use Athlete Mode if your FIT Index score is 60 or higher.

$$\text{FIT Index} = \text{Frequency} \times \text{Intensity} \times \text{Time}$$

Based on your FIT Index use the following mode when setting your personal data.

NORMAL = FIT Index < 60

ATHLETE = FIT Index ≥ 60

EXAMPLE: Running 5 times a week for 30 minutes.

$$\begin{array}{ccccccc} 75 & = & 5 & \times & 5 & \times & 3 \\ \text{(FIT Index)} & & \text{(Frequency)} & & \text{(Intensity)} & & \text{(Time)} \end{array}$$

FIT Index ≥ 60, ATHLETE mode

Number	Frequency of Exercise
5	Daily or almost daily
4	3 to 4 times per week
3	1 to 2 times per week
2	A few times per month
1	Less than once per month

1.6 Taking a Measurement

Intensity	Example Sports Activities	
5	Cycling - > 12 mph pace	Boxing
	Weightlifting - vigorous, powerlifting or bodybuilding	Football - competitive
	Rowing - moderate to vigorous	Handball, racquetball, or squash
	Rowing machines - moderate to vigorous effort	Ice hockey
	Aerobic dancing - high impact	Karate or kickboxing
	Step aerobics	Rockclimbing
	Running -> 5.0 mph	Rugby
	Rope jumping	Soccer - competitive
	Rollerblading (roller skating)	Tennis
	Ski machine	Swimming - competitive or lap
	Stairstepping	Speed skating - competitive
	Stationary cycling - moderate to vigorous effort	Skiing - cross-country
	Basketball - competitive	Skiing - downhill racing

Intensity	Example Sports Activities	
4	Cycling - < 12 mph pace	Rollerblading - leisurely
	Weightlifting - moderate effort	Basketball - shooting baskets
	Stationary cycling - light effort	Bowling
	Rowing - light effort	Golf
	Calisthenics	Gymnastics
	Stretching / Yoga	Horseback riding
	Rowing machines - light effort	Baseball
	Water aerobics or water exercise	Softball
	Aerobic dancing - low impact	Tai Chi
	Jogging - < 5.0 mph	Volleyball - competitive
	Walking - > 2.5 mph	Wrestling - competitive
	Swimming - leisurely	Ice skating - < 9 mph

Time	Duration
4	45 minutes or more
3	30 to 44 minutes
2	15 - 29 minutes
1	less than 15 minutes

1.7 Automatic Transfer of Readings

Requirements

To automatically transfer readings to your smartphone, the following conditions must be met.


■ Smartphone Requirements

- **Bluetooth communication must be enabled**
- **Location must be enabled**

To use automatic transfer, your smartphone's location setting must be enabled.

However, depending on the smartphone you use, you may not need to set the location.

- **The distance from the unit must be within 16.4 ft**

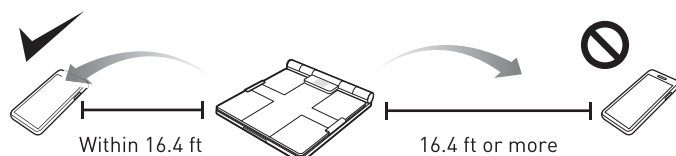- For iPhone, depending on the OS version, automatic transfer may not be possible when the app is running in the background. In this case, transfer the data when the app screen is displayed.

■ The Compatible OMRON App Requirements

- **This product must be registered**

After installing the compatible OMRON app, select this product when registering the device during initial setup.

- **"Automatic data transfer" in the device settings must be enabled**


For details, refer to the "Help" section in the compatible OMRON app.

Automatic Transfer Timing

The timing when the unit automatically transfers readings to your smartphone is as follows.


- **After determining your personal number**

<Personal number display>

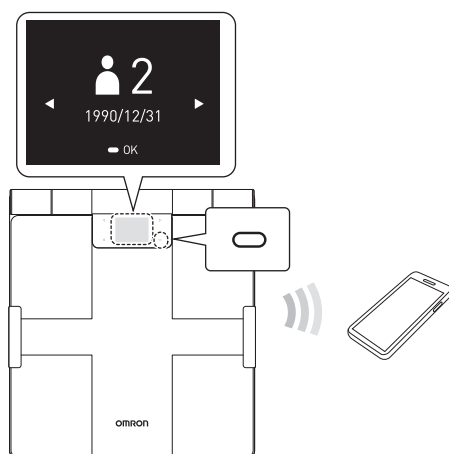- **The power will automatically turn off after approx. 16 seconds of inactivity, even if the personal number has not been identified for the measurement.**

<Personal number and body composition display>
(After 16 seconds)

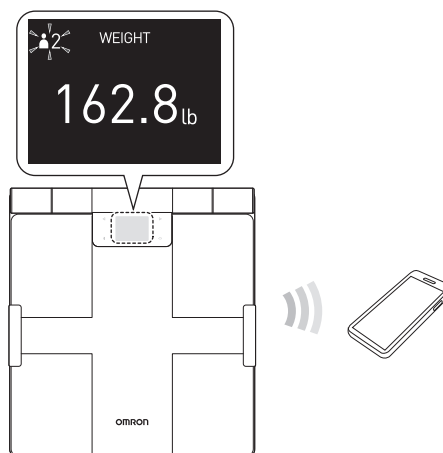Note that if the personal number is incorrect, the readings will be transferred to another registered user's smartphone.

Even if you do not have your smartphone with you when taking a measurement, if you bring your smartphone close to the unit (within 16.4 ft), the readings are automatically transferred.

1.8 Checking Your Readings Using the Compatible OMRON App

You can check the following readings on the compatible OMRON app.

- Body Weight/Body Fat Percentage/Body Fat Mass/Skeletal Muscle Percentage/Skeletal Muscle Mass/Skeletal Muscle Percentage Classification (Arms)/Skeletal Muscle Percentage Classification (Trunk)/Skeletal Muscle Percentage Classification (Legs)/Visceral Fat Level/Resting Metabolism/BMI/Subcutaneous Fat Percentage Classification/Subcutaneous Fat Percentage Classification (Arms)/Subcutaneous Fat Percentage Classification (Trunk)/Subcutaneous Fat Percentage Classification (Legs)/Body Water Percentage/Bone Mass/Fat Free Mass

For more details, refer to sub-section 8.1.

1.9 Measuring Weight Only

When measuring weight only, the readings are not transferred to the smartphone.


It will also not be saved in the unit's memory.

If you want your weight metrics to transfer to the compatible OMRON app, you must take a full-body measurement.

1. Press the power button to turn on the power.

- The weight measurement display will appear.

2. Touch and select “” mode.

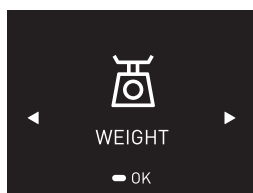3. Touch .

4. Stand still on the unit.

- Leave the display unit in the display unit holder.
- The unit will start measuring your weight.

5. Check the readings.

- After the unit makes a “beeping” sound, the display will show your weight.

6. When the measurement is completed, step off from the unit and turn off the power.

- Even if you forget to turn off the power, it will automatically turn off after approx. 10 seconds.

You can measure your weight simply by stepping on the unit, even with the display unit stored.


1.10 Measuring without Recording (Guest Mode)

When measuring without recording the readings, select “Guest” for the personal number and set the personal data (date of birth, gender, height, mode) before measuring. The readings will not be transferred to your smartphone. It will also not be saved in the unit’s memory.

1. Press the power button to turn on the power.

- The weight measurement display will appear.

2. Touch to select “Guest”.

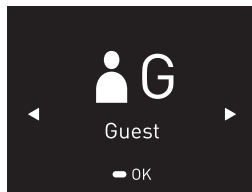3. Touch to confirm.

4. Set up your personal data.

	Set the year of your date of birth • Touch  to change the year. Each time you touch the button, the value will change in increments of one year, and if you continue to touch the button, the value will change in increments of 10 years. • Touch  to set the year.
--	---

	Set the month of your date of birth • Touch  to change the month. Each time you touch the button, the value will change in increments of one month, and if you continue to touch the button, the value will continue to change in one-month increments, but will increase at a faster rate. • Touch  to set the month.
--	---

	Set the day of your date of birth • Touch  to change the day. Each time you touch the button, the value will change in increments of one day, and if you continue to touch the button, the value will change in increments of 10 days. • Touch  to set the day.
--	--

	Set your gender • Touch  to select male or female. • Touch  to set the gender.
--	---

	Set your height • Touch  to change the height. The setting range is 3' 4" to 6' 6 1/2". Each time you touch the button, the value will change in 1/4" increments, and if you continue to touch the button, the value will change in 1" increments. • Touch  to set your height.
--	--

	Set the mode • Touch  to select normal or athlete. • Touch  to set the mode.
--	---

- Once the settings are complete, the configured information will be displayed, followed by “0.0 lb” (weight measurement).

--	--	--

5. Take a measurement.

- Refer to sub-section 1.6.

Note

- The measurement results are shown as 18 indicators on the display.


1.11 Selecting the Personal Number Manually Before Measurement

If your personal number is often not displayed on the unit display when you take a measurement, follow the procedure below to select your personal number and take a measurement.

1. Press the power button to turn on the power.

- The weight measurement display will appear.

2. Touch to select your personal number.

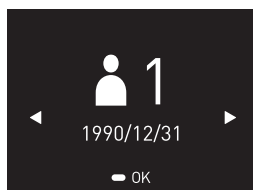3. Touch to confirm your personal number.

- The weight measurement display appears.

4. Take a measurement.

- Refer to sub-section 1.6.

2. Settings

2.1 User Settings/Display Settings/Deleting Data

Use this operation when you want to change your date of birth (year, month, day), gender, height, and mode in your personal data, display settings, or delete personal data on the unit.

You can also change your personal data using the compatible OMRON app. For details, refer to the “Help” section in the compatible OMRON app.

User settings


The readings will not be deleted.

1. Press the power button to turn on the power.

- The weight measurement display will appear.


2. Touch to select your personal number.

- Select the personal number you want to change or delete.

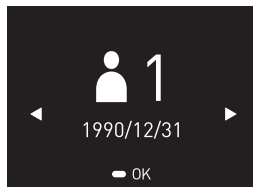
3. Touch to confirm.

- The weight measurement display appears.

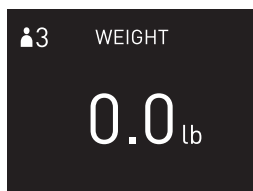4. Touch again.

5. Touch to select “User settings”, then touch to confirm.

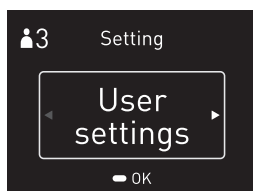- Even if you change the personal data, the readings will not be deleted. When using a personal number that has been used by another person, delete the personal data and then register your personal data again.

(If you use the compatible OMRON app, refer to sub-section 1.4. If you do not use the compatible OMRON app, refer to sub-section 2.5.)

- Changes made on the main unit or the compatible OMRON app will be reflected during communication.

6. Set up your personal data.

	Set the year of your date of birth • Touch  to change the year. Each time you touch the button, the value will change in increments of one year, and if you continue to touch the button, the value will change in increments of 10 years. • Touch  to set the year.
--	---

	Set the month of your date of birth • Touch  to change the month. Each time you touch the button, the value will change in increments of one month, and if you continue to touch the button, the value will continue to change in one-month increments, but will increase at a faster rate. • Touch  to set the month.
--	---

	Set the day of your date of birth • Touch  to change the day. Each time you touch the button, the value will change in increments of one day, and if you continue to touch the button, the value will change in increments of 10 days. • Touch  to set the day.
--	--

	Set your gender • Touch  to select male or female. • Touch  to set the gender.
--	---


2.1 User Settings/Display Settings/Deleting Data

Set your height


- Touch   to change the height.
The setting range is 3' 4" to 6' 6 1/2".
Each time you touch the button, the value will change in 1/4" increments, and if you continue to touch the button, the value will change in 1" increments.
- Touch  to set your height.

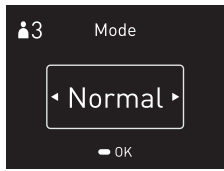
Set the mode

- Touch   to select normal or athlete.
- Touch  to set the mode.

- Once the settings are complete, the configured information will be displayed, followed by "0.0 lb" (weight measurement).

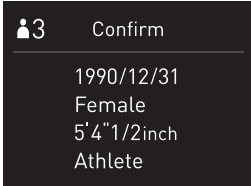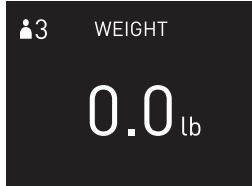7. Take a measurement.

- Refer to sub-section 1.6.


Display settings

You can select either the 4 indicators or the 18 indicators.

1. Select your personal number.

- Refer to steps 1 to 4 of the "User settings" in this section.


2. Touch to select "Display settings", then touch to confirm.

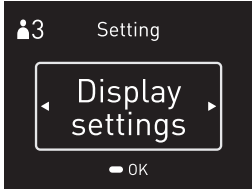
- Changes made on the main unit or the compatible OMRON app will be reflected during communication.


3. Set display indicators.


- Touch   to change indicators.
- Touch  to confirm.
- After 3 seconds, the "0.0 lb" is displayed.

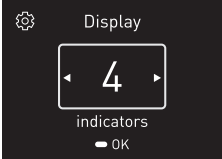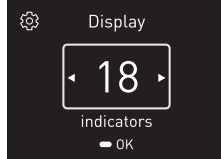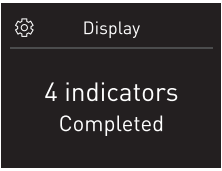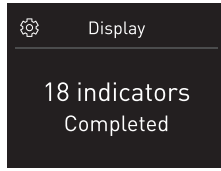
Deleting data

1. Select your personal number.


- Refer to steps 1 to 4 of the "User settings" in this section.


2. Touch to select "Delete data", then touch to confirm.

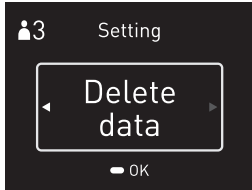- Changes made on the main unit or the compatible OMRON app will be reflected during communication.

3. Touch and hold the for more than 3 seconds.

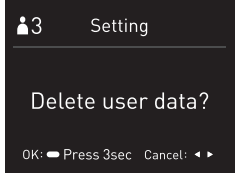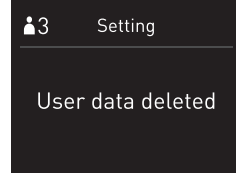- To cancel the setting, touch .
- After 10 seconds, the unit will automatically turn off.

2.2 Replacing the Batteries

When to Replace the Batteries

When the  and "A1" are displayed, the battery is depleted.

Replace all four at the same time with new AA alkaline batteries of the same brand.

Replace the Batteries

1. Remove the depleted batteries from the unit.

2. Insert new batteries.


- Refer to sub-section 1.3.

3. Perform the initial settings.

■ If you are using the compatible OMRON app

Touch and hold the .

- Communication with the smartphone will begin. Initial setting information (date and time) will be transferred to the unit.
- The following will be displayed on the unit's display. When the transfer is complete, the unit will make a "beeping" sound. The unit will automatically turn off after 10 seconds.

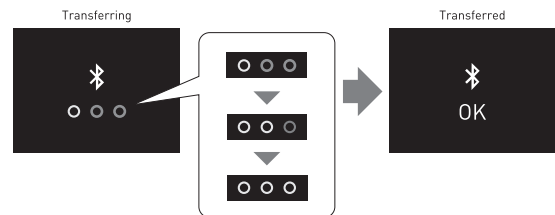■ If you are not using the compatible OMRON app
Perform the initial settings (date and time) on the unit. Refer to sub-section 2.4.

- Do not use manganese batteries or rechargeable batteries as they may significantly reduce the number of measurements or cause malfunction.
- Dispose of the replaced batteries according to the guidelines of your local regulation.
- When the battery is low, the  will flash. Prepare new batteries.
- Do not replace the battery while your readings are being transferred to your smartphone. This may result in the incorrect operation of your unit and failure to transfer your readings.


2.3 Restoring the Unit to Factory Settings

If you want to restore your unit to the factory settings, follow the steps below.

This will delete all registration information (date and time settings, personal data, communication settings) and all readings stored in the unit's memory.

1. Touch and hold the and for more than 3 seconds.

2. Touch and hold the and again for more than 3 seconds.

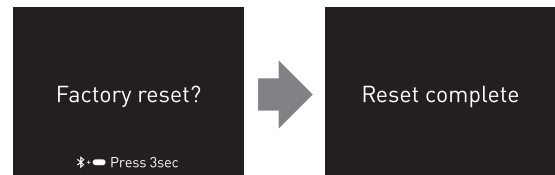
- After 2 seconds, the unit will automatically turn off.
- To cancel restoring, press the power button to turn off the power.

2.4 Performing Initial Settings on the Unit

You can configure the initial settings on the main unit immediately after purchase or after replacing its batteries.

1. Press the power button to turn on the power.

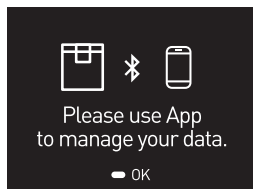2. Touch .

3. Set the measurement unit, date, and time.

	Set the measurement unit • Touch  to select lb or kg. This will also automatically change inch to cm. • Touch  to set the measurement unit.
--	--

	Set the year • Touch  to change the year. Each time you touch the button, the value will change in increments of one year, and if you continue to touch the button, the value will change in increments of 10 years. • Touch  to set the year.
--	---

	Set the month • Touch  to change the month. Each time you touch the button, the value will change in increments of one month, and if you continue to touch the button, the value will continue to change in one-month increments, but will increase at a faster rate. • Touch  to set the month.
--	---

	Set the day • Touch  to change the day. Each time you touch the button, the value will change in increments of one day, and if you continue to touch the button, the value will change in increments of 10 days. • Touch  to set the day.
--	--


	Set the hour • Touch  to change the hour. Each time you touch the button, the value will change in increments of one hour, and if you continue to touch the button, the value will change in increments of 10 hours. • Touch  to set the hour.
--	---

	Set the minute • Touch  to change the minute. Each time you touch the button, the value will change in increments of one minute, and if you continue to touch the button, the value will change in increments of 10 minutes. • Touch  to set the minute.
--	---

- Once the settings are complete, the configured information will be displayed, then the power will automatically turn off.

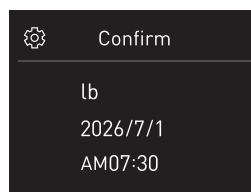- If you make a mistake while setting the numbers, turn off the power, then start over from step 1.
- To change the values after setting the date and time, follow the steps below.
 - (1) Remove the batteries
 - (2) Wait for more than 20 seconds, then insert the batteries again
 - (3) Repeat from step 1
- If you do not perform any operation for more than 3 minutes during the setup, the unit will automatically turn off. Start over from step 1.

2.5 Registering Personal Data on the Unit

You can also register your personal data (personal number, date of birth, gender, height, mode) on the unit itself.

The maximum number of personal data entries that can be registered is 4.


1. Press the power button to turn on the power.

- The weight measurement display will appear.

2. Touch to select the personal number you want to register.

Example:

If you select personal number "3"

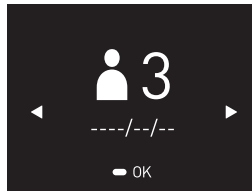3. Touch to set.

	Set the year of your date of birth • Touch  to change the year. Each time you touch the button, the value will change in increments of one year, and if you continue to touch the button, the value will change in increments of 10 years. • Touch  to set the year.
--	--

	Set the month of your date of birth • Touch  to change the month. Each time you touch the button, the value will change in increments of one month, and if you continue to touch the button, the value will continue to change in one-month increments, but will increase at a faster rate. • Touch  to set the month.
--	---


	Set the day of your date of birth • Touch  to change the day. Each time you touch the button, the value will change in increments of one day, and if you continue to touch the button, the value will change in increments of 10 days. • Touch  to set the day.
--	--

	Set your gender • Touch  to select male or female. • Touch  to set the gender.
--	---

	Set your height • Touch  to change the height. The setting range is 3' 4" to 6' 6 1/2". Each time you touch the button, the value will change in 1/4" increments, and if you continue to touch the button, the value will change in 1" increments. • Touch  to set your height.
--	--

	Set the mode • Touch  to select normal or athlete. • Touch  to set the mode.
--	---

- Once the settings are complete, the configured information will be displayed, followed by "0.0 lb" (weight measurement).

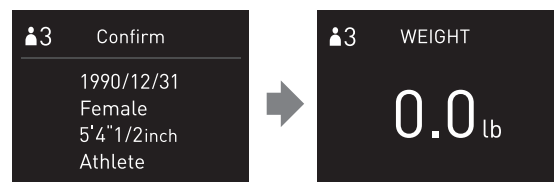4. Take a measurement.






- Refer to sub-section 1.6.
- Your weight and body composition data will be saved to estimate your personal number at the next measurement. The readings are saved in the unit's memory.




This will complete the registration of your personal data.

- To register other personal data, follow the procedure from step 1.
- To change or delete your personal data, refer to sub-section 2.1.

3. Error Messages and Troubleshooting

3.1 Error Messages

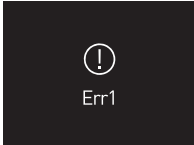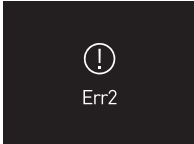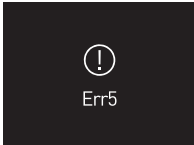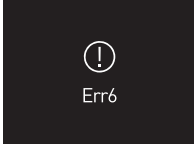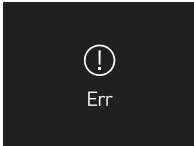Error Display	Causes and Solutions
	Got off the unit while your body composition was being measured. - Do not step off the unit until the measurement of body composition is complete. You moved during the measurement. - Stand still during measurement.
	The palms or soles are not in firm contact with the electrodes. - Place the palms or soles in close contact with the electrodes when taking measurements. Refer to sub-section 1.5. The palms or soles were dry. - Slightly dampen the palms or soles with a wet towel and try again. Measurements are taken with socks or stockings. - Measure with bare feet.
	Internal processing failed. Remove the batteries for more than 20 seconds, then install the batteries again and perform the initial settings. Then, touch .
	
	You moved the unit before "0.0 lb" was displayed. - Do not move the unit until the measurement is complete. Moved your body during weight measurement. - Stand still during measurement. Body weight exceeded 330.0 lb.
	Communication is not established properly. Check the display on your smartphone and follow the instructions in the compatible OMRON app, or refer to the "Help" section in the compatible OMRON app.

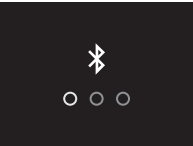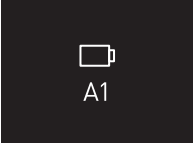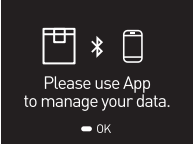Error Display	Causes and Solutions
	Personal data cannot be identified. When your personal data is not displayed and a "---/--/--" is shown instead, touch  and select the set personal number.
	Continued to touch . This is displayed when registering this product in the compatible OMRON app. When registering, follow the instructions in the compatible OMRON app. If you are not registering, you can cancel the operation by pressing the power button for more than 2 seconds.
	Touched . This is displayed when transferring readings to your smartphone. When transferring readings, follow the instructions in the compatible OMRON app. If you do not transfer readings, you can cancel the operation by pressing the power button for more than 2 seconds.
	More than 24 sets of readings are not transferred.  will disappear when you transfer the readings to the compatible OMRON app. 30 sets of readings are not transferred.  will disappear when you transfer the readings to the compatible OMRON app. If you continue measuring without transferring the data, the oldest readings will be deleted.
	Batteries are running low. Replace all four at the same time with new AA alkaline batteries of the same brand.
	Batteries are depleted. Replace all four at the same time with new AA alkaline batteries of the same brand.
	The initial settings were not performed after replacing the batteries. - Open the compatible OMRON app and touch  to communicate with your smartphone. - To perform the initial settings on the unit, refer to sub-section 2.4.

3.2 Troubleshooting

In case any of the below problems occur during use, first check that no other electrical device is within 12 inches (30 cm). If the problem persists, please refer to the table below.

Problem	Causes and Solutions
Nothing is displayed when the unit is turned on.	No batteries are installed. - Install the batteries. Battery $\oplus$ $\ominus$ polarities are installed in the wrong direction. - Install the batteries in the correct direction.
Even if you step on the unit, nothing is displayed.	The batteries are depleted. Replace all four at the same time with new AA alkaline batteries of the same brand.
The compatible OMRON app only displays weight and BMI.	The person has not yet reached the age for measurement. For children who are 5 years old or younger, only weight and BMI values are displayed. If the target person is 17 years old or younger, the visceral fat level, skeletal muscle percentage classification and subcutaneous fat classification will not be displayed.
The compatible OMRON app does not display the visceral fat level, skeletal muscle percentage classification and subcutaneous fat classification.	The date of birth or date setting is incorrect. - Open the compatible OMRON app and touch  to communicate with your smartphone. - To change your personal data on the unit, refer to sub-section 2.1. - If the date is incorrect, remove the batteries, wait for approx. 20 seconds, then install them again and perform the initial settings.

Problem	Causes and Solutions
Some readings are displayed as "----".	The registered values, weight, and body composition values were outside the measurable range. - Check if your personal data (date of birth, gender, height, mode) set in the compatible OMRON app is correct. - Even if the date of birth, gender, height, and mode are correct, "----" will be displayed if they are out of the display range or supported age range. Refer to section 5.
The correct personal number is not identified when you step on the unit.	Your body weight has changed significantly since the last measurement. People with similar builds are registered. Another person's readings were recorded. Select your personal number before taking measurements. Refer to sub-section 1.11.
Only weight, body fat percentage, skeletal muscle percentage, and BMI are measured.	- You can only check your weight, body fat percentage, skeletal muscle percentage and BMI on the unit. You can check indicators other than weight, body fat percentage, skeletal muscle percentage and BMI on the compatible OMRON app. Install the compatible OMRON app. - You can change the display indicators in sub-section 2.1, "Display settings".

3.2 Troubleshooting

Problem	Causes and Solutions
High weight and body composition measurements	Your posture is wrong during measurement. Take measurements with a correct posture. Refer to sub-section 1.5.
Low measured values	Measurements are taken on soft surfaces such as carpets or on uneven floor surfaces. Measure on a hard, flat surface. Refer to sub-section 1.5, "Install the Unit".
The values are extremely different for each measurement.	The soles of your feet and body are cold, impairing blood circulation. - Warm up your body to get blood circulation back to normal before taking a measurement. The soles of your feet are dry. - Slightly dampen the soles of your feet with a wet towel and then measure. The "0 lb calibration" is not implemented correctly. - Implement the "0 lb calibration". Refer to sub-section 1.5, "Implement the '0 lb Calibration'". Your date of birth, gender, height, or mode is not set correctly. - Check if your personal data set in the compatible OMRON app is correct.

Problem	Causes and Solutions
None of the touch buttons respond.	You touched a button when you stepped on the unit. - Do not operate the touch buttons while standing on the unit. You touched two or more buttons simultaneously. - Touch one button at a time. Your fingertips are dry. - Slightly dampen your fingertips with a wet towel before operating the touch buttons. No batteries are installed. - Install the batteries. Refer to sub-section 1.3.
Cannot transfer data to your smartphone.	Refer to the "Help" section in the compatible OMRON app.

- Rarely, an error could occur because of a person's physical constitution, making measurement impossible.
- All readings will be deleted if this product malfunctions or is repaired.

If the above corrections are not effective, contact us at 1-800-634-4350.

4. Maintenance


4.1 Maintenance

- Always keep the unit clean.
- If the device is visibly soiled or has been used by a person with a skin disease or foot infection, clean the device using a soft cloth moistened with water or a neutral detergent. Ensure that the cloth is well-wrung. Wipe thoroughly, ensuring no excess moisture enters the device, and then dry it completely with a clean cloth. (Otherwise, a bacterial infection can be caused.)
- Do not use benzine, thinner, or other volatile solvents to wipe the unit. (It may cause discoloration or damage.)
- Do not wash the unit with water. (It may cause the unit to malfunction.)
- You can use cleaning alcohol to clean the electrodes, but do not use it on other parts of the unit.

4.2 Storage

- Store the display unit in the main unit as shown.
- When storing the display unit, be careful not to crease the cord.

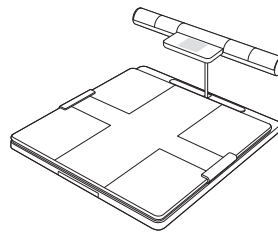Do not store the unit in the following locations. It may cause the unit to malfunction.


- Places where water may get into the unit.
- Places with high temperatures, high humidity, direct sunlight, dust, or places with high salinity.
- Places with sloping surfaces, vibrations, or shocks.
- Places where chemicals are stored or where corrosive gas is present.

Do not store the unit in the following conditions.


The readings displayed may differ from the actual results. Implement the "0 lb calibration" before taking a measurement.

Refer to sub-section 1.5, "Implement the '0 lb Calibration'".


- The unit is leaning against a wall or other objects

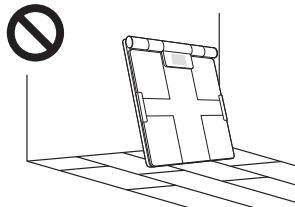- The unit is placed upside down

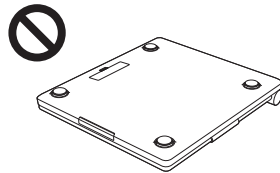- The unit is placed on an object

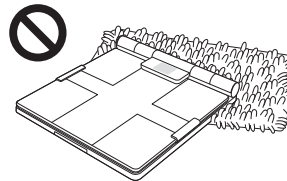- An object is placed on the unit

Disposal

- Remove the batteries from the unit and dispose of the unit.
- Dispose of the unit and batteries according to applicable local regulations. Unlawful disposal may cause environmental pollution.

5. Specifications

5.1 Specifications

Product name	VIVA 360+ Smart Scale
Model	BCM-750 REF HBF-H122T1-ZBK
Display*	
Body weight	6.6 to 330.0 lb, increments of 0.2 lb (3.0 to 150.0 kg, increments of 0.1 kg)
BMI	7.0 to 90.0, increments of 0.1
BMI classification	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels
Body fat percentage	Measuring range : 5.0 to 60.0 %, increments of 0.1 % Age range : 6 to 80 years old
Body fat percentage classification	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels Age range : 6 to 80 years old
Body fat mass	Measuring range : 0.4 to 198.0 lb, increments of 0.2 lb (0.2 to 90.0 kg, increments of 0.1 kg) Age range : 6 to 80 years old
Skeletal muscle percentage (whole body)	Measuring range : 5.0 to 50.0 %, increments of 0.1 % Age range : 6 to 80 years old
Skeletal muscle percentage classification (whole body)	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels Age range : 18 to 80 years old
Skeletal muscle mass (whole body)	Measuring range : 0.4 to 165.0 lb, increments of 0.2 lb (0.2 to 75.0 kg, increments of 0.1 kg) Age range : 6 to 80 years old
Resting metabolism	Measuring range : 385 to 3999 kcal, increments of 1 kcal Age range : 6 to 80 years old
Visceral fat level	Measuring range : 30 levels, increments of 1 level Age range : 18 to 80 years old
Visceral fat level classification	0 (Normal) / + (High) / ++ (Very High) 3 levels Age range : 18 to 80 years old
Skeletal muscle percentage classification (segmental)	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels Age range : 18 to 80 years old
Subcutaneous fat percentage classification (whole body)	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels Age range : 18 to 80 years old
Subcutaneous fat percentage classification (segmental)	- (Low) / 0 (Normal) / + (High) / ++ (Very High) 4 levels Age range : 18 to 80 years old
Body water percentage	Measuring range : 20.0 to 80.0 %, increments of 0.1 % Age range : 6 to 80 years old
Body water percentage classification	- (Low) / 0 (Normal) / + (High) 3 levels Age range : 6 to 80 years old
Bone mass	Measuring range : 1.10 to 17.64 lb, increments of 0.2 lb (0.5 to 8.0 kg, increments of 0.1 kg) Age range : 6 to 80 years old
Bone mass classification	- (Low) / 0 (Normal) / + (High) 3 levels Age range : 6 to 80 years old
Fat free mass	Measuring range : 2.8 to 314.0 lb, increments of 0.2 lb (1.2 to 142.5 kg, increments of 0.1 kg) Age range : 6 to 80 years old
Memory	30 readings (Total 120 readings for 4 people)
Transmission protocol	Bluetooth®
Wireless communication	Frequency range : 2.4 GHz (2400 - 2483.5 MHz) Modulation : GFSK Effective radiated power : < 20 dBm

Weight accuracy	6.6 to 88.2 lb : ± 0.88 lb (3.0 to 40.0 kg : ± 0.4 kg) 88.2 to 330.0 lb : ± 1 % (40.0 to 150.0 kg : ± 1 %)
Durable period	5 years or the time when reaching 7300 times of use.
IP Classification	IP21* *IP classification is degrees of protection provided by IEC 60529. This unit is protected against solid foreign objects of diameter 12.5 mm such as a finger and greater. This unit is protected against the ingress of vertically falling water drops.
Power supply	4 AA alkaline batteries (LR6)
Battery life	Approx. 6 months (When using AA alkaline batteries with four measurements per day under OMRON's test conditions.)
Operating conditions	+41 to +104 °F (+5 to +40 °C), 30 to 85 %RH (non-condensing), 860 - 1060 hPa
Storage / Transport conditions	-4 to +140 °F (-20 to +60 °C), 10 to 95 %RH (non-condensing)
Weight	Approx. 5.41 lbs (2.5 kg) (including batteries)
External dimensions	Approx. 12.0" (W) x 2.1" (H) x 13.1" (D) (304 x 53.5 x 332 mm)
Contents	Main unit, 4 AA alkaline batteries (LR6), instruction manual
Applied part	Type BF (grip electrodes, foot electrodes)
Protection against electric shock	Internally powered ME equipment
Operation mode	Continuous operation

* All indicators other than weight (lb) are computed by calculation. All calculated values, except BMI, are estimated using our proprietary formula.

- Subject to technical modification without prior notice.
- This OMRON product is produced under the strict quality system of OMRON HEALTHCARE Co., Ltd., Japan.
- Communication with a smart device is paired and encrypted according to Bluetooth Low Energy specifications. Pairing requires user interaction.
- Operation mode classification complies with IEC 60601-1.

Trademarks

- The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by OMRON HEALTHCARE Co., Ltd. is under license.
- Other brand names and product names mentioned in this document are trademarks or registered trademarks of the respective companies.

5.2 Important Information regarding Electromagnetic Compatibility (EMC)

This device conforms to IEC 60601-1-2:2014+A1:2020 Electromagnetic Compatibility (EMC) standard.
Further documentation in accordance with this EMC standard is available at:
<https://healthcare.omron.com/electro-magnetic-compatibility>.
Refer to the EMC information for this device on the website.

6. Limited Warranty

Your BCM-750 VIVA 360+ Smart Scale, excluding the batteries, is warranted to be free from defects in materials and workmanship appearing within 1 year from the date of purchase, when used in accordance with the instructions provided with the unit. The above warranty extends only to the original retail purchaser, and only to products purchased from an Omron authorized seller who is subject to and follows Omron's quality control standards, unless otherwise prohibited by law. We will, at our option, replace without charge any unit covered by the above warranty. Replacement is our only responsibility and your only remedy under the above warranty.

To obtain warranty service contact Customer Service by calling **1-800-634-4350** for the address of the inspection center and the return shipping and handling fee that may apply.

Enclose the original printed receipt. Include a letter, with your name, address, phone number, and description of the specific problem. Pack the product carefully to prevent damage in transit. Because of possible loss in transit, we recommend insuring the product with return receipt requested.

THE FOREGOING IS THE SOLE WARRANTY PROVIDED BY OMRON IN CONNECTION WITH THIS PRODUCT, AND OMRON HEREBY DISCLAIMS ANY OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NON-INFRINGEMENT. IMPLIED WARRANTIES AND OTHER TERMS THAT MAY BE IMPOSED BY LAW, IF ANY, ARE LIMITED IN DURATION TO THE PERIOD OF THE ABOVE EXPRESS WARRANTY.

OMRON SHALL NOT BE LIABLE FOR LOSS OF USE OR ANY OTHER SPECIAL, INCIDENTAL, CONSEQUENTIAL OR INDIRECT COSTS, EXPENSES OR DAMAGES.

SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

This warranty provides you with specific legal rights, and you may have other rights that vary by jurisdiction. Because of special local requirements, some of the above limitations and exclusions may not apply to you.

FOR CUSTOMER SERVICE

Visit our web site at:

OmronHealthcare.com

Call toll free:

1-800-634-4350

7. FCC/ISED Statement

FCC Compliance statement

CAUTION

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note

- This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.

This transmitter must not be co-located or operated in conjunction with any other antenna or transmitter.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment and meets the FCC radio frequency (RF) Exposure Guidelines as this equipment has very low levels of RF energy.

ISED Canada compliance statement

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with ISED radiation exposure limits set forth for an uncontrolled environment and meets the RSS-102 of the ISED radio frequency (RF) Exposure rules as this equipment has very low levels of RF energy.

8. Information on Body Composition

8.1 Information on Body Composition

Principle of body composition calculation

Body fat has low electric conductivity

The unit measures the body fat percentage by the Bioelectrical Impedance (BI) method. Muscles, blood vessels and bones are body tissues with a high water content that conducts electricity easily. Body fat is tissue that has little electric conductivity. The unit sends an extremely weak electrical current of 50 kHz and less than 500 µA through your body to determine the amount of fat tissue. This weak electrical current is not felt while operating the unit.

In order for the unit to determine body composition, it uses the electrical impedance, along with your height, weight, age and gender information to generate results based on OMRON's data of body composition.


Measures the whole body to avoid the influence of fluctuations

During the course of a day, the amount of water in the body tends to gradually shift to the lower limbs. This is why there is a tendency for the legs and ankles to swell in the evening or at night. The ratio of water in the upper body and lower body is different in the morning and evening, and this means that the electrical impedance of the body also varies. Since the unit uses electrodes for both hands and feet to take measurements, it can reduce the influence of these fluctuations on measurement results.

Recommended measurement times

Understanding the normal changes in your body fat percentage can help you in preventing or reducing obesity.

Being aware of the times when the body fat percentages shift within your own daily schedule will assist you in obtaining an accurate trending of your body fat. It is recommended to use this unit in the same environment and daily circumstances. [See chart]

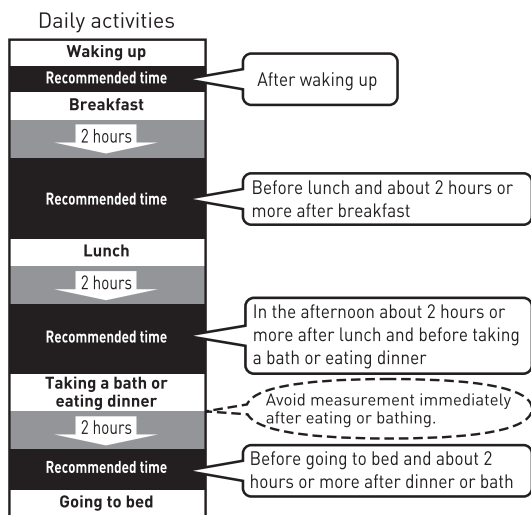Avoid taking measurements under the following conditions:

- Immediately after vigorous exercise, after a bath or sauna.
- After drinking alcohol or a large amount of water, after a meal (about 2 hours).

If a measurement is taken under these physical conditions, the calculated body composition may differ significantly from the actual one because the water content in the body is changing.

What is BMI (Body Mass Index)?

BMI uses the following simple formula to indicate the ratio between weight and height of a person.

$$\text{BMI} = \text{weight (lb)} / \text{height (inches)} / \text{height (inches)} \times 703$$

The unit uses the height information stored in your personal profile number or when entering information in the Guest Mode to calculate your BMI classification.

If the fat level revealed by BMI is higher than the international standard, there is an increased likelihood of common diseases. However, not all types of fat can be revealed by BMI.

8.1 Information on Body Composition

What is Body Fat Percentage?

Body fat percentage refers to the amount of body fat mass in regards to the total body weight expressed as a percentage.

$$\text{Body fat percentage (\%)} = \{\text{Body fat mass (lb)} / \text{Body weight (lb)}\} \times 100$$


The unit uses BI method to estimate your body fat percentage.

Depending on where fat is distributed in the body, it is classified as visceral fat or subcutaneous fat.

What is Visceral Fat Level?

Visceral fat = fat surrounding internal organs

Too much visceral fat is thought to be closely linked to increased levels of fat in the bloodstream, which can lead to common diseases such as hyperlipidemia and diabetes, which impairs the ability of insulin to transfer energy from the bloodstream to the cells for use. In order to prevent or improve conditions of common diseases, it is important to try and reduce visceral fat levels to an acceptable level. People with high visceral fat levels tend to have large stomachs. However, this is not always the case and high visceral fat levels can lead to being metabolically obese. Metabolically obese [visceral obesity with normal weight] represents fat levels that are higher than average, even if a person's weight is at or below the standard for their height.

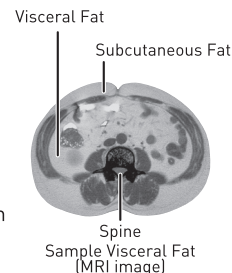
What is Subcutaneous Fat?

Subcutaneous fat = fat below the skin

Subcutaneous fat not only accumulates around the stomach but also around the upper arms, hips and thighs, and can cause a distortion of the body's proportions.

Although not directly linked to increased risk of disease, it is thought to increase pressure on the heart and other complications. Subcutaneous fat is not displayed in this unit, but is included in the body fat percentage.

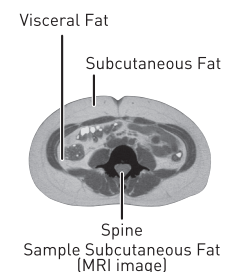What is Skeletal Muscle?

Muscle is divided into two types, muscle in internal organs, such as the heart, and skeletal muscle attached to bones that is used to move the body. Skeletal muscle can be increased through exercise and other activity.

Increasing the ratio of skeletal muscle means that the body can burn energy more easily, which means that it is less likely to turn to fat, and makes it easier to lead an energetic lifestyle.

What is Resting Metabolism?

Regardless of your activity level, a minimum level of caloric intake is required to sustain the body's everyday functions. Known as the resting metabolism, this indicates how many calories you need to ingest in order to provide enough energy for your body to function.

The reason calculated results may differ from actual body fat percentage

The body fat percentage measured by this unit may significantly differ from the actual body fat percentage in the following situations:

Elderly people (over 81 years old) / People with a fever / Body builders or highly trained athletes / Patients undergoing dialysis / Patients with osteoporosis who have very low bone density / Pregnant women / People with swelling

These differences may be related to changing ratios of body fluid and/or body composition.

8.2 Measurement Results Interpretation Diagrams

Interpreting the Body Fat Percentage Result

The table below is based on research by HD McCarthy et al, in the International Journal of Obesity, Vol. 30, 2006, and by Gallagher et al., American Journal of Clinical Nutrition, Vol. 72, Sept. 2000, and classified into four levels by Omron Healthcare.

Gender	Age	- (Low)	0 (Normal)	+ (High)	++ (Very High)
Female	6	< 13.8 %	13.8 - 24.9 %	25.0 - 27.0 %	≥ 27.1 %
	7	< 14.4 %	14.4 - 27.0 %	27.1 - 29.6 %	≥ 29.7 %
	8	< 15.1 %	15.1 - 29.1 %	29.2 - 31.9 %	≥ 32.0 %
	9	< 15.8 %	15.8 - 30.8 %	30.9 - 33.8 %	≥ 33.9 %
	10	< 16.1 %	16.1 - 32.2 %	32.3 - 35.2 %	≥ 35.3 %
	11	< 16.3 %	16.3 - 33.1 %	33.2 - 36.0 %	≥ 36.1 %
	12	< 16.4 %	16.4 - 33.5 %	33.6 - 36.3 %	≥ 36.4 %
	13	< 16.4 %	16.4 - 33.8 %	33.9 - 36.5 %	≥ 36.6 %
	14	< 16.3 %	16.3 - 34.0 %	34.1 - 36.7 %	≥ 36.8 %
	15	< 16.1 %	16.1 - 34.2 %	34.3 - 36.9 %	≥ 37.0 %
	16	< 15.8 %	15.8 - 34.5 %	34.6 - 37.1 %	≥ 37.2 %
	17	< 15.4 %	15.4 - 34.7 %	34.8 - 37.3 %	≥ 37.4 %
	18 - 39	< 21.0 %	21.0 - 32.9 %	33.0 - 38.9 %	≥ 39.0 %
	40 - 59	< 23.0 %	23.0 - 33.9 %	34.0 - 39.9 %	≥ 40.0 %
	60 - 80	< 24.0 %	24.0 - 35.9 %	36.0 - 41.9 %	≥ 42.0 %
Male	6	< 11.8 %	11.8 - 21.7 %	21.8 - 23.7 %	≥ 23.8 %
	7	< 12.1 %	12.1 - 23.2 %	23.3 - 25.5 %	≥ 25.6 %
	8	< 12.4 %	12.4 - 24.8 %	24.9 - 27.7 %	≥ 27.8 %
	9	< 12.6 %	12.6 - 26.5 %	26.6 - 30.0 %	≥ 30.1 %
	10	< 12.8 %	12.8 - 27.9 %	28.0 - 31.8 %	≥ 31.9 %
	11	< 12.6 %	12.6 - 28.5 %	28.6 - 32.6 %	≥ 32.7 %
	12	< 12.3 %	12.3 - 28.2 %	28.3 - 32.4 %	≥ 32.5 %
	13	< 11.6 %	11.6 - 27.5 %	27.6 - 31.3 %	≥ 31.4 %
	14	< 11.1 %	11.1 - 26.4 %	26.5 - 30.0 %	≥ 30.1 %
	15	< 10.8 %	10.8 - 25.4 %	25.5 - 28.7 %	≥ 28.8 %
	16	< 10.4 %	10.4 - 24.7 %	24.8 - 27.7 %	≥ 27.8 %
	17	< 10.1 %	10.1 - 24.2 %	24.3 - 26.8 %	≥ 26.9 %
	18 - 39	< 8.0 %	8.0 - 19.9 %	20.0 - 24.9 %	≥ 25.0 %
	40 - 59	< 11.0 %	11.0 - 21.9 %	22.0 - 27.9 %	≥ 28.0 %
	60 - 80	< 13.0 %	13.0 - 24.9 %	25.0 - 29.9 %	≥ 30.0 %

8.2 Measurement Results Interpretation Diagrams

Interpreting the Visceral Fat Level Result

Visceral Fat Level	Level Classification
1 - 9	0 (Normal)
10 - 14	+ (High)
15 - 30	++ (Very High)

Classified based on OMRON HEALTHCARE Co., Ltd. study data.

Interpreting the Skeletal Muscle Percentage Result

Gender	Age	- (Low)	0 (Normal)	+ (High)	++ (Very High)
Female	18 - 39	< 24.3 %	24.3 - 30.3 %	30.4 - 35.3 %	≥ 35.4 %
	40 - 59	< 24.1 %	24.1 - 30.1 %	30.2 - 35.1 %	≥ 35.2 %
	60 - 80	< 23.9 %	23.9 - 29.9 %	30.0 - 34.9 %	≥ 35.0 %
Male	18 - 39	< 33.3 %	33.3 - 39.3 %	39.4 - 44.0 %	≥ 44.1 %
	40 - 59	< 33.1 %	33.1 - 39.1 %	39.2 - 43.8 %	≥ 43.9 %
	60 - 80	< 32.9 %	32.9 - 38.9 %	39.0 - 43.6 %	≥ 43.7 %

Classified based on OMRON HEALTHCARE Co., Ltd. study data.

Interpreting the BMI Result

BMI	BMI (Designation by the WHO)
< 18.5	- (Underweight)
18.5 - 24.9	0 (Normal)
25 - 29.9	+ (Overweight)
≥ 30	++ (Obese)

Source: Values for obesity judgment proposed by WHO, the World Health Organization.

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© 2026 OMRON HEALTHCARE, INC.
Classé d'après des données d'étude d'OMRON HEALTHCARE Co., Ltd.