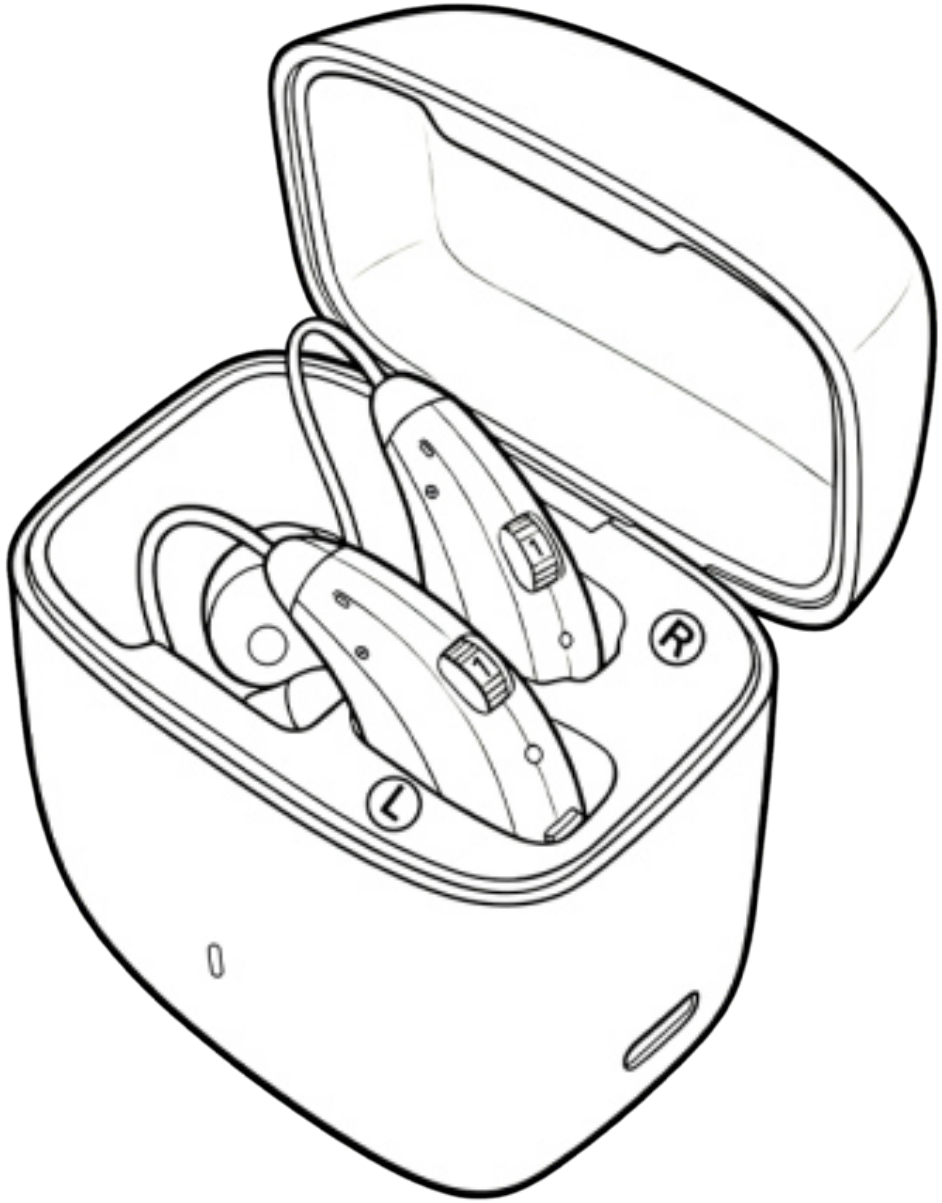
VHP-1735i

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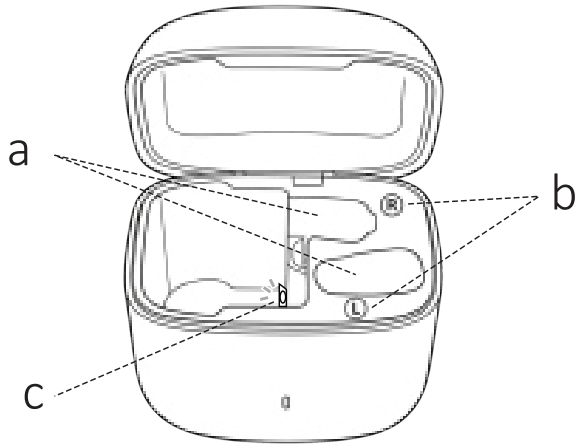
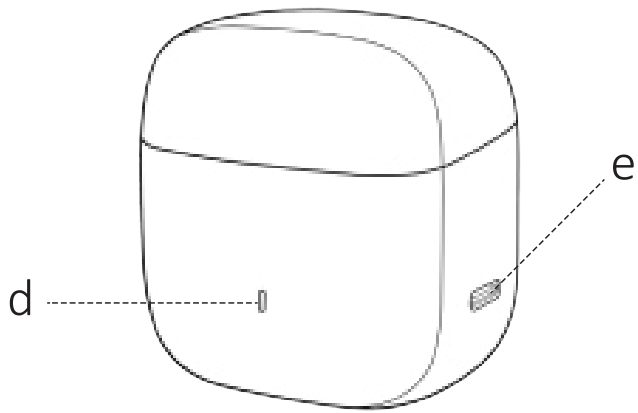
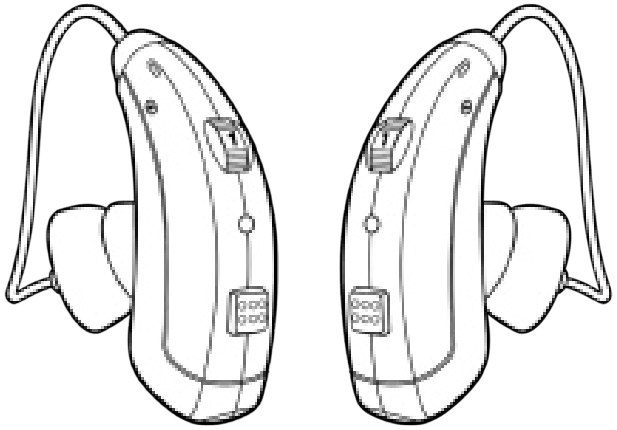
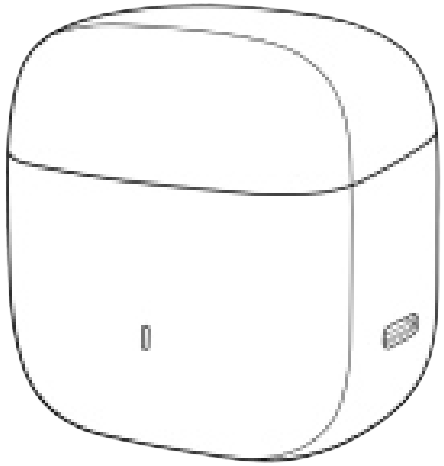
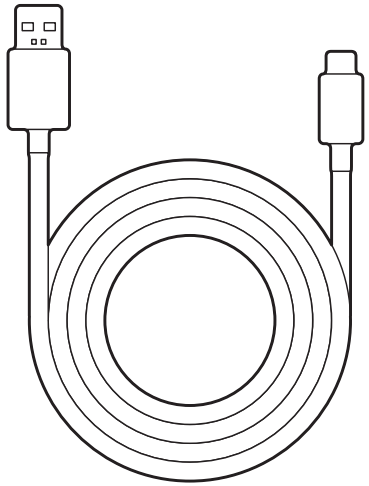
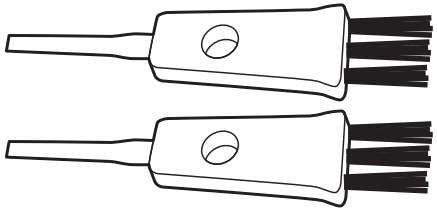
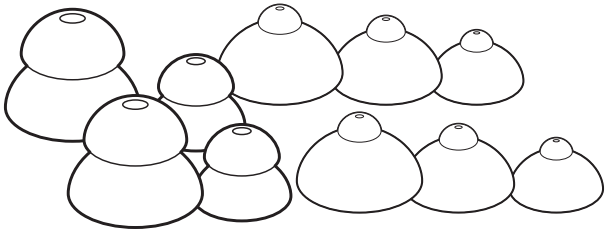
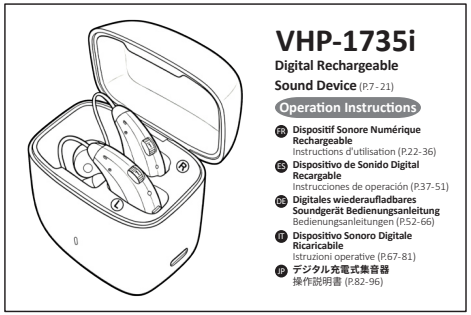
x 1

Eng Package Contents

Device with Eartip ×2
Charging Case ×1
Type-C Charging Cable x 1
Cleaning Brush x 2
Extra Eartips ×10
Operation Instruction x 1

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
First-Time Activation (Exit Shipping Mode)


Eng

To prevent battery drain during shipping, insulating tabs cover the charging contacts on both devices. Remove the tabs and activate each device before first use.

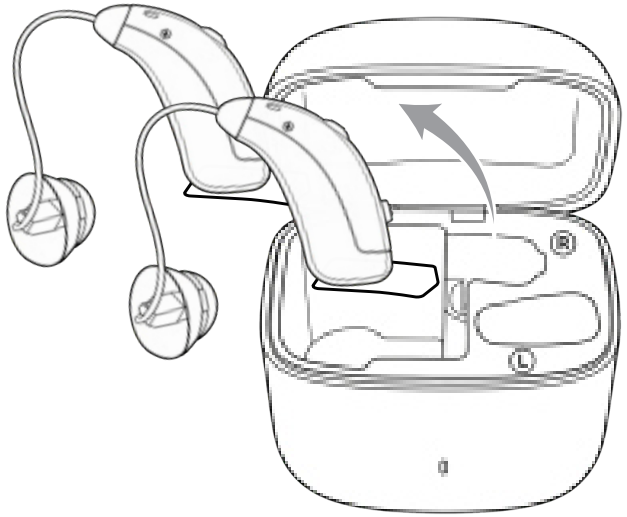
1. Remove Both Devices

Take both devices out of the charging case.

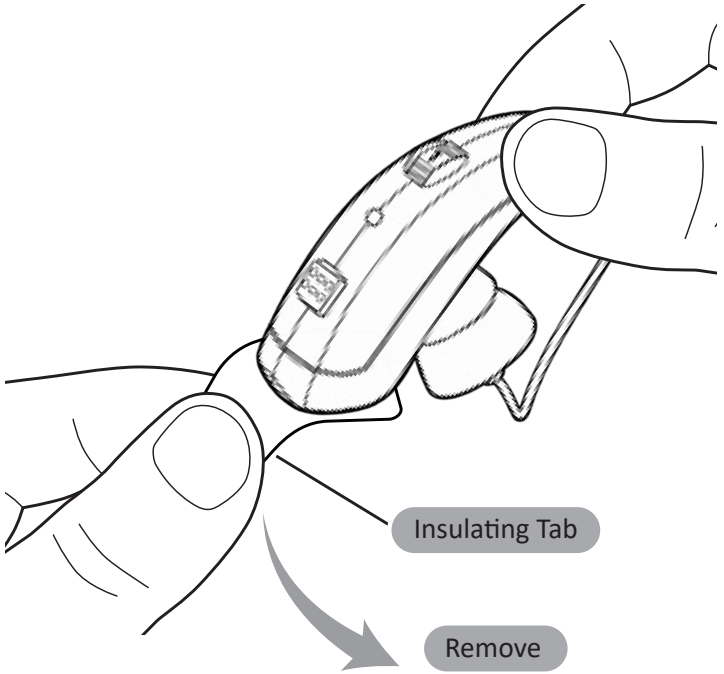
2.Remove the Insulating Tabs

Peel off the insulating tabs from the charging contacts on both devices.

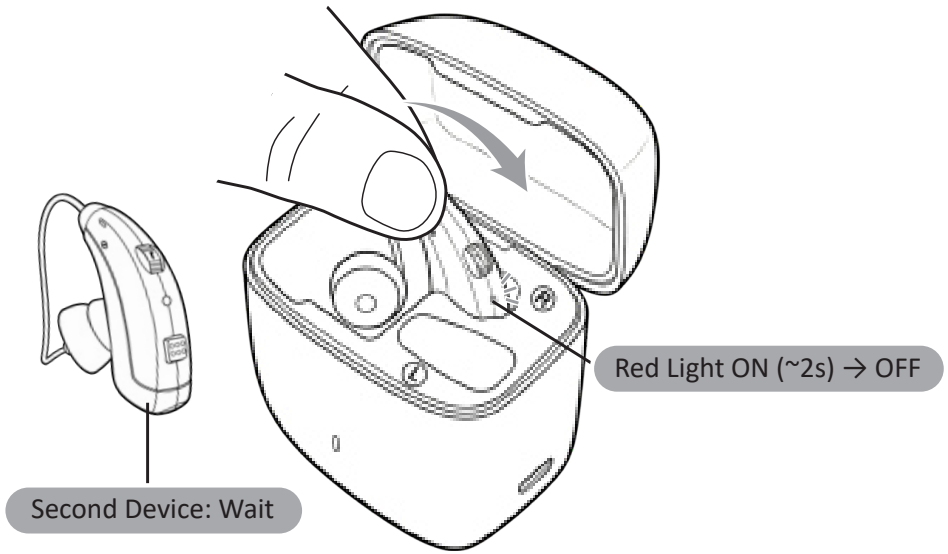3. Activate the First Device

Place one device into the charging case. A dim red light will turn on and then off after about 2 seconds. The device is now activated.

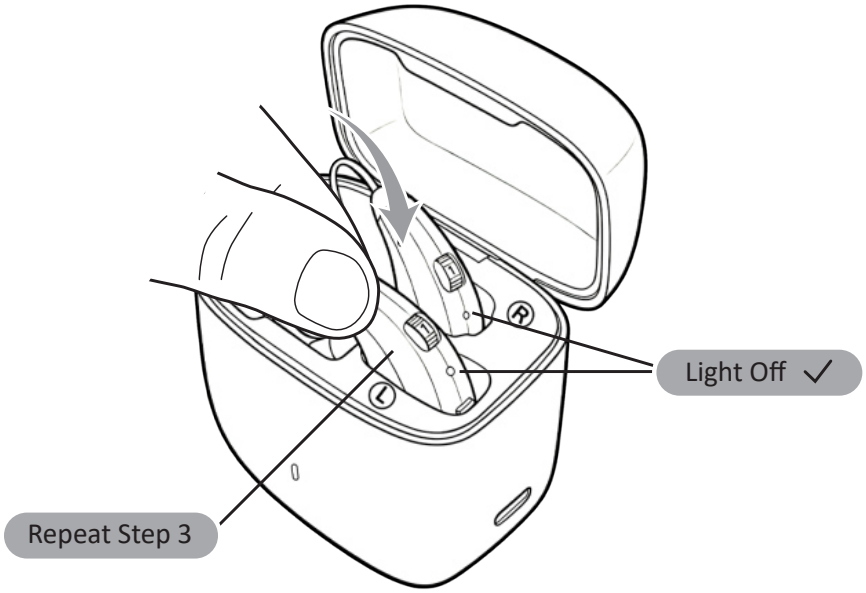4. Activate the Second Device


After the first device is activated, repeat Step 3 with the second device.

Important Notes

- **Do not activate both devices at the same time.** Activate them one at a time.
- If activation fails, remove the device and place it back into the charging case.
- Once both devices are activated, follow the manual to start using them.


First-time L / R Devices Syncing

Eng

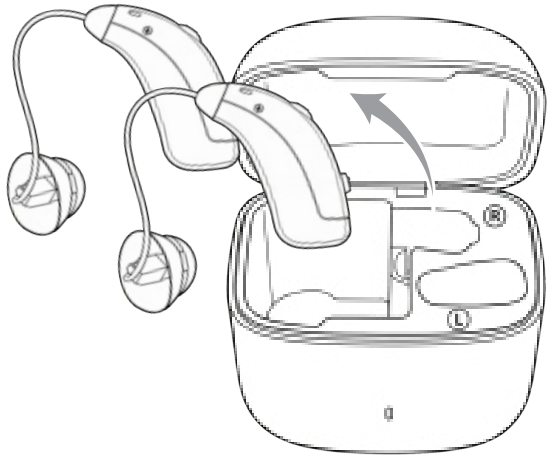
1.Power On Both Devices

Remove both the Left (L) and Right (R) devices from the charging case at the same time. The devices will power on automatically.

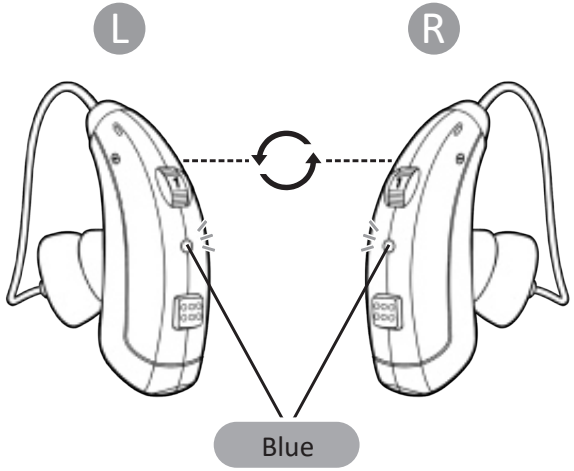2.Wait for Syncing

The blue indicators on both devices will flash while the Left and Right devices are syncing automatically.

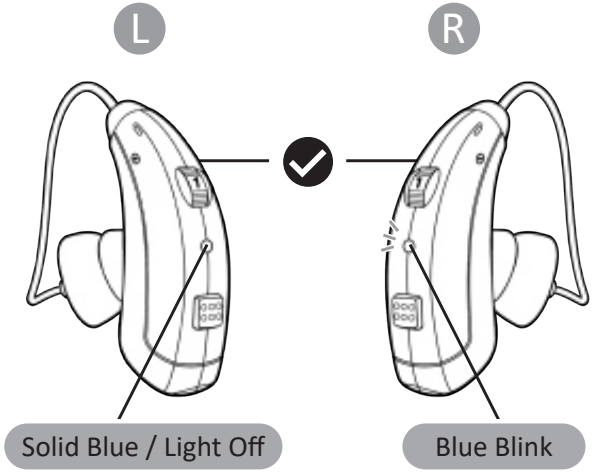
3.Syncing Complete

Syncing is complete when one device shows a solid blue light or the light turns off, while the blue indicator on the other device continues to flash.

Note: Syncing is not required for single-ear use. Simply power on one device.


Wireless Connection Operations

Eng

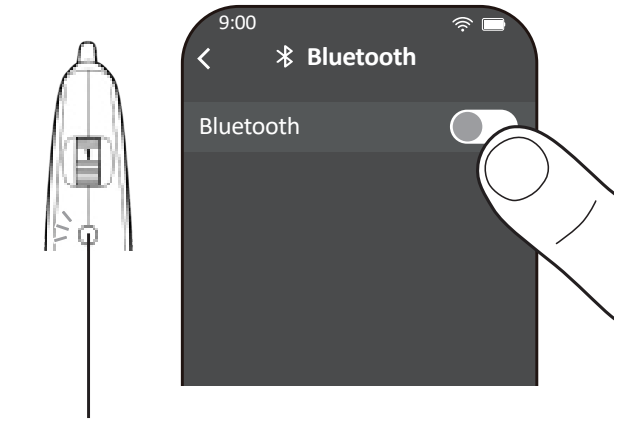
4.Turn On Wireless Connection

Open the wireless settings on your phone and turn on the wireless connection function.The blue indicator on the device will continue to flash while waiting for connection.

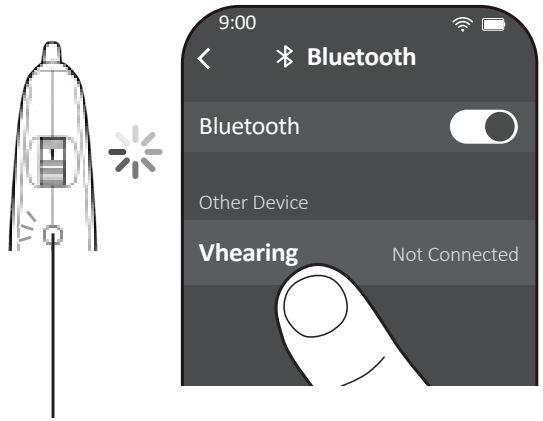5.Select “Vhearing”

Wait for “Vhearing” to appear in the available device list, then tap “Vhearing” to connect.

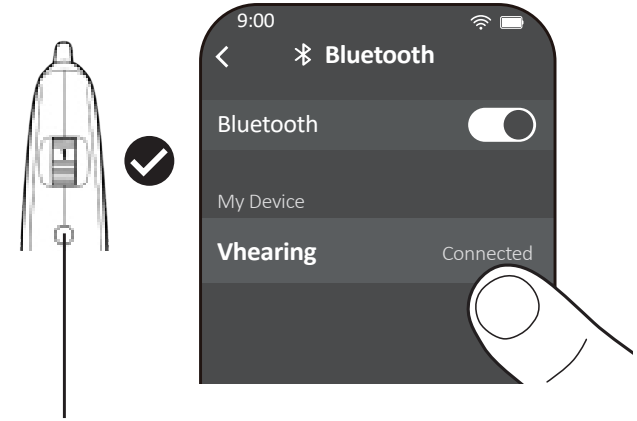6.Connection Complete

The wireless connection is successful when the blue indicator turns off.

A “CONNECTED” voice prompt will sound.

Note: Upon successful connection, you will hear a "Connected" voice prompt.
If the wireless connection is interrupted or turned off from phone, a "Disconnected" prompt will sound.

Connection Tips

1. Automatic Reconnection

After initial pairing, if “Vhearing” is saved on your phone, the devices will automatically reconnect when the phone’s wireless function is on.

2. Wireless Mode Control

Press and hold the Mode button for 3 seconds to turn off wireless connection and enter standard device mode.Press and hold it again for 15 seconds to turn wireless connection back on.

3. Switch Connected Device

To connect to another phone or device, first disconnect or remove “Vhearing” from the current device.
Place the devices in the charging case, then restart the connection process from Step 1.

4. Sleep Mode Before Connection

If not paired or connected within approximately 3 minutes, the wireless function will enter sleep mode.
Restart the devices to reconnect.

5. Sleep Mode After Connection

After connection, if no signal is received for approximately 6 minutes, the wireless function will enter sleep mode.Restart the devices to reconnect.

6. Auto Return After Playac


Approximately 6 seconds after wireless playback stops, the devices will automatically return to standard device mode.

7.Connection Range

Keep the connected device within approximately 6 m.
Walls, obstacles, and the wireless environment may affect connection performance.


Device Operation

Eng

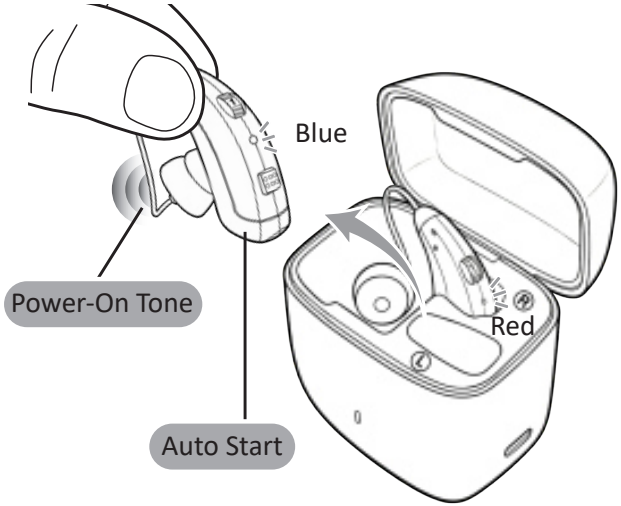1.Auto Power On


Take the device out of the charging case. It will automatically power on, the blue indicator light will turn on, and a power-on tone will sound.

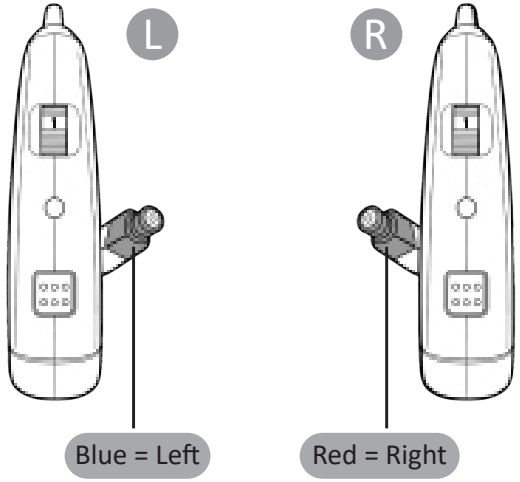2.Identify Left & Right


Check the color marking on the receiver before use.

Blue indicates the left ear (L), and red indicates the right ear (R).

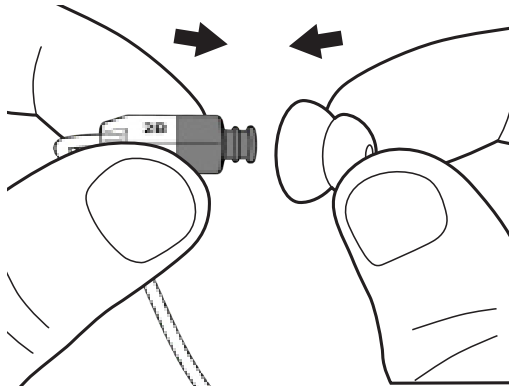3.Choose Your Eartip


Choose an eartip that fits comfortably. Replace the eartip as shown in Fig.3 if needed.

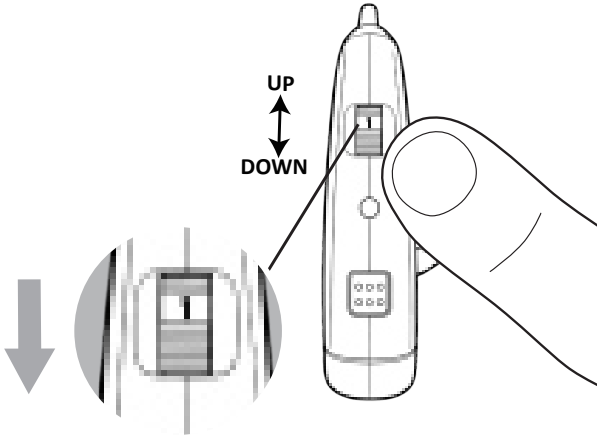4.Lower the Volume

Before wearing the device, adjust the volume to the lowest level.

Start at a low volume and increase it gradually after the device is properly fitted.

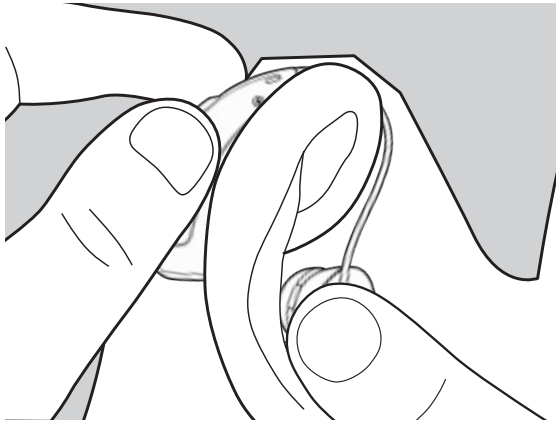
5. Wear the Device

Gently place the eartipin

your ear until it feels secure and comfortable. If whistling occurs,check the fit and reposition the

device.

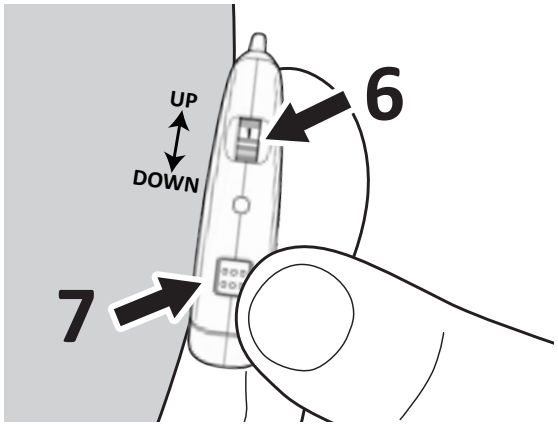
6.Adjust the Volume

Press the Volume Button to adjust the sound gradually until you reach a comfortable listening level.

7.Select the Mode

Press the Mode Button to switch between available listening modes. Select the mode that feels most comfortable for your current environment.

Mode Button Operation

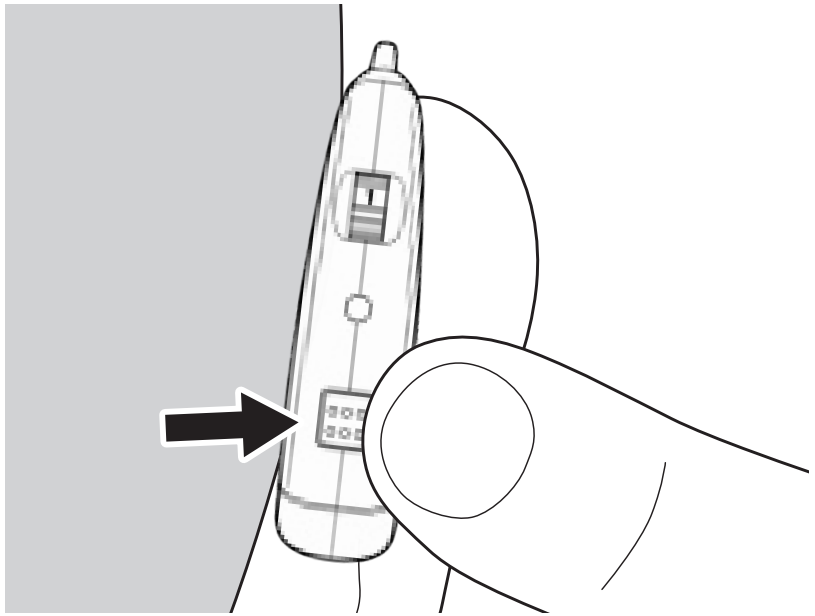Long Press — Wireless Connection

Press and hold the Mode button for 3 seconds to turn off wireless connection and enter Standard Device Mode.

A “POWER OFF” voice prompt will sound.

Press and hold the Mode button again for 15 seconds to turn wireless connection back on.

A “CONNECTED” voice prompt will sound.

Short Press — Mode Switch

Short press the Mode button to cycle through the 4 modes.

Mode 1 = 1 beep

Mode 2 = 2 beeps

Mode 3 = 3 beeps

Mode 4 = 4 beeps

Mode Guide


Eng

The device is equipped with 4 preset modes for different situations and environments. Press the "Ⓜ" button to switch between modes. The device will emit 1–4 beeps to indicate the selected mode.

Quiet Indoor Listening/Conversation Mode (MODE 1 = 1 beep)

Suitable for quiet indoor settings, conversations, or watching TV. This mode provides a balanced sound experience in quiet surroundings.

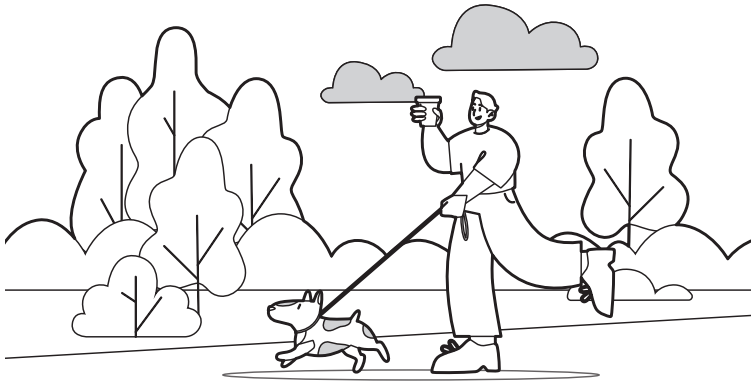
Relaxing Mode (MODE 2 = 2 beeps)

Suitable for quieter environments such as parks, reading areas, or libraries. This mode provides a comfortable sound experience in calm surroundings.

Noisy Social & Restaurant Mode (MODE 3 = 3 beeps)

Suitable for noisy social gatherings, restaurants, and group conversations. This mode helps reduce distracting background noise, such as surrounding chatter and dining sounds, for a more comfortable listening experience.

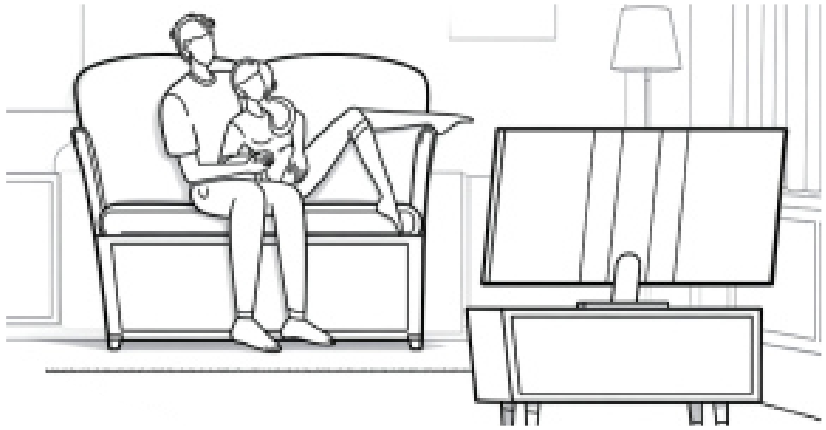TV & Video Mode (MODE 4 = 4 beeps)

Designed for watching television at home, including TV shows, movies, and news. This mode provides a comfortable listening experience, helping you enjoy your favorite programmes and follow conversations more easily.

Adaptive Noise Reduction


The device detects surrounding noise in real time and automatically reduces excessive background noise. This process takes approximately **9 seconds** after switching modes or turning on the device.

Note: *The scenarios above are for reference only. Select the mode based on your environment and preference.*

Charging Case Operation


Eng

Charging Your Devices with A **Portable Charging Case.**

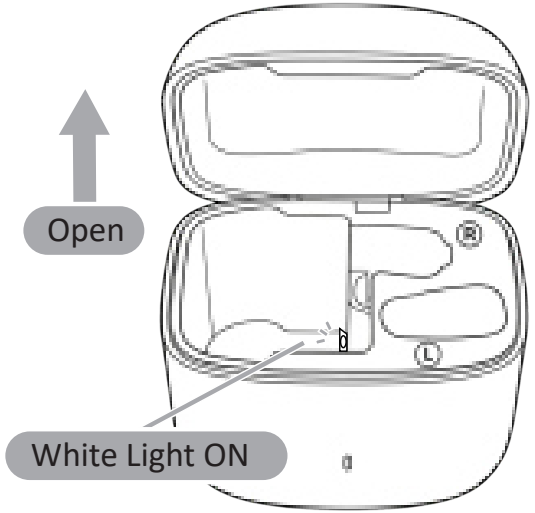1.Open the Charging Case

Open the case. The white light inside the case will turn on automatically.

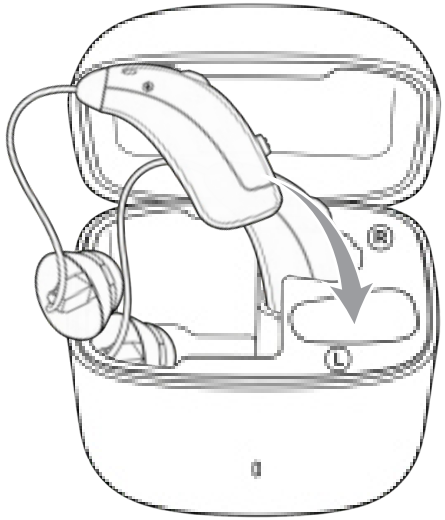
2.Insert the Devices

Place the devices into the charging case before charging. They will automatically power off once inserted.


3.Charge the Devices

The Red indicators remain on while the devices are charging. The Red indicators turn off when the devices are fully charged.

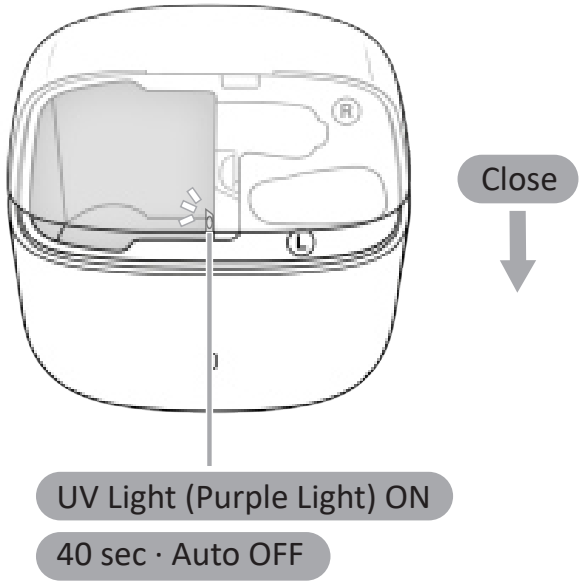4.Close the Charging Case


After closing the case, the white light turns off and the UV light (purple light) turns on automatically. The UV light operates for approximately 40 seconds and then turns off automatically.

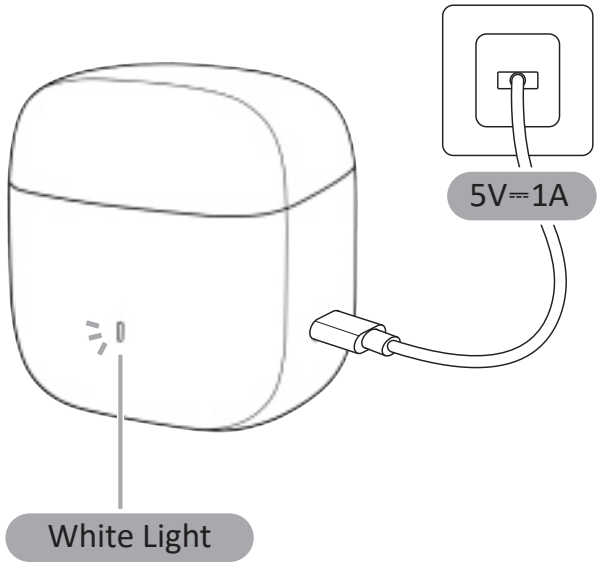5.Charge the Charging Case

Connect the charging case to a 5V⎓1A power source using the supplied Type-C cable.

Do not use fast-charging power sources or power banks.

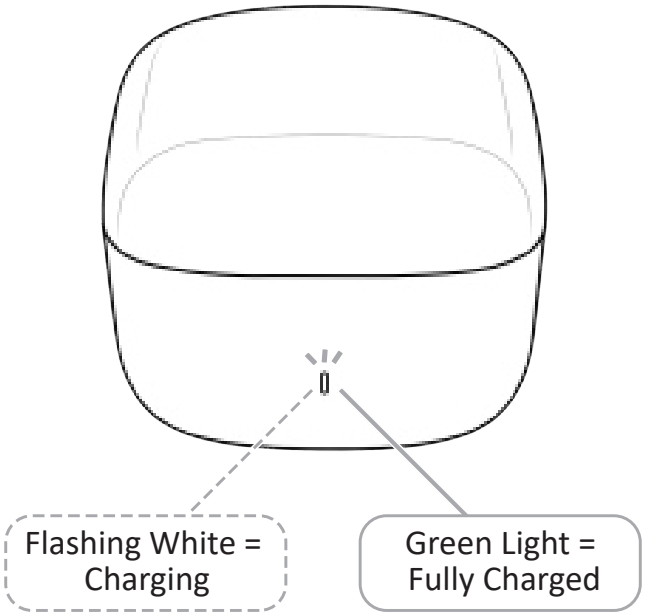6.Check the Charging Status

The Charging Case Status Indicator flashes white while the case is charging and turns green when fully charged.

A fully charged case can charge both devices up to 3 times.

Troubleshooting

Problem	Possible solution
Sound Device does not work	• Charge the device • Make sure the device is turned ON • Increase the volume
Sound is weak	• Increase the volume • Recharge the device • Clean or exchange the ear dome
Device emits howling & whistling sound	• If the device whistles while in the charging case, keep it in the case and charge it for at least 5 minutes. The whistling will stop automatically. • Reinsert the ear dome until it fits securely • Decrease the volume • Clean the ear dome • Use different ear dome size
Sound is distorted	• Decrease the volume • Clean or exchange the ear dome
Device not charging	• Clean the charge contact points with a cloth • Do not use portable charger or PowerBank to charge the device

Eng

Warning

- It is important to increase the volume control very gradually to avoid sudden increase in sound.
- Children can't use this product. If you have any problem, please consult audiologist.
- Please use the required adaptor to charge to prevent safety issue. Its parameters is

INPUT: 100V~240V, AC 50/60Hz, 0.2-0.25A
OUTPUT: 5V ⎓ 1A

Attention

1. This product is not waterproof, therefore the user must take care to prevent water from entering the product. If water enter the product, place the product in a warm place to dry.
2. Use it only when situations require amplification of sounds, such as movies, television and the like. it can't enlarge the extremely sounds.
3. Turn OFF the device when not use for a long time.
4. Do not use in or near an X ray or MRI machine.

Storage

Store in a dry place with good ventilation and related humidity no higher than 80%, and free of corrosive gas.

WARNING

CHOKING HAZARD - Small parts Not for children under 3 years.

FCC Compliance Statement

Eng

Warning: This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1.This device may not cause harmful interference.

2.This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC Class B Digital Device Statement

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy. If not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception (which can be determined by turning the equipment off and on), users are encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE: This device and its antenna(s) must not be co-located or operate in conjunction with any other antenna or transmitter.

The device has been evaluated to meet general RF exposure requirements. The device can be used under portable exposure conditions without restriction.

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